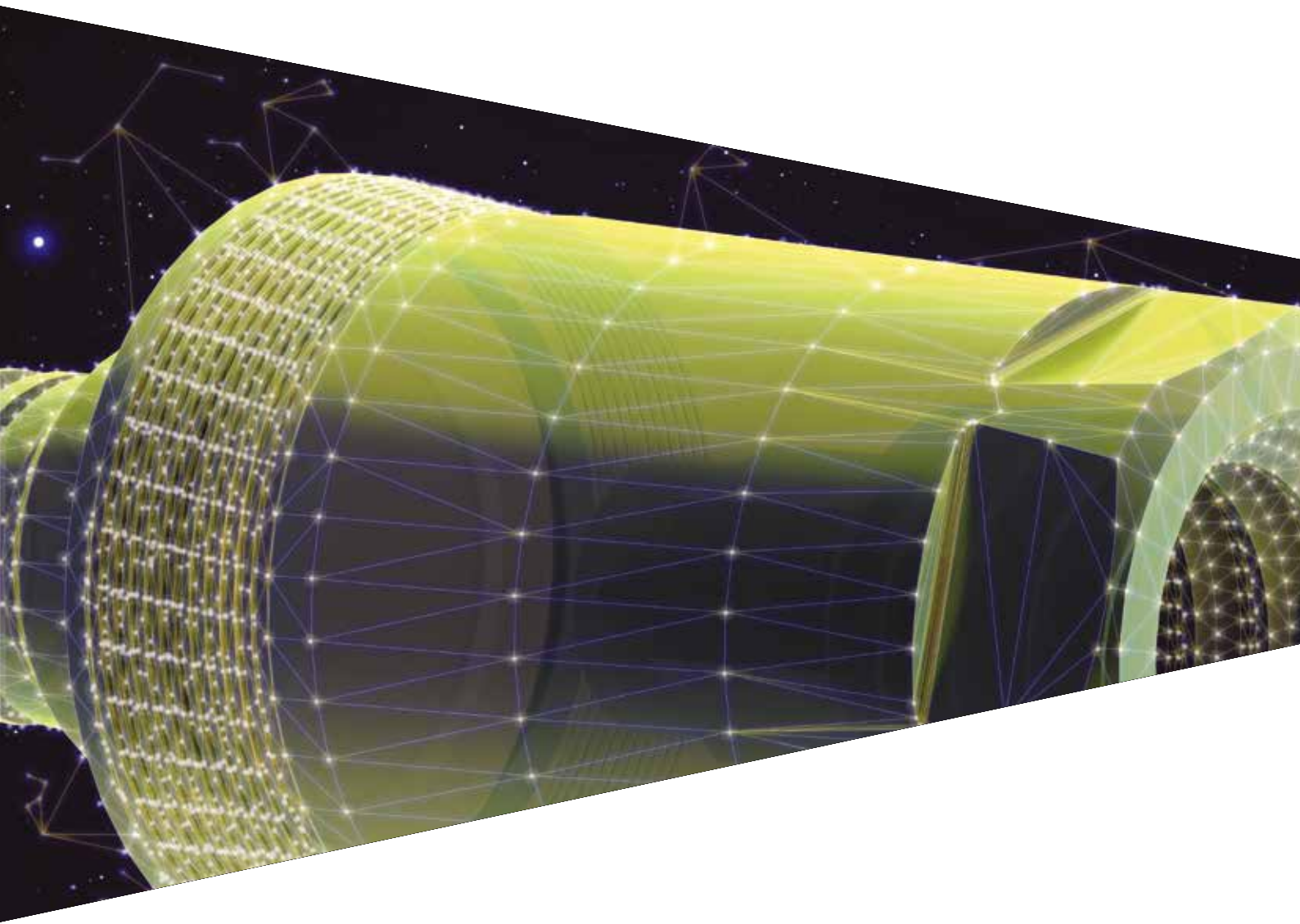


XA Prosthesis

For SHELTA implants



The principles of the XA concept

The tooth crosses an epithelium and the stabilization of the soft tissues defines its survival. The **soft tissues** around the tooth are **essential to maintain its function and the bone tissue** around it. Rehabilitations on implants are submitted to the same biological concepts related to the natural tooth. In order to obtain good functionality and preserve the tissues around the rehabilitation on implants, proper stabilization of soft tissues must be reached. The connective tissue guides the other tissues: the epithelium and the bone tissue.

XA prosthetic system was specifically meant to **guide and thicken the connective tissue**, improving its own biotype. The space obtained among the crown, the post, the soft tissue and the bone defines an area where a major quantity of collagen will be produced, thickening the tissues and preventing the bone loss (Chamber Concept, Degidi IJPRD 2013).



At the same time the presence of a **micro-threading** at the basis of the post **stimulates the production of collagen** (Contact Guidance Concept, Brunette IJOMI 1998 & Guillem Martí COIR 2012). The combination of the Platform Switching and the **conical structure** of the posts facilitates the **stabilization** of the circular fibers of the connective tissue at a more coronal level compared to a standard rehabilitation (Rodríguez, Vela IJOMI 2011; Rodríguez, Vela IJPRD 2016). Circular fibers stabilization is what determines a better bone level around the implant, which will improve with the passing of time. The peculiar feather-edge morphology of XA posts allows the possibility of inserting the post directly at the moment of exposure of the implant in the oral cavity, avoiding the unnecessary and noxious connections and disconnections during the prosthetic phase.

Clinical benefits in the use of XA posts are clear both from an aesthetic point of view, with **healthy and thick tissues**, and from the functional point of view, since the new bone tissue around the implant integrates itself permanently and **gives support and stability to the whole rehabilitation**.



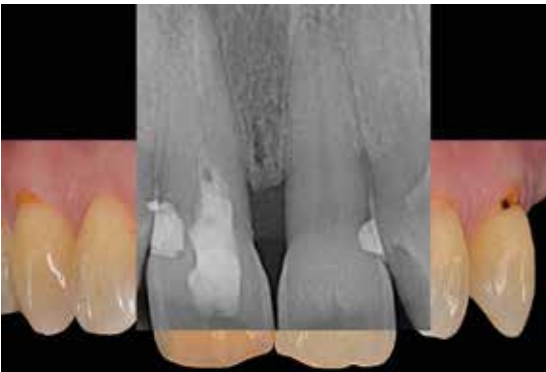
Clinical case

By courtesy of Dr. Xavier Vela Nebot, Dr. Xavier Rodríguez Ciurana and Mr. Javier Pérez López (laboratory technician)

The patient presents a failed of the endodontic retreatment at the incisor 1.1 and interproximal caries of the incisor 2.1. The substitution of incisor 1.1 with Shelta implant 4.25 mm and XA post for screw-retained prosthesis is planned, exploiting the micro-threading to optimize the aesthetic final result. Incisor 2.1 is rehabilitated using B.O.P.T. technique with a cap on natural post, prepared without any margin.



Frontal view: notice the gingival recession of both the incisors.



Initial case: radiographic and clinical occlusal view. The incisor 1.1 presents the failed of the endodontic retreatment, discouraging new treatments. The incisor 2.1 presents interproximal caries.



Preparation of the post on natural tooth 2.1.



Removal of the margin of the natural post with a conical incision up to the crestal bone level with the B.O.P.T. technique.



Extraction of the incisor 1.1.



Insertion of the Shelta implant with 4.25 mm in subcrestal position, exploiting the palatal wall of the socket to obtain a favorable emergence for the realization of a screw-retained prosthesis.



Insertion of the XA post for screw-retained prosthesis: the micro-threaded section remains at a subcrestal level.



Positioning of the Conoweld cap for conometric prosthesis.
Realization of a resin prosthesis incorporating the Conoweld cap and marking of the profile of the sulcus and the gingival margin.



Finalization of the temporary prosthesis: vestibular and palatal profiles (mesial and distal view).



Insertion of the graft material in the soft tissues and of a collagen sponge to preserve the horizontal dimension.



Double fastening of the temporary prosthesis with chlorhexadine gel on the natural post and with the Conoweld cap on the XA post.



Follow up after 4 months: the process of the gingival stabilization is promoted by the conical profiles of both the rehabilitations. The formation of the papilla is started, which is occupying progressively the space between the two crowns.



Detail of the healing of the tissues, where it can be noticed a unidirectional cells' replacement flow up to the coronal part. This phenomenon is the result of the change of the rehabilitation and it will gradually throw the surplus graft material out of the soft tissue.



Taking the impression on the XA post with its transfer, which avoids the disconnection of the post from the implant and facilitates the stabilization of soft and hard tissues.



Reproduction of the clinical situation in the gypsum model with removable gingiva.



Definitive rehabilitation in metal-ceramic for screwing on the XA post.



Comparison between the cap on natural tooth and the screw-retained rehabilitation on the XA post.



Follow up after 6 months: soft tissues appear healthy and continue their growth in coronal direction, they occupy progressively the interproximal space intended for the papilla. Gingival margins present the same trend of stabilization at a coronal level avoiding gingival recession.

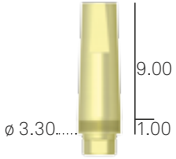
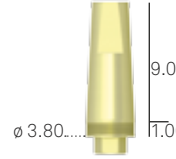
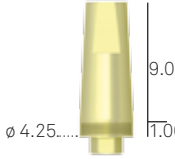
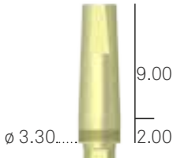
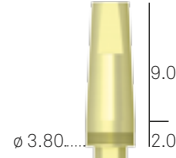
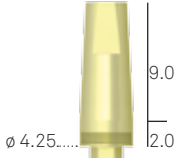


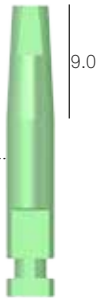

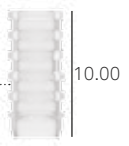
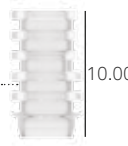
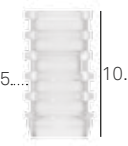
XA posts for cemented prosthesis

These posts, made of Gr. 5 titanium, allow cemented prosthesis to be fitted using the **One-Abutment-One-Time** technique, in which the XA post is left screwed into the patient's mouth, reproducing its form and position in the laboratory model with great precision, using a transfer and a dedicated analog.

Posts are available with diameters of 3.30, 3.80 and 4.25 to permit **Platform Switching** on the three platforms of the Shelta implant.



XA posts for cemented prosthesis

description	post \varnothing 3.30	post \varnothing 3.80 mm	post \varnothing 4.25
for implants	Premium 3.30 - 3.80 Kohno 3.80 Shelta 3.80	Premium 3.80 Kohno 3.80 Shelta 3.80 - 4.25 - 5.00	Shelta 4.25 - 5.00
Premade XA post Repositionable Transgingival H. 1.00 mm	SH-MD-F-330-1 	SH-MD-F-380-1 	SH-MD-F-425-1 
Premade XA post Repositionable Transgingival H. 2.00 mm	SH-MD-F-330-2 	SH-MD-F-380-2 	SH-MD-F-425-2 
Single pack Pack of 10 pieces Fixation screw with conical support	L-VMS-180 L-VMS-180-10 	Use L-VMS-180	Use L-VMS-180
Analog for premade XA posts for cemented prosthesis	SH-ANA-MD-F-330 	SH-ANA-MD-F-380 	SH-ANA-MD-F-425 
Transfer for premade XA posts for cemented prosthesis	SH-TRA-MD-F-330 	SH-TRA-MD-F-380 	SH-TRA-MD-F-425 

Recommended torque for definitive fastening of fixation screws: 20–25 Ncm.

IMPORTANT WARNING: The \varnothing 3.30 mm prosthetic components allow prosthetic Platform Switching with \varnothing 3.80 mm implants. It is recommended to use these posts exclusively for single crowns in front sectors (excluding premolars), and only as a support for multiple prostheses in distal sectors.

The \varnothing 3.80 mm prosthetic components are compatible with \varnothing 3.80 mm, \varnothing 4.25 mm and \varnothing 5.00 mm implants.

They do not allow prosthetic Platform Switching on \varnothing 3.80 mm implants; they allow prosthetic Platform Switching on \varnothing 4.25 mm and \varnothing 5.00 mm implants.

XA posts for screw-retained prosthesis

Intermediate XA abutments, with a feather-edge morphology, offer several prosthetic possibilities of exploiting the biological benefits of the XA concept, from temporary through to final prostheses. As with the line of cemented prostheses, abutments for screw-retained prostheses can be fitted using the **One-Abutment-One-Time** technique, leaving the abutments screwed into the patient's mouth and reproducing their form and position in the laboratory model using a suitable analog.

The only available diameter of 3.80 mm permits the greatest possible prosthetic simplification, and ensures Platform Switching, which is extremely favourable in reducing peri-implant bone reabsorption to a minimum. On posts with heights of 5.50 and 6.50 mm, an extremely safe **conometric prosthesis** can be fitted, by using Conoweld caps.

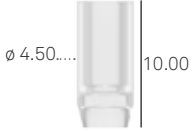

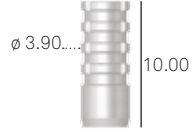
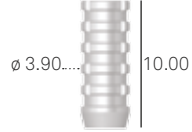




XA posts for screw-retained prosthesis

∅ prosthetic component	∅ 3.80	
for implants	Premium 3.80 Kohno 3.80 Shelta 3.80 - 4.25 - 5.00	
Intermediate XA abutment H. 4.50 mm	SH-ABU-F-TS-380-4	
Intermediate XA abutment H. 5.50 mm	SH-ABU-F-TS-380-5	
Intermediate XA abutment H. 6.50 mm	SH-ABU-F-TS-380-6	
Analog for intermediate XA abutments	SH-ANABU-F-380	
Transfer for intermediate XA abutments	SH-TRABU-F-380	
Transfer screw	SH-VTRABU-F-200	
Healing cap in PEEK for intermediate XA abutments	SH-CG-ABU-F-380	

Recommended torque for definitive fastening of directly screwed posts: 25–30 Ncm.





Accessories for screw-retained prosthesis on XA posts

∅ prosthetic component	∅ 3.90 - 4.50
for implants	Premium 3.80 Kohno 3.80 Shelta 3.80 -4.25- 5.00
Castable sleeve for XA abutments Repositionable Fixation screw included	SH-CCABU-F-380  ∅ 4.50... 10.00
Castable sleeve for XA abutments Non-repositionable Fixation screw included	SH-CCABU-F-380-ROT  ∅ 4.50... 10.00
Titanium sleeve for XA abutments Repositionable Fixation screw included	SH-CTABU-F-380  ∅ 3.90... 10.00
Titanium sleeve for XA abutments Non-repositionable Fixation screw included	SH-CTABU-F-380-ROT  ∅ 3.90... 10.00
Fixation screw for XA abutments	A-PLAIN-VP200  M2.0...
Final cap for Conoweld luting*	CAP-TS-DEF 

Recommended torque for definitive fastening of directly screwed structures on posts: 20–25 Ncm.

* the conometric cap CAP-TS-DEF can be used only with posts with H. 5.50 mm and 6.50 mm.

Prosthesis screwdrivers

description	code
Screwdriver for fixation screws, with hexagonal connector for ratchet or hand knob, short	HSM-20-EX  A short screwdriver with a hexagonal connector at the top and a 1.25 mm tip. The text "1.25 MM" is visible on the shaft.
Screwdriver for fixation screws, with hexagonal connector for ratchet or hand knob, long	HSML-20-EX  A medium-length screwdriver with a hexagonal connector at the top and a 1.25 mm tip. The text "1.25MM-L" is visible on the shaft.
Screwdriver for fixation screws, with hexagonal connector for ratchet or hand knob, extra-long	HSMXL-20-EX  An extra-long screwdriver with a hexagonal connector at the top and a 1.25 mm tip. The text "1.25 MM-XL" is visible on the shaft.
Driver for fixation screws, with right angle shank	HSM-20-CA  A screwdriver with a right-angle shank and a 1.25 mm tip. The text "HSM-20-CA" is visible on the shaft.

The clinical and scientific evolution behind the development of the basic concept of XA prosthesis



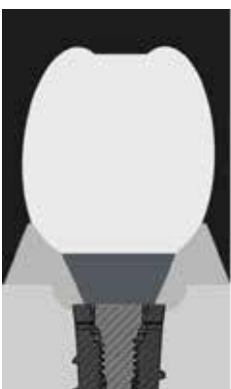
The main preoccupation of the practitioner is to limit the inexorable bone reabsorption and the apical migration of soft tissues to a minimum, as these generate many problems that are not solely aesthetic, but that can instead also be linked to the medium-term and long-term survival of the rehabilitation, due to peri-implantitis.

What we are talking about is a paradigm shift, given that with the use of **XA posts** and an adequate rehabilitation protocol, not only can apical migration be avoided, but coronal migration of the bone and of soft tissues can also be promoted over time. This implies both aesthetic improvements and greater protection against the onset of peri-implantitis.

We believe that the best way of demonstrating the evolution of the XA concept is to analyze the various protocols that we have adopted in our work in recent years, with the final aim of stabilizing tissues.



1. Implant without Platform Switching (PS) presenting an anatomical profile (divergent) that mimics the form of the tooth, and a protocol that envisions prosthesis disconnections. In this situation, tissues showed a tendency to **apical migration**. The inevitable gingival recession generated not only an aesthetic risk, but also the possibility of peri-implantitis.



2. Implant with PS presenting an anatomical profile (divergent) and a protocol that envisages prosthesis disconnections. The introduction of the concept of PS has allowed us to reduce peri-implant bone reabsorption, improving **tissue stability**. This concept allowed us to learn more about the relationship between the shape of posts and the stabilization of tissues by means of the circular fibres of connective tissue.



3. Implant with PS presenting a straight profile and a protocol that envisages prosthesis disconnections. The use of straight posts allowed us to increase tissue thickness and to improve results in a **predictable** way.



4. Implant with PS presenting a conical profile (convergent) of the XA posts and absence of disconnections. Finally, the introduction of **XA posts** for screw-retained and cemented prostheses allowed us not only to stabilize tissues, but also to promote their **coronal migration** over time.

The **tapering** of these posts promotes the short-term and long-term coronal migration of tissues. The absence of a margin (**marginless**) makes it possible to avoid the repeated disconnections and reconnections that are typical of conventional rehabilitation protocols, causing tissue destabilization.

The presence of **micro-threads** at the base of the posts promotes the alignment of fibroblasts, accelerating and increasing the production of collagen (contact guidance), a crucial factor for biological sealing.

The apparent **simplicity** of this concept satisfies a multitude of biological criteria that contributes in the search for the best possible tissue response. The use of the components is simple and facilitates the work of both the dentist and the dental technician. We took care to ensure that the posts were compatible with the different operating protocols (screw-retained or cemented prostheses), and that they could be used in all clinical situations (anterior and posterior sectors, immediate or deferred loading). For the first time, posts act as a guide for connective tissue and for the periosteum, in turn promoting bone growth. All these reasons lead us to believe that we are truly faced with a paradigm shift.

Insanity?

The definition of insanity is doing the same thing over and over and expecting different results.

A. Einstein

Dr. Xavier Vela Nebot
Dr. Xavier Rodríguez Ciurana

Centro BORG (Barcelona Osseointegration Research Group)

Bibliography on the XA principle

Vela X., Rodríguez X., Rodado C., Segalá M.

Benefits of an Implant Platform Modification Technique to Reduce Crestal Bone Resorption
Implant Dentistry / Volume 15, Number 3 2006, 313-320

Steigmann M., Monje A., Chan H., Wang H.

Emergence Profile Design Based on Implant Position in the Esthetic Zone
Int J Periodontics Restorative Dent 2014; 34: 559-563

Schoenbaum T.R., Chang Y., Klokkevold P.R., Snowden J.S.

Abutment Emergence Modification for Immediate Implant Provisional Restoration
Journal of Esthetic and Restorative Dentistry Vol 25, N.2, 103-107, 2013

Zucchelli G., Mazzotti C., Mounssif I., Marzadori M., Stefanini M.

Esthetic Treatment of Peri-implant Soft issue Defects: A Case Report of a Modified Surgical-Prosthetic Approach
Int J Periodontics Restorative Dent 2013, 33:327-335

Vela X., Méndez V., Rodríguez X., Segalá M., Tarnow D.

Crestal Bone Changes on Platform-Switched Implants and Adjacent Teeth When the Tooth-Implant Distance is Less Than 1.5 mm
Int J Periodontics Restorative Dent 2012; 32: 149-155

Guillem-Martí J., Delgado L., Pegueroles M., Herrero M., Gil F.J.

Fibroblast adhesion and activation onto micro-machined titanium surfaces
Clin. Oral Impl. Res. 00, 2012, 1-11

Caram S.J., Huynh-Ba G., Schoolfield J.D., Jones A.A., Cochran D.L., Belser U.C.

Biologic Width Around Different Configurations. A Radiographic Evaluation of the Effect of Horizontal Offset and Concave Abutment Profile in the Canine Mandible
Int J Oral Maxillofac Implants. 2014 Sep-Oct;29(5):1114-22.

Rodríguez X., Vela X., Méndez V., Segalà M., Calvo-Guirado J.L., Tarnow D.P.

The effect of abutment dis/reconnections on peri-implant bone resorption: a radiologic study of platform-switched and non-platform-switched implants placed in animals
Clin Oral Implants Res. 2013 Mar;24(3):305-11. Epub 2011 Oct 3.

Brunette D.M.

The effects of implant surface topography on the behavior of cells
Int J Oral Maxillofac Implants. 1988 Winter;3(4):231-46.

Romanos G.E.

Tissue preservation strategies for fostering long-term soft and hard tissue stability
Int J Periodontics Restorative Dent. 2015 May-Jun;35(3):363-71.

Rodríguez X., Vela X., Calvo-Guirado J.L., Nart J., Stappert C.F.

Effect of platform switching on collagen fiber orientation and bone resorption around dental implants: a preliminary histologic animal study
Int J Oral Maxillofac Implants. 2012 Sep-Oct;27(5):1116-22.





rev.09-16



Sweden & Martina S.p.A.

Via Veneto, 10
35020 Due Carrare (PD), Italy
Tel. +39.049.9124300
Fax +39.049.9124290
info@sweden-martina.com
www.sweden-martina.com

Sweden & Martina Mediterranea S.L.

Sorolla Center, Oficina 801
Avda Cortes Valencianas 58, 8pl
46015-Valencia, Spain
Tel. +34.96.3525895
Tel. 900993963
info.es@sweden-martina.com

Sweden & Martina Ltd

Basepoint Business Centre Crawley
Metcalf Way, Crawley,
West Sussex, RH11 7XX, UK
info.uk@sweden-martina.com

Sweden & Martina Inc.

c/o DCI Management
301 Pleasant St.
Abbottstown, 17301 PA, US
Toll free 844-8MARTINA
844-862-7846
info.us@sweden-martina.com

The implants, standard prosthetic components and surgical instruments contained in this catalogue are Medical devices and are manufactured by Sweden & Martina S.p.A. They conform to the ISO 9001 and ISO 13485 standards and are certified with the CE Mark (Class I) and CE 0476 mark (Class IIA and class IIB) in compliance with European Medical Device Directive No. 93/42 and European Directive No. 2007/47/CE.

We have met the good manufacturing standards (GMP) set forth by many countries worldwide, including the United States FDA.



Some products may not be regulatory/released for sale in all markets.
All trademarks herein are the property of Sweden & Martina S.p.A. unless otherwise indicated.
This material is intended for laboratories and clinicians and is not intended for patient distribution.
This material is not to be redistributed, duplicated, or disclosed without the express written consent of Sweden & Martina S.p.A.
For additional product information, including indications, contraindications, warnings, precautions, and potential adverse effects, see Sweden & Martina S.p.A. website.

The contents are updated at the time of publication. Check with the company for any subsequent updates.