


PREMIUM - SHELTA **POWER**

International Product Overview



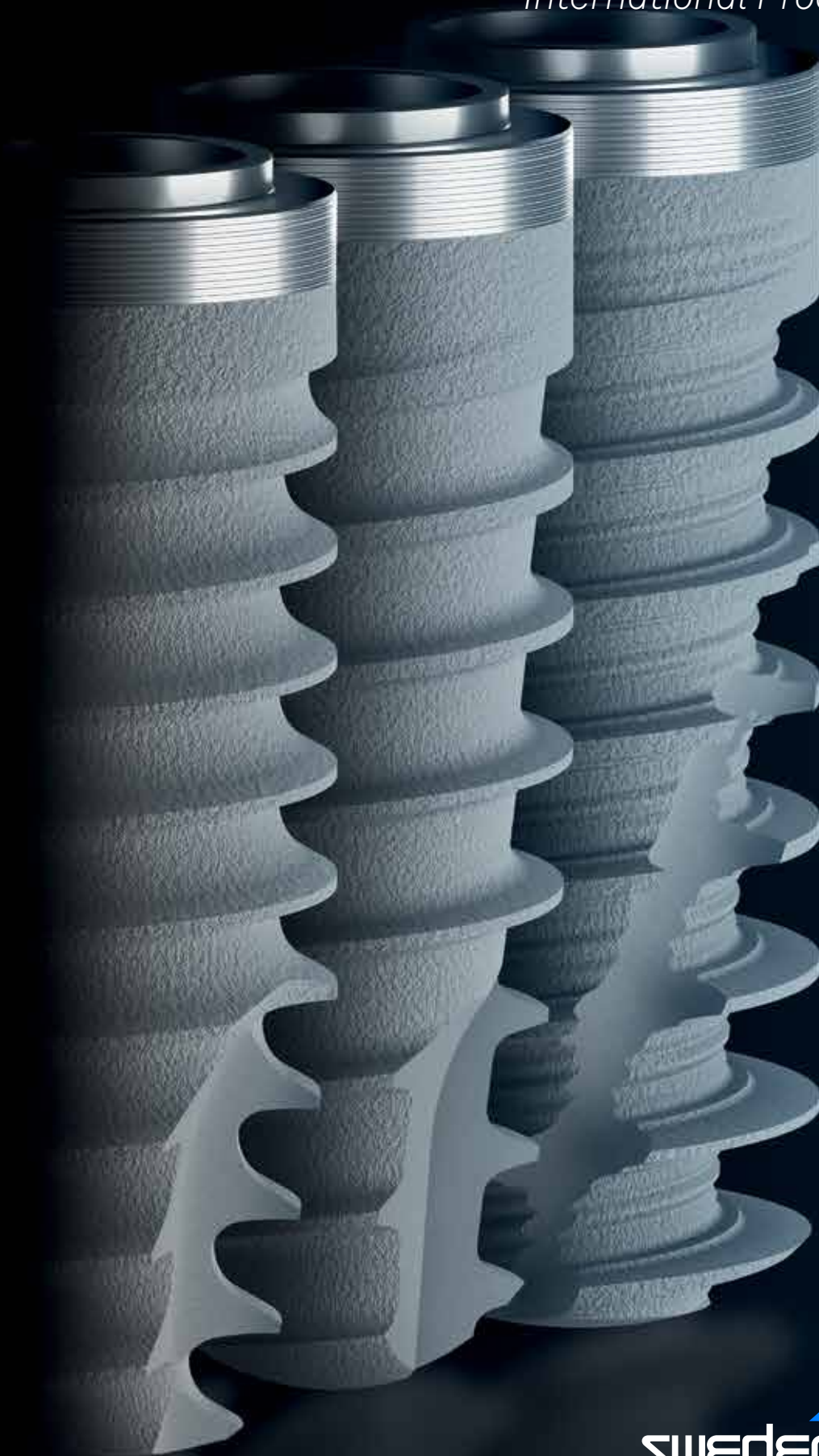

Sweden & Martina











Dental Arena

PREMIUM POWER AND SHELTA POWER












List of contents

4	Table of Premium Power and Shelta Power implants
6	Premium Power Slim and Premium Power
8	Shelta Power and Shelta SL Power
10	UTM and ZirTi: two surfaces in synergy
12	The advantages of Collex connection
14	Digital Workflow with Premium Power and Shelta Power
16	Ergonomics of Surgical Procedures
18	Surgical kits
24	Surgical kits for guided surgery
28	Complements for guided surgery
28	Optional surgical instruments
31	Depth notches
32	Prosthetic solutions for excellence
34	Implant platforms and prosthetic combinations
35	Healing abutments
36	REGISTRA Technique
37	Impression taking and model phase
38	Temporary posts
39	Preformed posts
40	Posts and components for vertical technique
41	Millable posts
42	XA prosthetic line
44	Castable posts
45	Solutions with angled screw hole
46	Prosthetic line PAD ^e (Aesthetic screw-retained disparallel prosthesis)
48	PLAIN B-Space prosthetic line
49	Conico Conometric prosthesis
50	Components for digital prosthesis
55	Components for chairside digital prosthesis
56	Solutions for overdenture
59	Clinical and laboratory prosthetic screws
60	Implant site preparation sequences

Table of Premium Power and Shelta Power implants

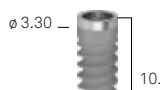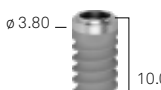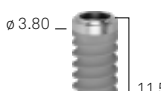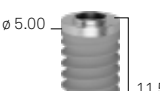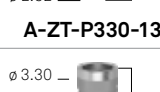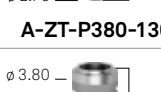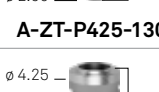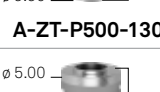
	Premium Slim*	Premium Power		
h fixture	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
7.00 mm	-	 A-ZT-P380-070	 A-ZT-P425-070	 A-ZT-P500-070
8.50 mm	 A-ZT-P330-085	 A-ZT-P380-085	 A-ZT-P425-085	 A-ZT-P500-085
10.00 mm	 A-ZT-P330-100	 A-ZT-P380-100	 A-ZT-P425-100	 A-ZT-P500-100
11.50 mm	 A-ZT-P330-115	 A-ZT-P380-115	 A-ZT-P425-115	 A-ZT-P500-115
13.00 mm	 A-ZT-P330-130	 A-ZT-P380-130	 A-ZT-P425-130	 A-ZT-P500-130
15.00 mm	 A-ZT-P330-150	 A-ZT-P380-150	 A-ZT-P425-150	 A-ZT-P500-150
18.00 mm	-	 A-ZT-P380-180	 A-ZT-P425-180	-
surgical cover screws	 A-VT-330	 A-VT-380	 SH-VT-425-BL	 SH-VT-500-VI

* Sweden & Martina dental implants are implantable medical devices indicated for the rehabilitation of patients with total or partial edentulism. They are intended for surgical placement in the mandibular or maxillary bone. Sweden & Martina implant systems are designed for both one-stage and two-stage surgical procedures. Due to their diameter, implants with a diameter of 3.30 mm or less are recommended in cases where bone thickness is limited. They may be used for the rehabilitation of single crowns replacing incisors and canines. Owing to their reduced dimensions, they are not suitable for the rehabilitation of single crowns in premolar or molar positions. In posterior regions, these implants may only be used in bar- or bridge-supported rehabilitations, in combination with larger-diameter implants.

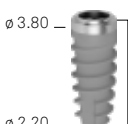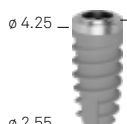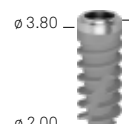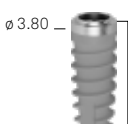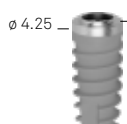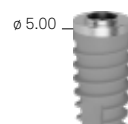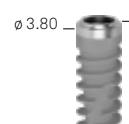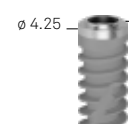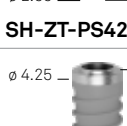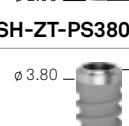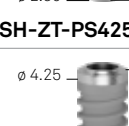
	Shelta Power			Shelta Power SL		
	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
	-	-	-	-	-	-
	 SH-ZT-PS380-085	 SH-ZT-PS425-085	 SH-ZT-PS500-085	 SH-ZT-PS380SL-085	 SH-ZT-PS425SL-085	 SH-ZT-PS500SL-085
	 SH-ZT-PS380-100	 SH-ZT-PS425-100	 SH-ZT-PS500-100	 SH-ZT-PS380SL-100	 SH-ZT-PS425SL-100	 SH-ZT-PS500SL-100
	 SH-ZT-PS380-115	 SH-ZT-PS425-115	 SH-ZT-PS500-115	 SH-ZT-PS380SL-115	 SH-ZT-PS425SL-115	 SH-ZT-PS500SL-115
	 SH-ZT-PS380-130	 SH-ZT-PS425-130	 SH-ZT-PS500-130	 SH-ZT-PS380SL-130	 SH-ZT-PS425SL-130	 SH-ZT-PS500SL-130
	 SH-ZT-PS380-150	 SH-ZT-PS425-150	 SH-ZT-PS500-150	 SH-ZT-PS380SL-150	 SH-ZT-PS425SL-150	 SH-ZT-PS500SL-150
	-	-	-	-	-	-
	 SH-VT-380-VE	 SH-VT-425-BL	 SH-VT-500-VI	 SH-VT-380-VE	 SH-VT-425-BL	 SH-VT-500-VI

Premium Power


SLIM Ø 3.30

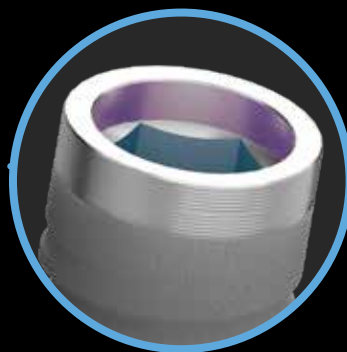


**COLLEX ONE
CONNECTION**
For a symplified
prosthetic approach

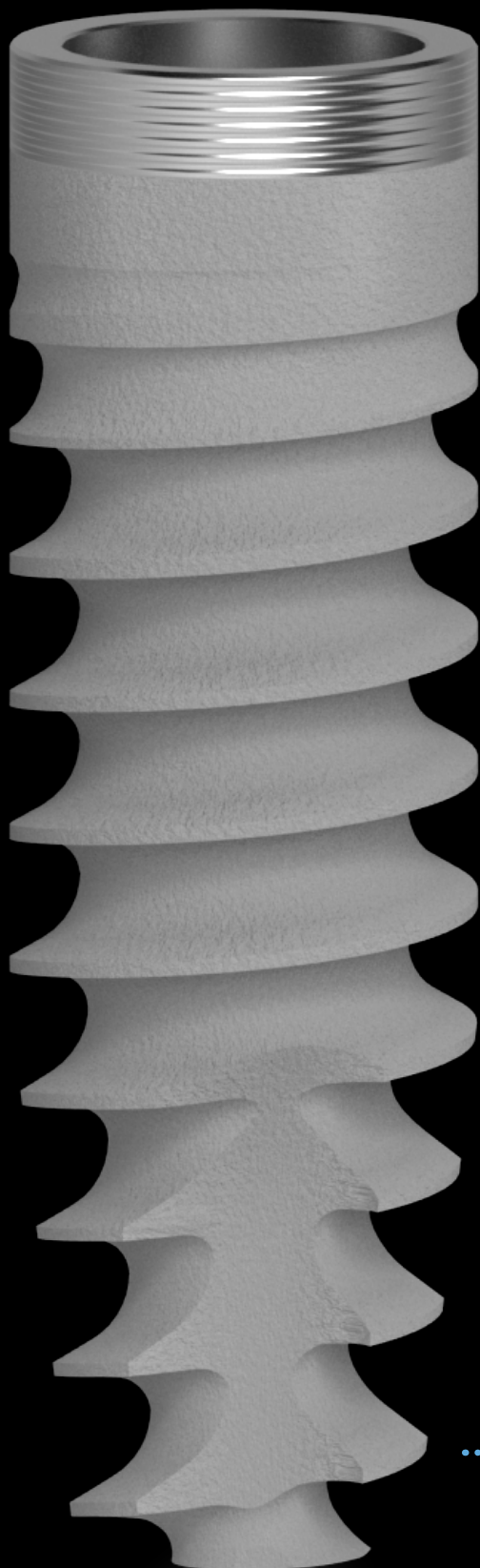


**PREMIUM SLIM POWER:
CORONAL MORPHOLOGY**
More stable in thin cortical bone

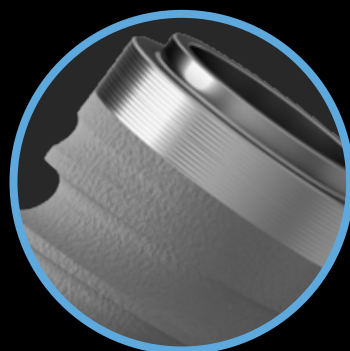

SLIM MORPHOLOGY
High performance core in
every bone type, to save grip in
reduced space

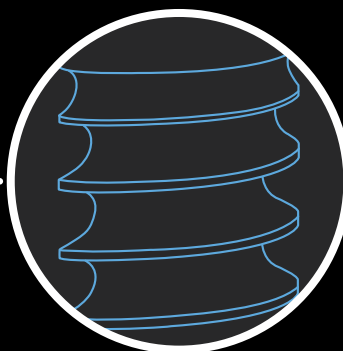


EASY INSERTION
Thanks to the
core/thread ratio,
specially studied

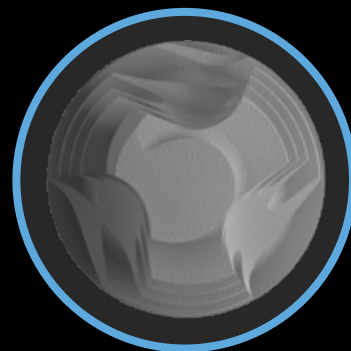


Premium Power

Ø 3.80, Ø 4.25, Ø 5.00

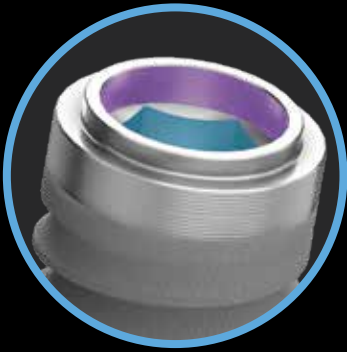


UTM NECK

The surface meant for the interface area between hard and soft tissues

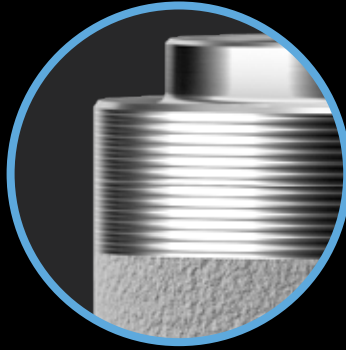


PREMIUM POWER: CORONAL MORPHOLOGY

Excellent performance even in hard cortical bone

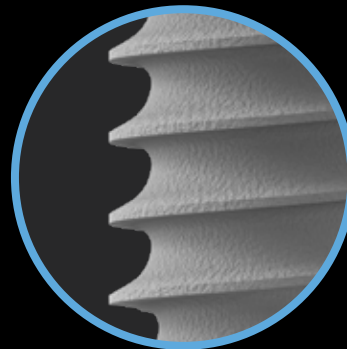


REVERSE BUTTRESS THREAD

High cutting capacity and penetration even in hard bone

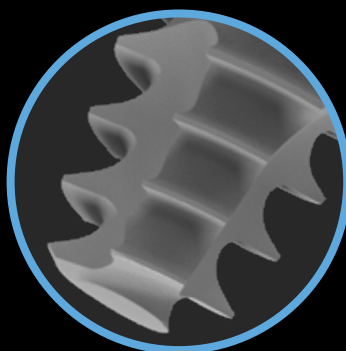


APICAL NOTCHES

Designed to facilitate insertion and clot decompression

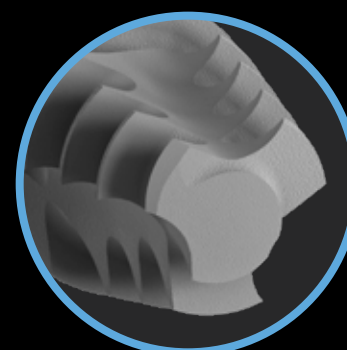


APEX

Tapered, with deeper thread throughout the length, it assures stability mainly in post-extraction cases

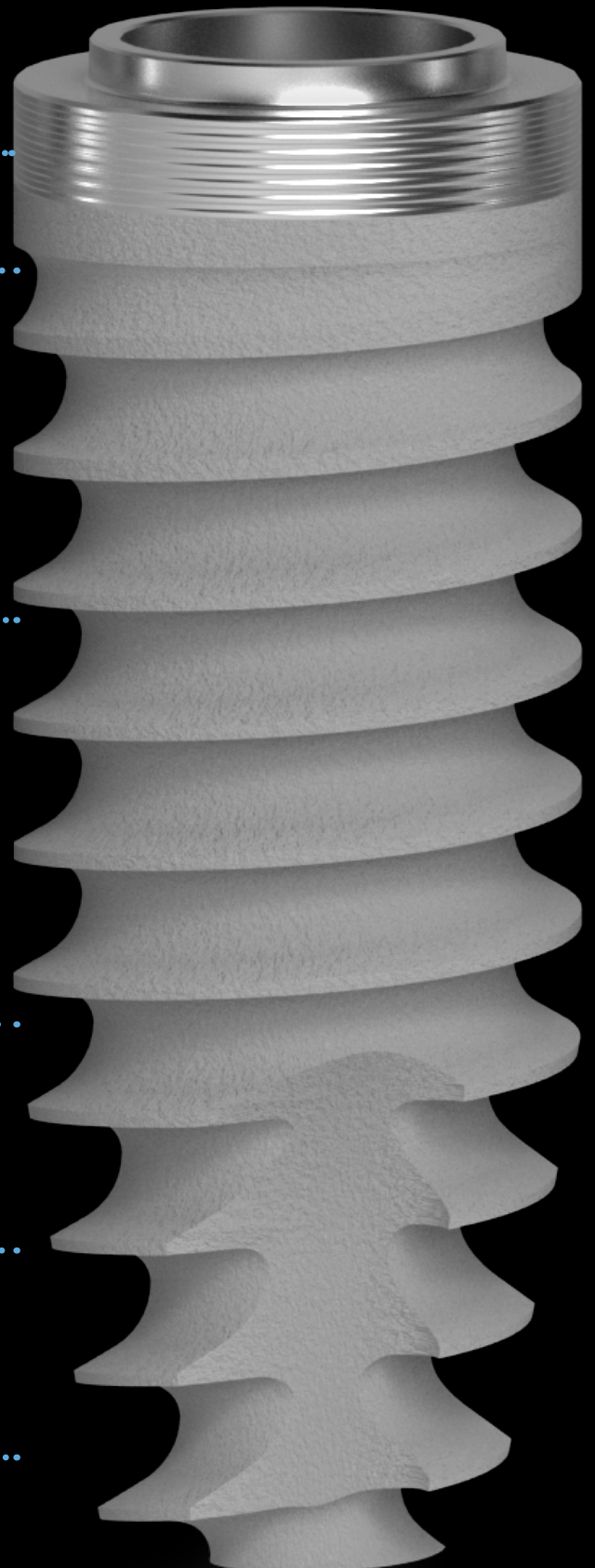

Shelta Power

Ø 3.80, Ø 4.25, Ø 5.00

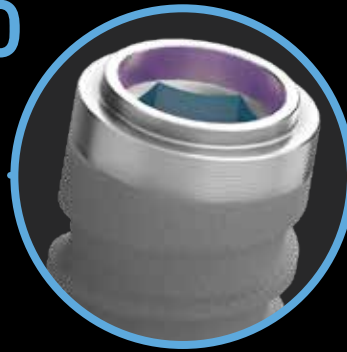


CONNECTION COLLEX ONE

common to all diameters for a simplified prosthetic approach



CORONAL MORPHOLOGY

which makes insertion precise and increases primary stability in the upper jaw

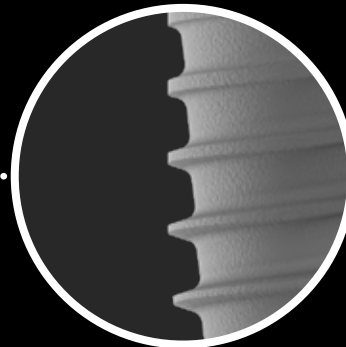


CONICAL CORE AND CONICAL THREAD

allows excellent sensitivity during insertion



THREADING


with progressive thickening to increase the condensation of bone poorly mineralised

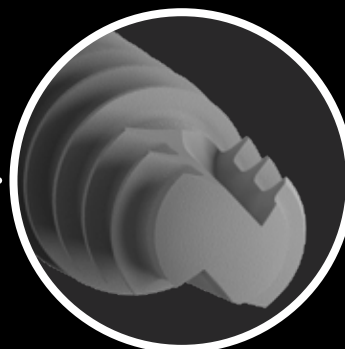

ROUNDED APEX


perfect for sinus lift

Ø 3.80, Ø 4.25, Ø 5.00

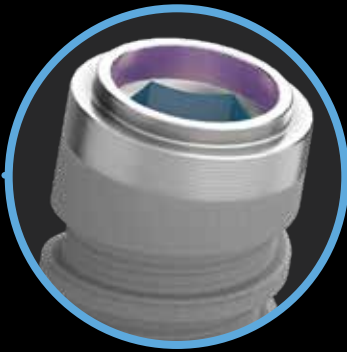


UTM NECK

The surface designed for the interface area between soft tissue and hard tissue

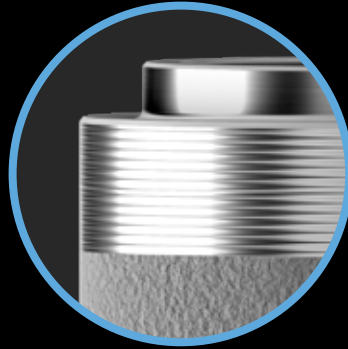


LARGE THREAD AND GROOVES IN THE CORE

for a greater primary stability to increase the contact surface when the bone is less dense

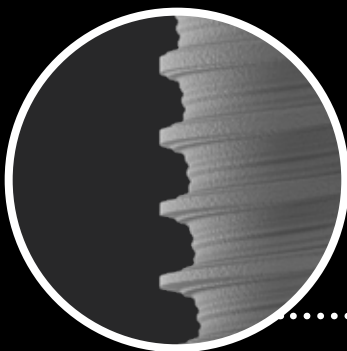


HELICAL APICAL NOTCH

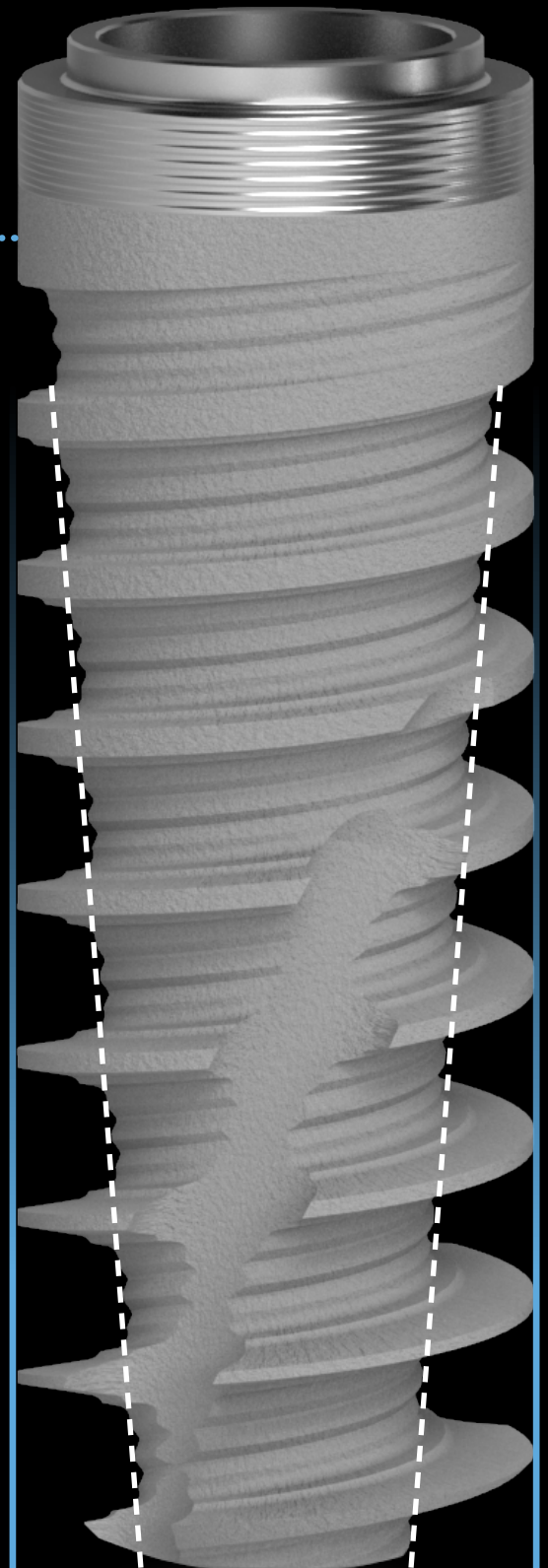

UTM and ZirTi: the synergy between two surfaces


UTM Surface: the substrate

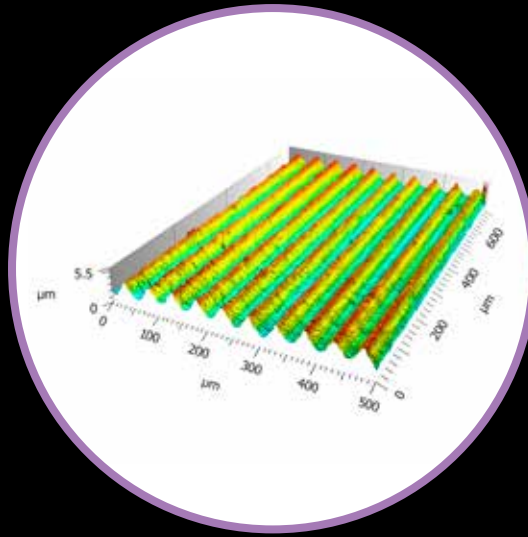

It's not a machined surface, but a regular micro-threading that provides the correct sizing to enhance cellular organization.

Sa: 0.6 μm

ZirTi surface: it's ideal for

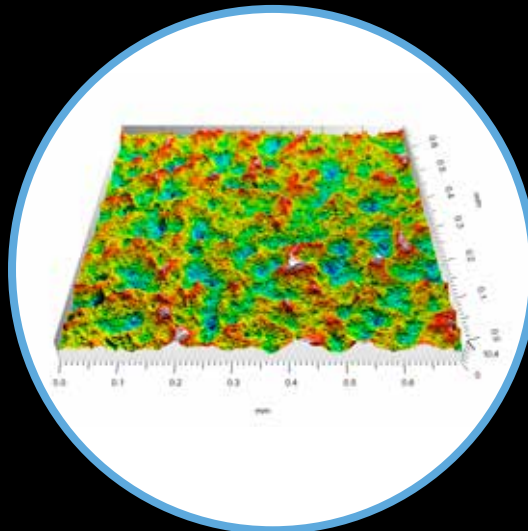


The ZirTi surface is obtained with a sandblasting process with zirconium oxide, followed by an acid etching, techniques that contribute to a topography to which the osteoblasts have shown to intimately adapt and to activate rapid proliferation.


Sa: 1.4-1.7 μm

Histology in humans at 3 months: even in bone D4, excellent mineralization is observed in contact with the ZirTi surface.

Image courtesy of Dr. Antonio Rocci and Prof. Marco Gargari.

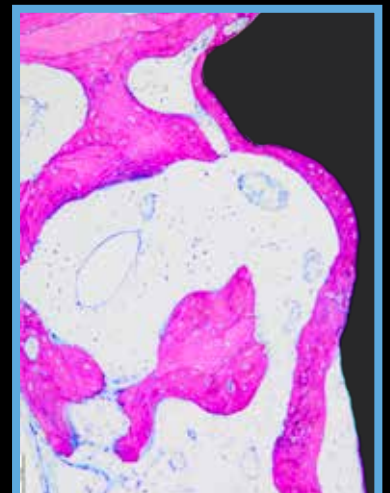


designed for soft tissues

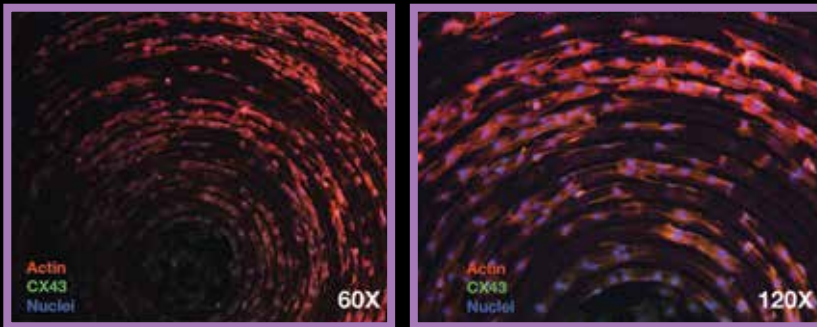

Arrangement of fibroblasts on UTM surface after 72 hours, observed with fluorescence microscopy - *in vitro* experiment.

The micro-thread on the implant neck provides a guide for the unidirectional movement of cells, with **the biological advantage of faster, organized, and oriented deposition of collagen fibers**. The resulting clinical benefit is a **quicker healing process** long-term maintenance of **healthy and stable tissues** over time.

osseointegration

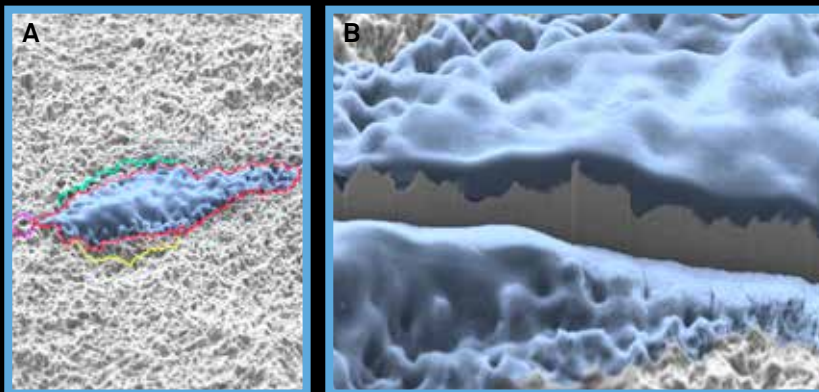


Enlargements of primary osteoblasts on the ZirTi surface: in images A and B, the **intimate adaptation of a cell on the peaks** that characterize the ZirTi can be observed, thanks to a section obtained with FIB (Focused Ion Beam).

Image courtesy of Prof. Macaluso, University of Parma, Department of Medicine and Surgery - Dentistry Center.



Bibliography on
ZirTi surface.

The advantages of the COLLEX UNIQUE CONNECTION

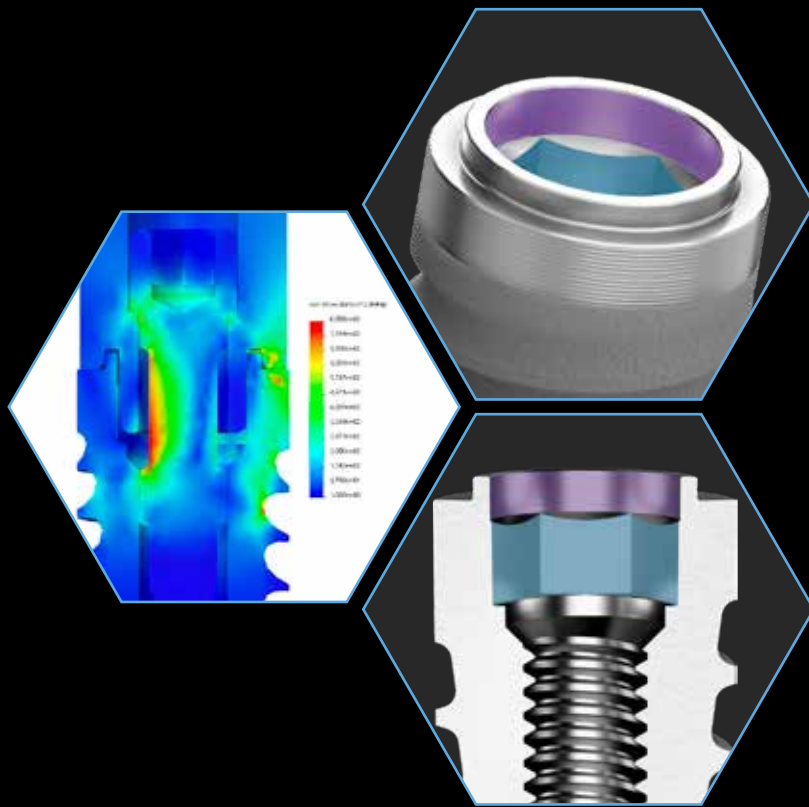

COLLEX ONE CONNECTION:

several studies have documented its robustness and prosthetic stability compared to other connections without a collar.





Bibliography
on Collex One
connection

INTERNAL HEXAGONAL


connection with a collar of the same size for all diameters.

**FOUR DIFFERENT
DIAMETERS** united by a
single connection.

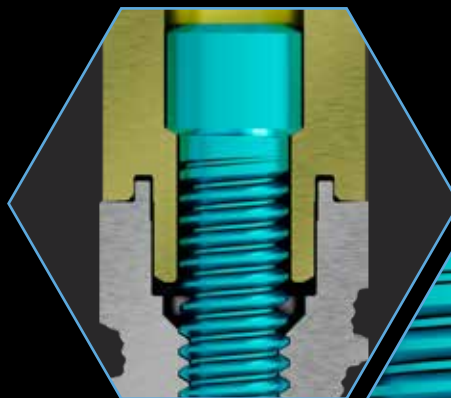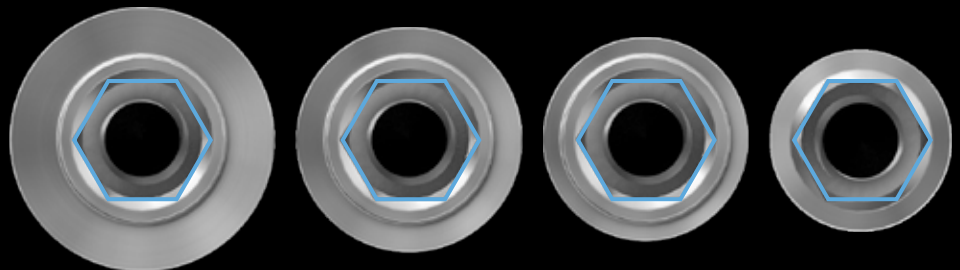

CONTRACONE SEAL:

- Special micromechanical machining that makes both surfaces conical for mutual fitting;
- Peripheral seal limiting bacterial infiltration;
- Documented in literature.

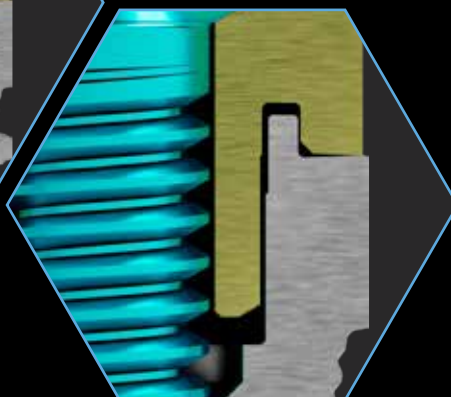



Shelta Power SL $\varnothing$ 5.00 mm
with $\varnothing$ 3.80 mm post



Shelta Power SL $\varnothing$ 5.00 mm
with $\varnothing$ 4.25 mm post




Shelta Power SL $\varnothing$ 5.00 mm
with $\varnothing$ 5.00 mm post

Possibility of prosthetic **PLATFORM SWITCHING** with a abutment of smaller diameter than the implant, or traditional prosthesis, with an abutment of diameter congruent with the connection platform.

The ability to mount components of different diameters on the same platform allows achieving **DIFFERENT MISMATCHING**, to obtain the biological and clinical benefits of platform switching without compromising mechanical strength.

$\varnothing$ implant	3.30 ●	3.80 ●	4.25 ●	5.00 ●
	3.30 ●	3.30 ●	3.30 ●	3.30 ●
$\varnothing$ compatible prosthetic component	-	3.80 ●	3.80 ●	3.80 ●
	-	-	4.25 ●	4.25 ●
	-	-	-	5.00 ●

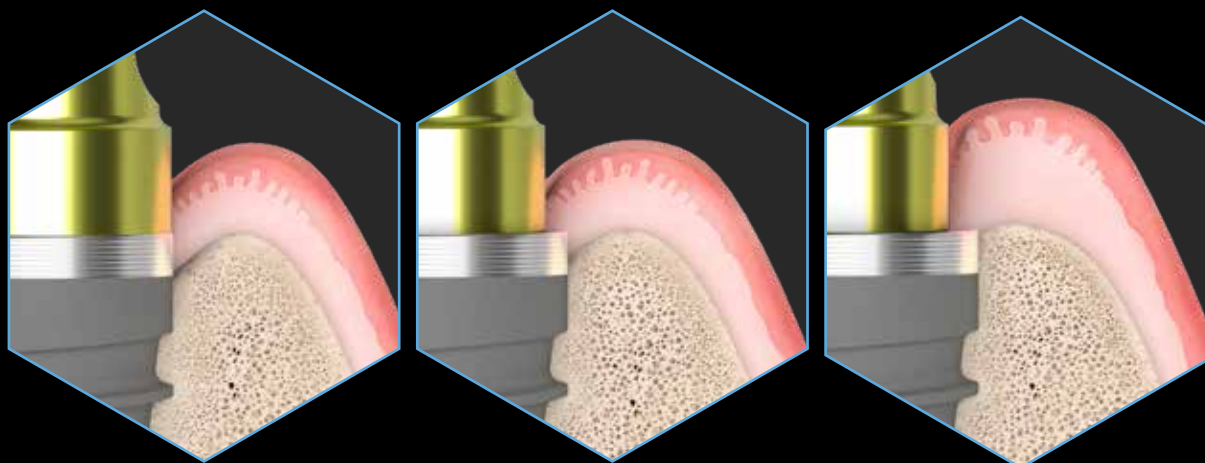

As demonstrated by the literature, there is a correlation between the extent of Platform Switching (mismatching) and the preservation of marginal bone dimensions: the greater the mismatching, the greater the volumes of hard and soft tissue preserved around the implant.





DIGITAL WORKFLOW

WITH PREMIUM POWER AND SHELTA POWER:


everything you need to save time, enhance treatment predictability, and boost the productivity of both the clinic and the laboratory



Guided surgery

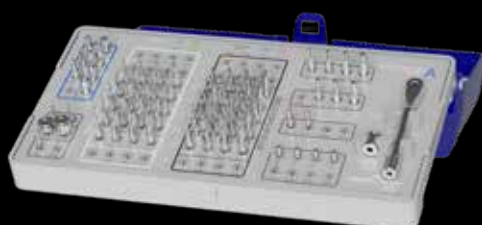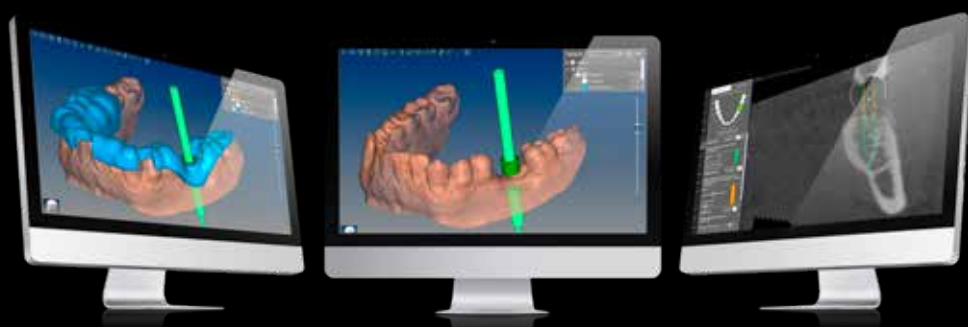

ECHOPLAN PRO A

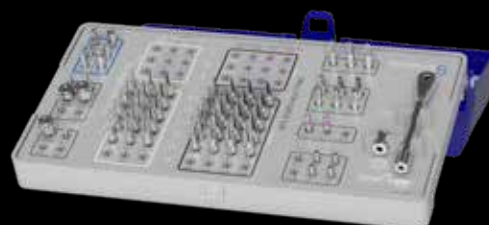

ECHOPLAN PRO S

Complete kits dedicated to guided surgery: the drills with integrated stops allow for time-saving and achieving predictable results. For the detailed content of the guided surgery kits compatible with Premium Power and Shelta Power, refer to page 25 and following.




DRILLS WITH INTEGRATED STOP

The uniqueness of the Echoplan PRO drills lies in their not requiring attachment to depth stops, as they come with an integrated vertical stop plane within the drills morphology. Additionally, the grooves on the guide body of the drill allow for increased physiological flow inside the preparation and reduce friction between the drill and the guide sleeve. This results in optimized surgical times as well as safe site preparation.



Data acquisition

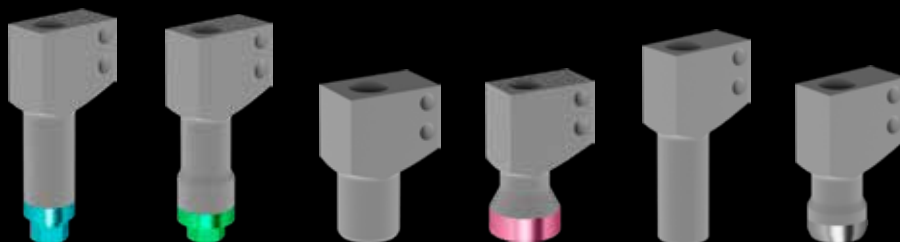


Capability to acquire the impression at the implant platform or abutment level, even for special prosthetic lines such as Plain B-Space, and XA.




New digital analogs: they allow for a unique and stable repositioning in 3D-printed models, thanks to their integration into CAD libraries.



CAD DESIGN

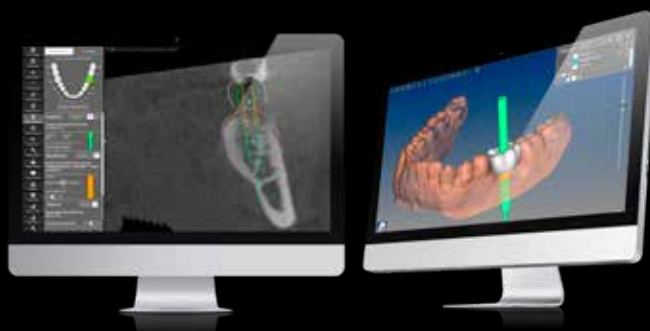

A wide range of T-Connect options for single and multiple rehabilitations, with cementation cones of different heights and customization, featuring straight and angled screw holes, allow finding the most suitable support for each milled structure. For the details of the T-Connects available for the Collex One connection, refer to page 53 and following.

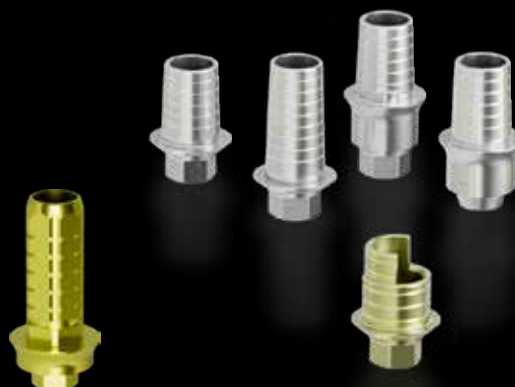


Ergonomics of Surgical Procedures

The **SURGICAL KITS** are comprehensive, **compact, ergonomic**, and contain all the necessary tools for the surgery of **Premium Power** and **Shelta Power** implants.

The color codes, also present on the surgical trays, guide the correct sequence of instrument usage.



Both the drills and the stops are characterized by a color code that guides operators in tool selection.

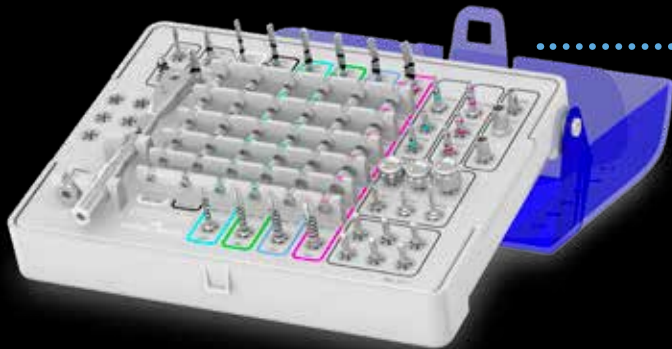

The **SURGICAL TRAYS** are made of **Radel**, a polymer that maintains its characteristics even after numerous cycles of washing and autoclave sterilization.

CYLINDRICAL or **CONICAL DRILLS**, based on different preparations, are supplied complete with depth stops.

The drills are made of **stainless steel** for surgical use and maintain their characteristics unchanged over time and through multiple uses.

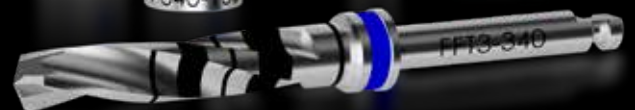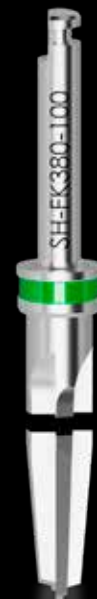

The drills are characterized by the presence of clearly visible depth notches.

The drills are **highly sharp, centered, and perform well**. The morphology of the conical drills allows for the collection of valuable bone for regeneration needs.





The kits include a convenient **RATCHET**, which can be used both as a fixed wrench and as a dynamometric wrench for controlling the implant insertion torque and tightening of prosthetic screws.

- The torque can be adjusted from 10 to 70 Ncm;
- The torque adjustment is highly accurate because each individual ratchet is individually calibrated.

Included is a wide range of **PROSTHETIC DRIVERS** for digital use, for the dynamometric ratchet and for micromotor.

The range is suitable for different oral access needs, from extra-short tools to longer instruments for Toronto-type prosthetic rehabilitations.

The slightly tapered coupling between the driver and the screw allows for adequate retention, ensuring secure transport of the screw in the oral cavity.



The surgical procedure involves a mountless approach that ensures visibility and ergonomics with the use of the patented **EASY INSERT** drivers:

- it can be easily removed from implant wells even in the presence of high insertion torque;
- the dodecagonal design prevents deformation of the connection edges, ensuring the maintenance of precision and prosthetic stability;
- available in the kit in various configurations, both for use with a micromotor and the ratchet;
- the instruments have a hexagon on the shaft, aligned with the connection hexagon, to facilitate rotational orientation of the fixture.

Premium Power kit

Complete surgical kit containing all the necessary tools for the insertion of Premium Power implants (code ZPREMIUM-ONE-INT).
All instruments can be purchased individually as spare parts.

1

Precision drill
FS-230



2

Intermediate drill*
◦ FG-200/280XS
• FG-330/425XS



3

Pilot drill*
• FPT3-200-LXS



4

Intermediate drill $\varnothing$ 2.50 mm*
◦ FFT3-250-LXS



20

Adapter for instruments
with digital knob
AVV3-MAN-DG



19

Torque control ratchet
CRI5-KIT



18

Countersink drills
• A-FCC-380
• A-FCC-425
• A-FCC-500



17

Stops for cylindrical
drills FFT3-425-LXS

- STOP4-425-070
- STOP4-425-085
- STOP4-425-100
- STOP4-425-115
- STOP4-425-130
- STOP4-425-150



16

Stops for cylindrical
drills FFT3-340-LXS

- STOP4-340-070
- STOP4-340-085
- STOP4-340-100
- STOP4-340-115
- STOP4-340-130
- STOP4-340-150



15

Stops for cylindrical drills
FFT3-300-LXS

- STOP4-300-070
- STOP4-300-085
- STOP4-300-100
- STOP4-300-115
- STOP4-300-130
- STOP4-300-150



14

Stops for cylindrical drills
FFT3-280-LXS

- STOP4-280-070
- STOP4-280-085
- STOP4-280-100
- STOP4-280-115
- STOP4-280-130
- STOP4-280-150



13

Stops for intermediate
drill $\varnothing$ 2.50 mm


- STOP4-250-070
- STOP4-250-085
- STOP4-250-100
- STOP4-250-115
- STOP4-250-130
- STOP4-250-150



12

Stops for pilot drill

- STOP4-200-070
- STOP4-200-085
- STOP4-200-100
- STOP4-200-115
- STOP4-200-130
- STOP4-200-150

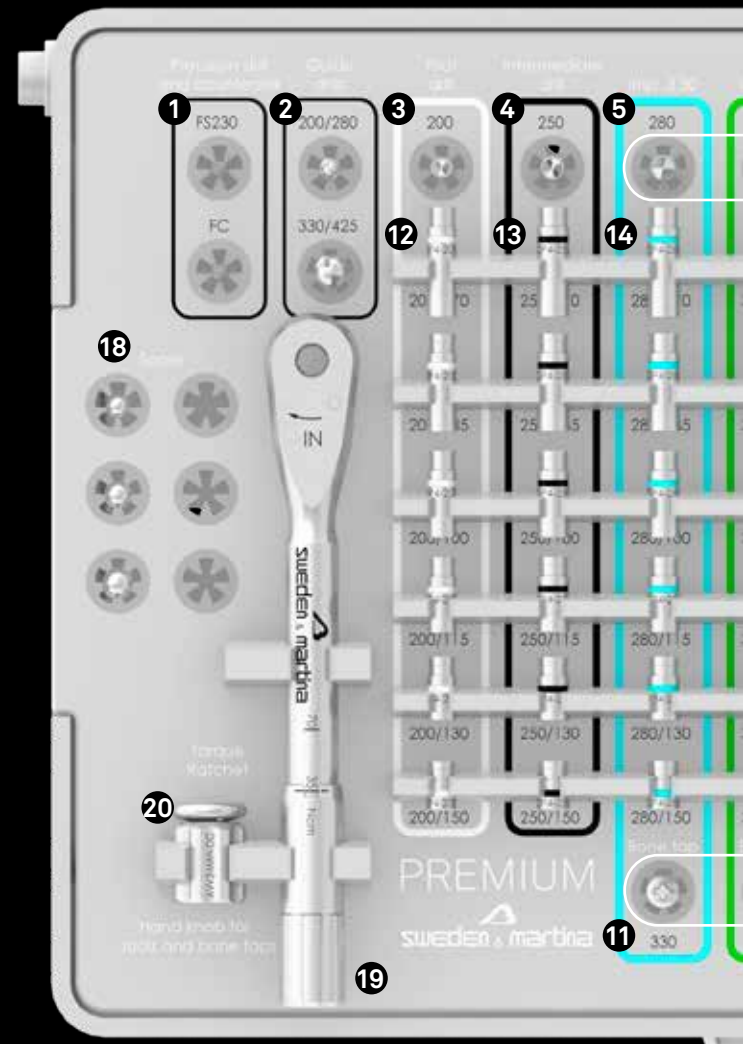

AS-TRAY-INT instrument carrying case in Radel.



18 *Refer to page 31 for warnings regarding instrument measurements

5



Final cylindrical drill*

- FFT3-280-LXS for implants $\varnothing$ 3.30 mm
- FFT3-300-LXS for implants $\varnothing$ 3.80 mm
- FFT3-340-LXS for implants $\varnothing$ 4.25 mm
- FFT3-425-LXS for implants $\varnothing$ 5.00 mm

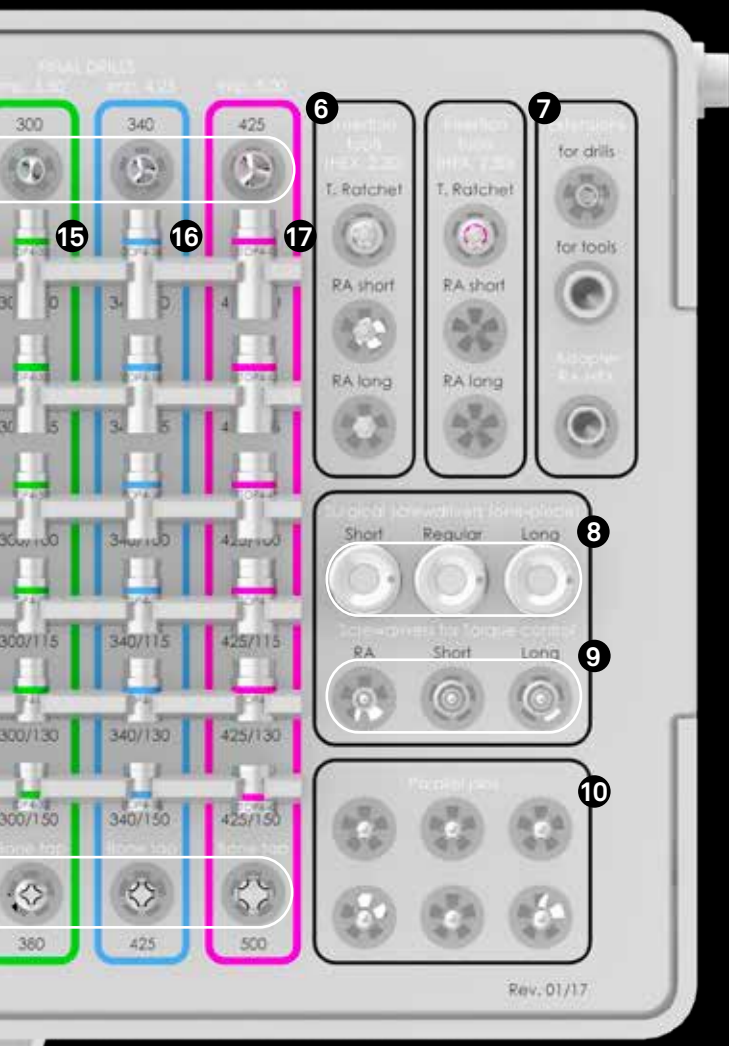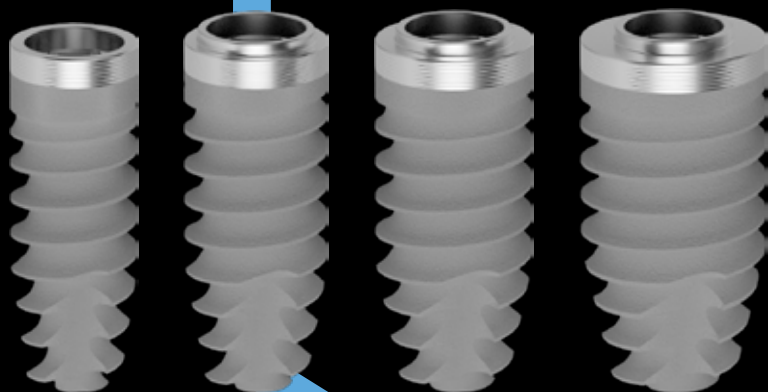

11



Bone taps*

- SM-MS-330 for implants $\varnothing$ 3.30 mm
- A-MS-380 for implants $\varnothing$ 3.80 mm
- A-MS-425 for implants $\varnothing$ 4.25 mm
- A-MS-500 for implants $\varnothing$ 5.00 mm

10



Parallelism pins PPTL-2-28

9

Prosthetic screwdrivers HSM-20-CA



8

Digital screwdrivers HSMXS-20-DG




6

Easy Insert Drivers for implant insertion EASYC4-EX230-CA



7

Extensions and adapters PROF-CAL3 Extension for surgical drills

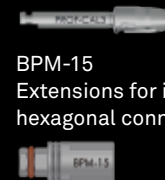

X-ray templates included in the kit, reorderable separately with the codes AS-L100 (real dimensions), AS-L120 (dimensions increased by 20%), AS-L130 (dimensions increased by 30%).

Shelta Power kit

Complete surgical kit containing all the necessary tools for the insertion of Shelta Power implants (code ZSHELTA-INT).
All instruments can be purchased individually as spare parts.

- 1** Precision drill
FS-230
- 2** Pilot drill
• CSR-FP-200
- 3** Intermediate drill $\varnothing$ 2.50 mm
• SE-FK250
- 4** Final cylindrical drills and stops
• FFT3-280-LXS
• CSR-STOP4-280-085
• CSR-STOP4-280-100
• CSR-STOP4-280-115
• CSR-STOP4-280-130
• CSR-STOP4-280-150

- 22** Bone taper
• SM-MS-330-CA for Premium Power implants $\varnothing$ 3.30 mm

- 20** Torque control ratchet
CRI5-KIT

- 21** Bone taps
• SH-MS-380-CA for implants $\varnothing$ 3.80 mm
• SH-MS-425-CA for implants $\varnothing$ 4.25 mm
• SH-MS-500-CA for implants $\varnothing$ 5.00 mm

- 18** Digital screwdrivers
HSMXS-20-DG

HSM-20-DG



HSMML-20-DG



- 19** Adapters and torque control ratchet
AVV-CA-DG-EX



- 17** Prosthetic screwdrivers
HSM-20-CA



HSM-20-EX



HSMML-20-EX



HSMXL-20-EX



- 16** Extensions
PROF-CAL3
Extension for surgical drills



BPM-15

Extension for instruments with hexagonal connector for torque-control ratchet



- 15** Parallelism pins
PPTL-2-28




- 14** Easy Insert Hex. 2.30 mm
EASYC4-EX230-CA



EASYL4-EX230-CA



EASY4-EX230-EX

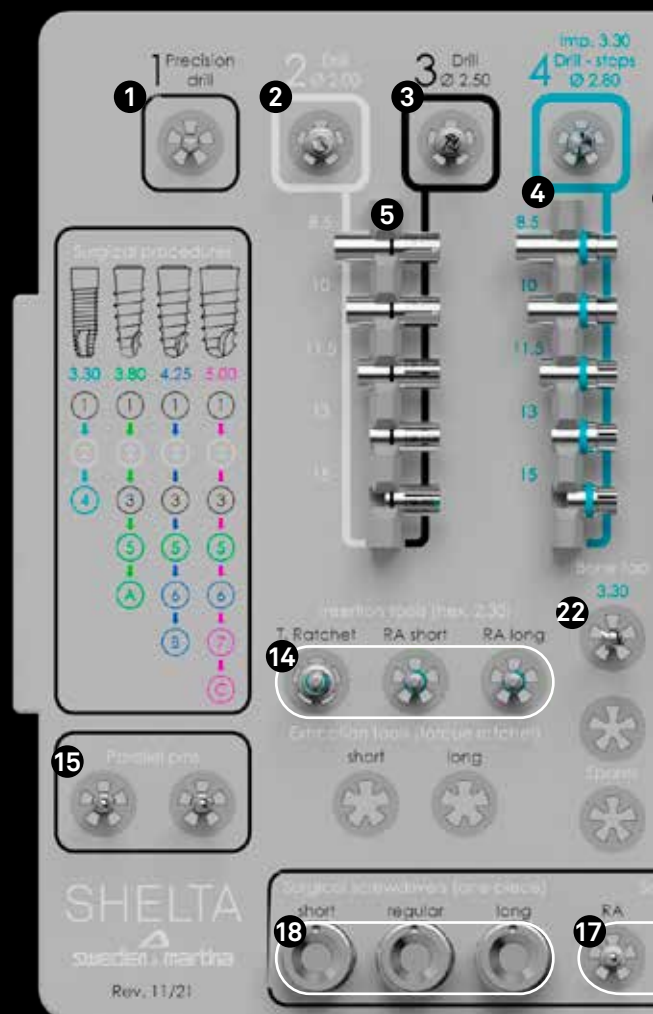

SH-TRAY-INT instrument carrying case in Radel.

*Refer to page 31 for warnings regarding instrument measurements.

5



Stops for pilot and intermediate drills



- SH-STOP4-FK380
- SH-STOP4-FK425
- SH-STOP4-FK500

6



Stop for final conical drills

- SH-FK380-085
- SH-FK380-100
- SH-FK380-115
- SH-FK380-130
- SH-FK380-150

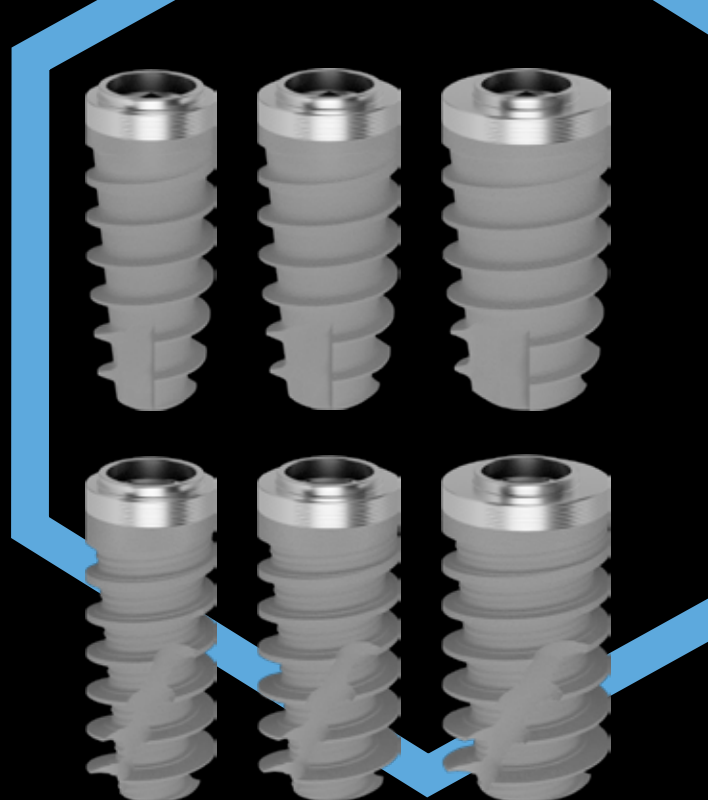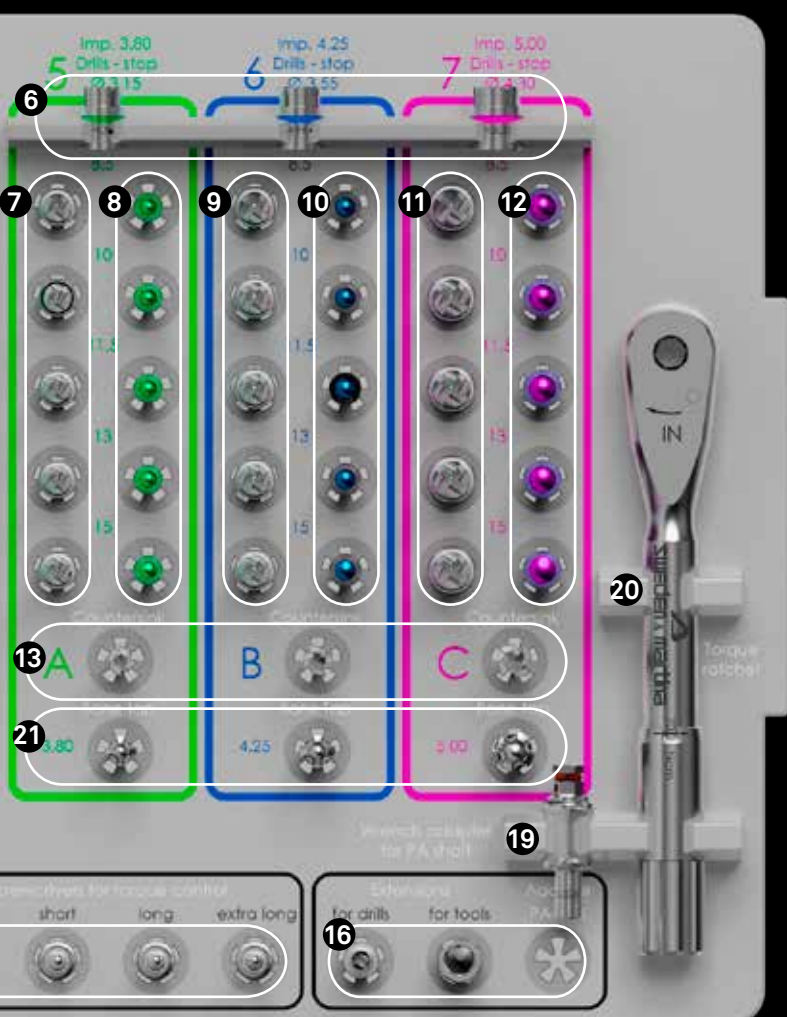

7



Final conical drills for ϕ 3.80 mm implants

- SH-FK380-085
- SH-FK380-100
- SH-FK380-115
- SH-FK380-130
- SH-FK380-150

8



Replies

- SH-380-085-RP
- SH-380-100-RP
- SH-380-115-RP
- SH-380-130-RP
- SH-380-150-RP

9



Final conical drills for ϕ 4.25 mm implants

- SH-FK425-085
- SH-FK425-100
- SH-FK425-115
- SH-FK425-130
- SH-FK425-150

13



Countersink drills

- A-FCC-380
- A-FCC-425
- A-FCC-500

12



Replies

- SH-500-085-RP
- SH-500-100-RP
- SH-500-115-RP
- SH-500-130-RP
- SH-500-150-RP

11



Final conical drills for ϕ 5.00 mm implants

- SH-FK500-085
- SH-FK500-100
- SH-FK500-115
- SH-FK500-130
- SH-FK500-150

10



Replies





- SH-425-085-RP
- SH-425-100-RP
- SH-425-115-RP
- SH-425-130-RP
- SH-425-150-RP

X-ray templates included in the kit, reorderable separately with the codes SH-L100 (real dimensions), SH-L120 (dimensions increased by 20%), SH-L130 (dimensions increased by 30%).

ATTENTION: For USA a different version of Shelta surgical kit is available: please contact your sales rep for details.

Combined Premium Shelta Power kit

Complete surgical kit containing all the necessary instruments for the insertion of both Premium Power and Shelta Power implants (code ZPRESH-INT).

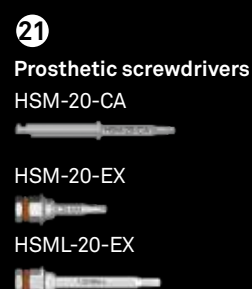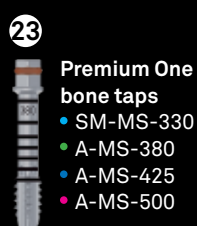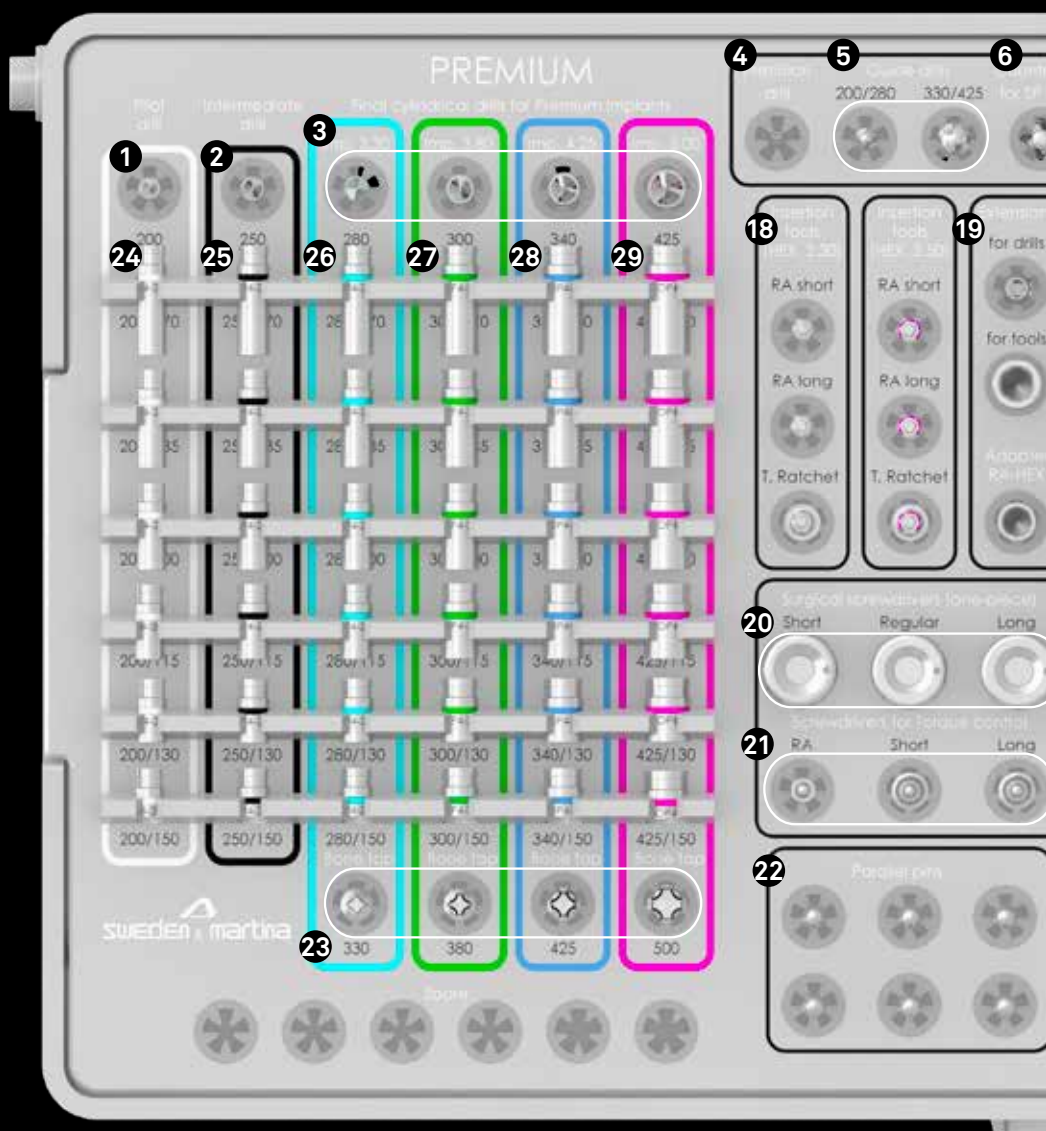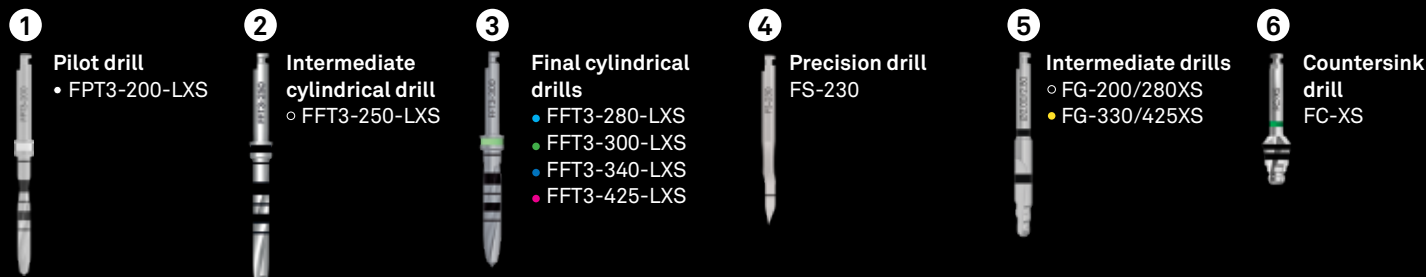

ASH-TRAY-INT Radel instrument grommetless tray.

22 *All the instruments in the Combi Kit are also available in the two kits: Premium Power on pages 18 and 19, and Shelta Power on pages 20 and 21.

7

Intermediate conical drill
 ○ SE-FK250



8

Stops for conical drills
 ● SH-STOP4-FK380
 ● SH-STOP4-FK425
 ● SH-STOP4-FK500



9

Stops for intermediate conical drill
 ● CSR-STOP-2028-085
 ● CSR-STOP-2028-100
 ● CSR-STOP-2028-115
 ● CSR-STOP-2028-130
 ● CSR-STOP-2028-150



10

Final conical drills
 ● SH-FK380-085
 ● SH-FK380-100
 ● SH-FK380-115
 ● SH-FK380-130
 ● SH-FK380-150



11

Replies
 ● SH-380-085-RP
 ● SH-380-100-RP
 ● SH-380-115-RP
 ● SH-380-130-RP
 ● SH-380-150-RP



12

Final conical drills
 ● SH-FK425-085
 ● SH-FK425-100
 ● SH-FK425-115
 ● SH-FK425-130
 ● SH-FK425-150



13

Replies
 ● SH-425-085-RP
 ● SH-425-100-RP
 ● SH-425-115-RP
 ● SH-425-130-RP
 ● SH-425-150-RP



14

Final conical drills
 ● SH-FK500-085
 ● SH-FK500-100
 ● SH-FK500-115
 ● SH-FK500-130
 ● SH-FK500-150




15

Replies
 ● SH-500-085-RP
 ● SH-500-100-RP
 ● SH-500-115-RP
 ● SH-500-130-RP
 ● SH-500-150-RP



16

Shelta Bone taps
 ● SH-MS-380-CA
 ● SH-MS-425-CA
 ● SH-MS-500-CA

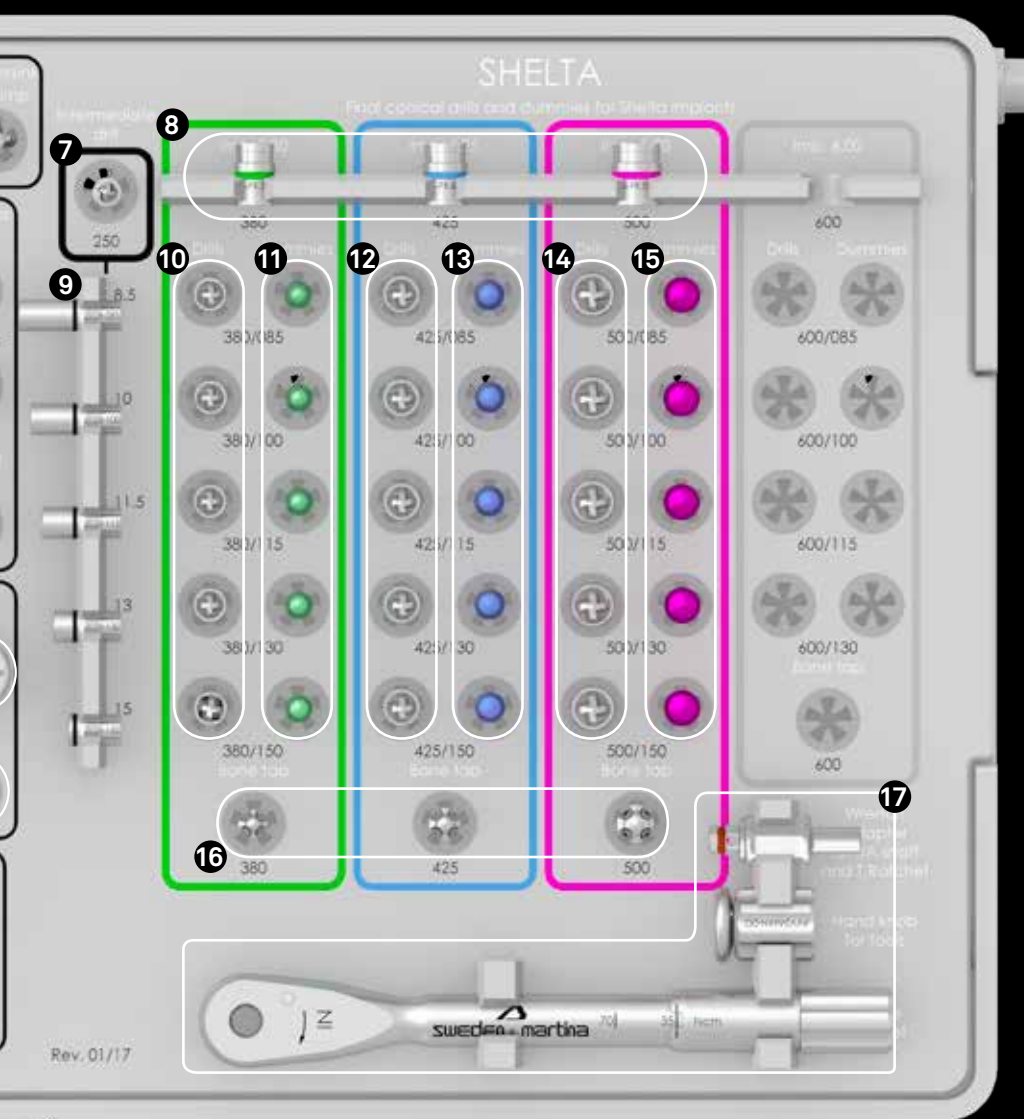

20

Digital screwdriver
 HSMXS-20-DG



HSM-20-DG



HSML-20-DG



19

Extensions and adapter
 PROF-CAL3



BPM-15



B-AVV-CA3



18

Easy Insert Hex. 2.30 mm
 EASYC4-EX230-CA



EASYL4-EX230-CA




EASY4-EX230-EX



17

Adapters and torque control ratchet
 AVV-CA-DG-EX

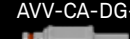

AVV3-MAN-DG




CRI5-KIT










X-ray templates included in the kit

Surgical kit for guided surgery with Premium One Power implants

Complete surgical kit containing all the necessary instruments for guided surgery with Premium Power implants (code ZGS-PRO-A-INT). In the Tray there are a number of instrument slots referring to different systematics. These instruments are therefore not included in the standard kit.

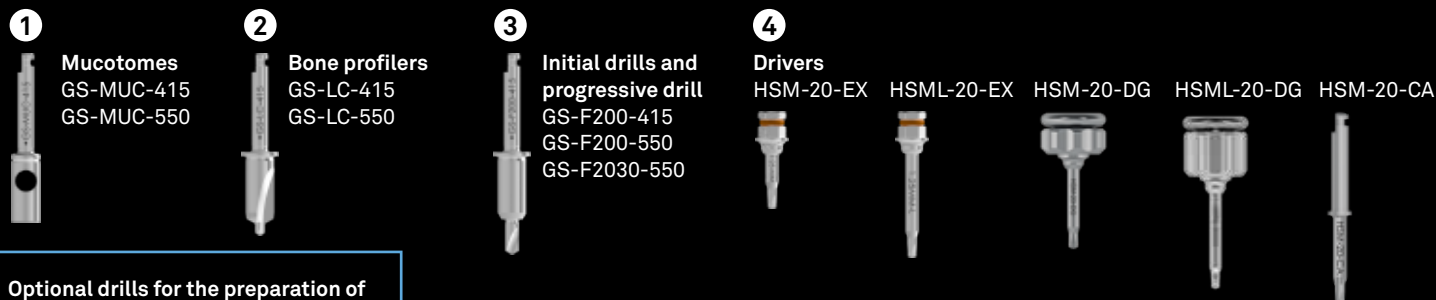

Optional drills for the preparation of implants h 18.00 mm.

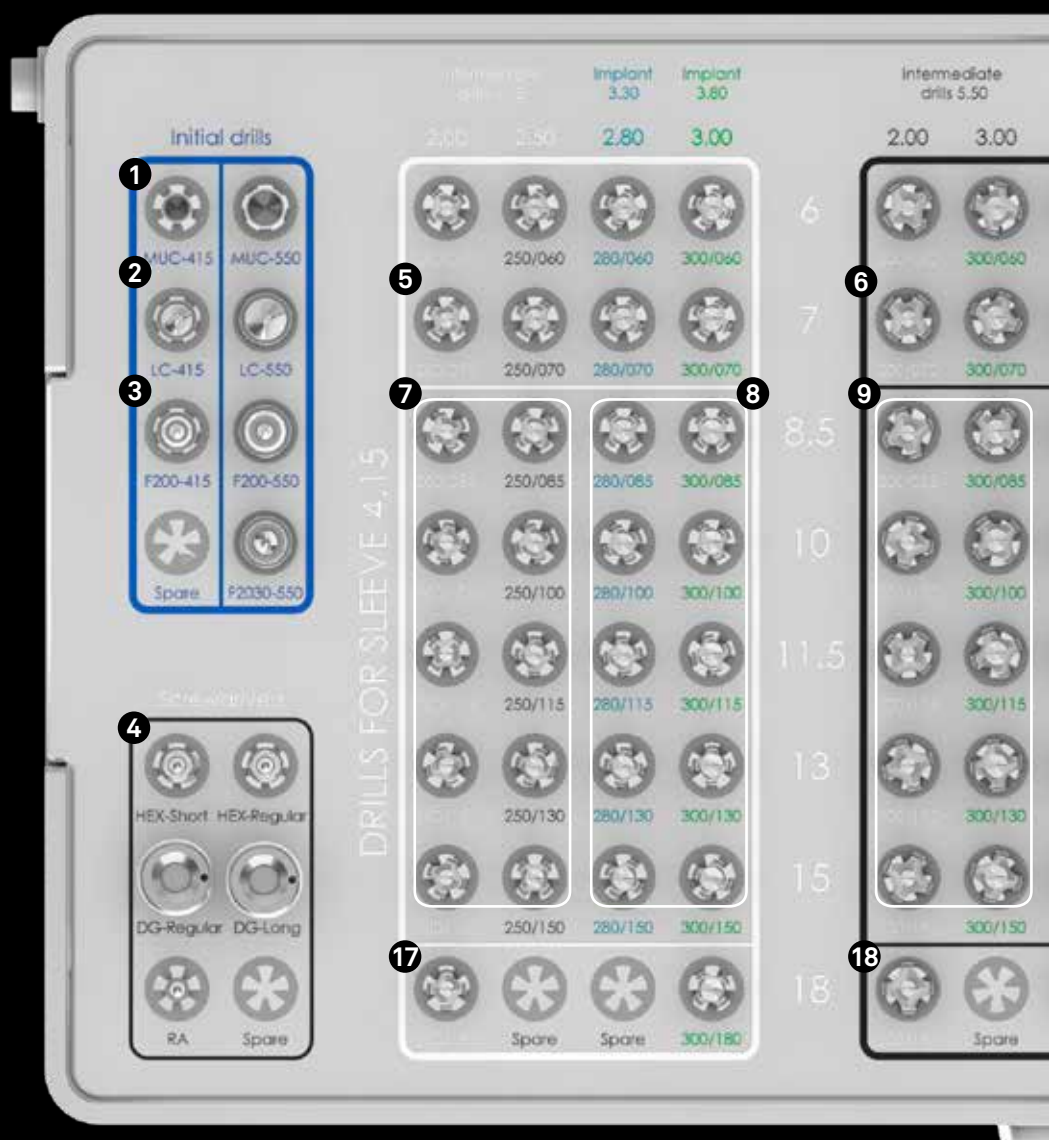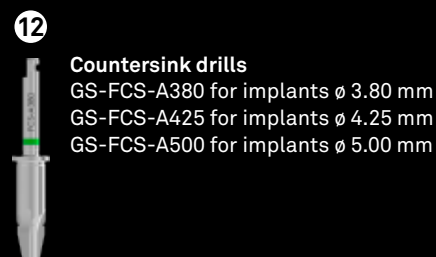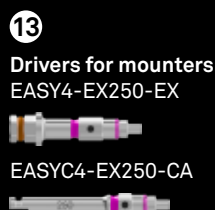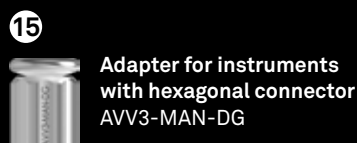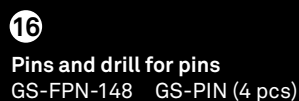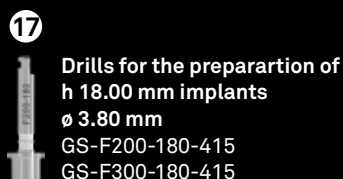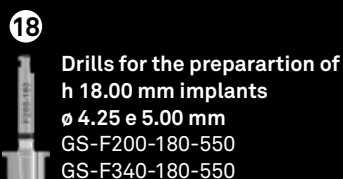

Optional drills for preparing 7.00 mm high implants can be ordered also cumulatively as a supplementary set with the code GS-PROA-INTEGRA-070.

5



Optional drills for h 7.00 mm
 $\varnothing$ 3.80 mm implants
 GS-F200-070-415
 GS-F250-070-415
 GS-F280-070-415
 GS-F300-070-415

6




Optional drills for h 7.00 mm
 $\varnothing$ 4.25 and 5.00 mm implants
 GS-F200-070-550
 GS-F300-070-550
 GS-F340-070-550
 GS-F425-070-550

7



Intermediate drills
 for $\varnothing$ 4.15 mm sleeve
 GS-F200-085-415
 GS-F200-100-415
 GS-F200-115-415
 GS-F200-130-415
 GS-F200-150-415
 GS-F250-085-415
 GS-F250-100-415
 GS-F250-115-415
 GS-F250-130-415
 GS-F250-150-415

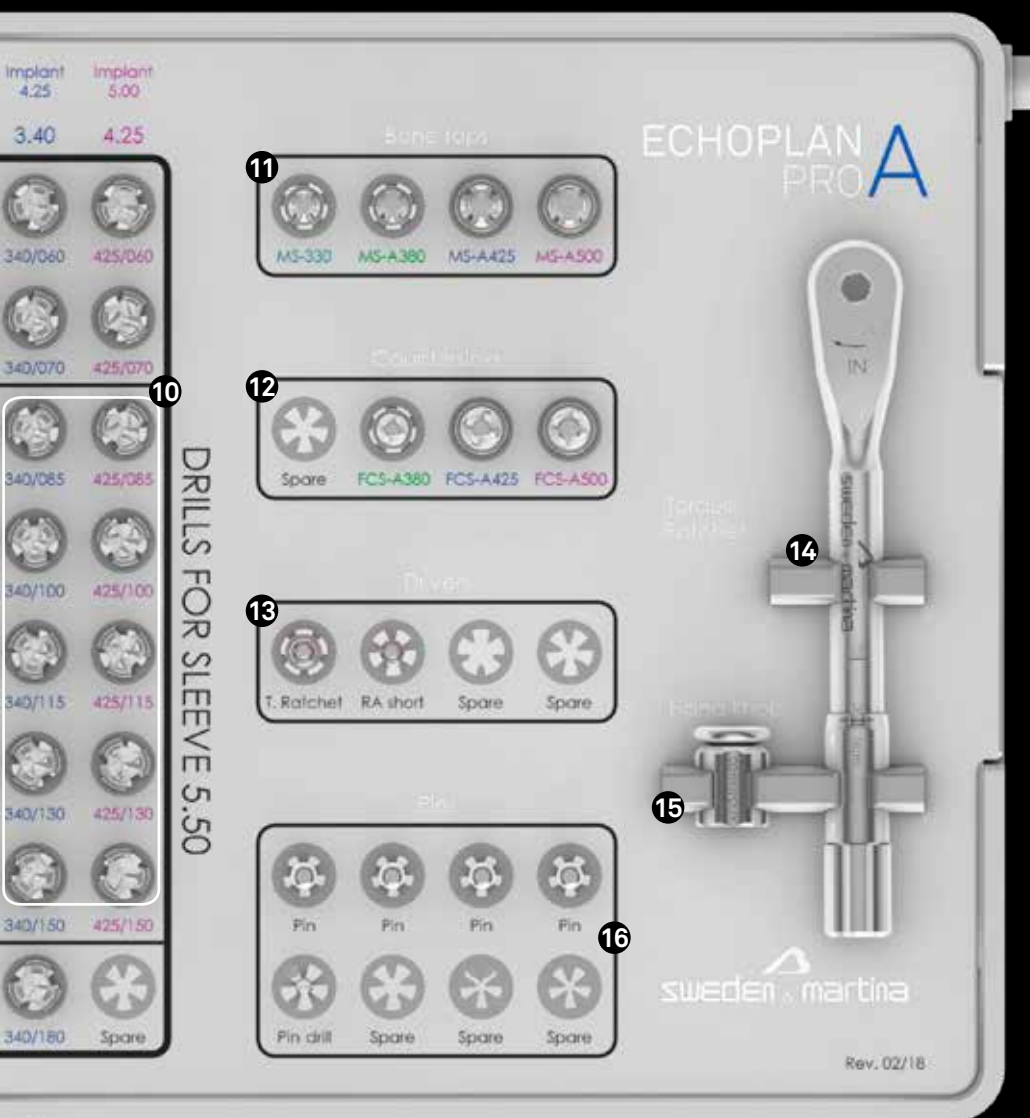

8



Final drills for
 $\varnothing$ 4.15 mm sleeve
 GS-F280-085-415
 GS-F280-100-415
 GS-F280-115-415
 GS-F280-130-415
 GS-F280-150-415
 GS-F300-085-415
 GS-F300-100-415
 GS-F300-115-415
 GS-F300-130-415
 GS-F300-150-415

9



Intermediate drills for
 $\varnothing$ 5.50 mm sleeve
 GS-F200-085-550
 GS-F200-100-550
 GS-F200-115-550
 GS-F200-130-550
 GS-F200-150-550
 GS-F300-085-550
 GS-F300-100-550
 GS-F300-115-550
 GS-F300-130-550
 GS-F300-150-550

11



Bone taps
 GS-SM-MS-330 for implants $\varnothing$ 3.30 mm
 GS-MS-A380 for implants $\varnothing$ 3.80 mm
 GS-MS-A425 for implants $\varnothing$ 4.25 mm
 GS-MS-A500 for implants $\varnothing$ 5.00 mm

10



Final drills for $\varnothing$ 4.25 mm implants
 with $\varnothing$ 5.50 mm sleeve
 GS-F340-085-550
 GS-F340-100-550
 GS-F340-115-550
 GS-F340-130-550
 GS-F340-150-550

Final drills for $\varnothing$ 5.00 mm implants
 with $\varnothing$ 5.50 mm sleeve
 GS-F425-085-550
 GS-F425-100-550
 GS-F425-115-550
 GS-F425-130-550
 GS-F425-150-550

Surgical kit for guided surgery with Shelta Power and Shelta SL Power implants

Complete surgical kit containing all the necessary instruments for guided surgery with Shelta Power implants (code ZGS-PRO-S-INT). In the Tray there are a number of instrument slots referring to different systematics. These instruments are therefore not included in the standard kit.

1



Mucotomes
GS-MUC-415
GS-MUC-550

2



Bone profilers
GS-LC-415
GS-LC-550

3



Initial drills and progressive drill
GS-F200-415
GS-F200-550

16




Adapter for instruments with hexagonal connector
AVV3-MAN-DG

15



Torque control ratchet
CRI5-KIT

14

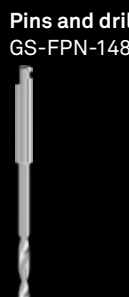

Pins and drill for pins
GS-FPN-148 GS-PIN (4 pcs)

13



Drivers for mounters
EASY4-EX250-EX EASYC4-EX250-CA

12




Countersink drills
GS-FCS-SH380
GS-FCS-SH425
GS-FCS-SH500
GS-FCS-SE380
GS-FCS-SE425
GS-FCS-SE500

11



Bone taps
GS-MS-S380
GS-MS-S425
GS-MS-S500

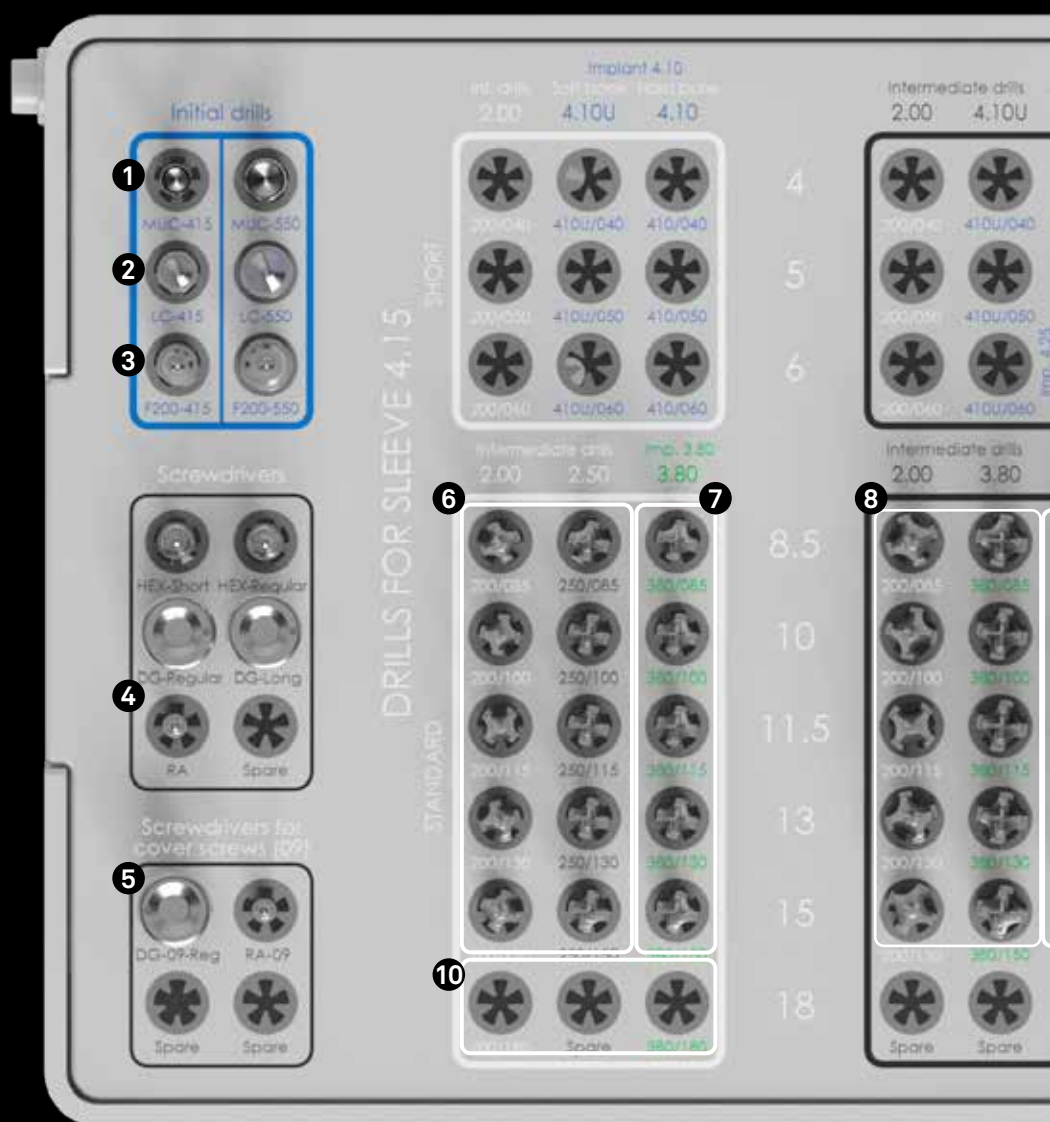

Guided surgery accessories

Guide sleeves and mounter for guided surgery

implants	ø	sleeve	mounter	screw	easy mounter
Premium Power Slim	3.30	 GS-B415-EX-6* GS-B415-EX-SET*	 GS-MOU-A330	 GS-VTMOU-180**	 GS-EASY-AS-415-CA  GS-EASY-AS-415-EX
Premium Power and Shelta Power	3.80	 GS-B415-EX-6* GS-B415-EX-SET*	 GS-MOU-A380	 GS-VTMOU-180**	 GS-EASY-AS-415-CA  GS-EASY-AS-415-EX
Premium Power and Shelta Power	4.25 5.00	 GS-B550-EX-6* GS-B550-EX-SET*	 GS-MOU-A380SP	 GS-VTMOU-180**	 GS-EASY-AS-550-CA  GS-EASY-AS-550-EX

*The indexed guide sleeves are available in packs of 6 pieces, and in sets of 100 pieces, not for sale individually.

**The fixation screw is included with tmounter but can also be ordered individually as a spare part.

Optional surgical instruments

Intermediate drills



FG-200/280XS

Long precision drill



FSL-230

Cylinder for the construction of radiological guide



DIM

*in titanium for guiding the precision drill

Drills for the preparation of Premium Power implants with a height of 7.00 mm*

pilot drill



FPS-200

intermediate Shorty drill



FGS-200/300



FFS-300



FFS-340



STOPPS4-200-070

-



STOPPS4-300-070



STOPPS4-340-070

*All preparation drills are included in the ZSHORTY-INT drill kit, which also contains drills for the preparation of shorter implants (h 5.00 mm and h 6.00 mm) from other systems.

Drills for distal sectors

ø 2.00



FPT5-200-LXS

ø 2.80



FFT5-280-LXS

ø 2.90



FFT5-290-LXS

ø 3.00



FFT5-300-LXS

ø 3.20



FFT5-320-LXS

ø 3.30



FFT5-330-LXS

ø 3.40



FFT5-340-LXS

ø 3.60



FFT5-360-LXS

ø 4.25



FFT5-425-LXS

ø 4.45



FFT5-445-LXS

Screwdrivers not included in the kit



HSMXL-20-EX

with hexagonal connector, extra long



HSMXS-20-CA

with right angle shank, extra short



AVV-CAP-TIT-1

instrument for maintaining the titanium cap for ball attachments



BASCC-EX

for ball attachments, with hexagonal connector



AVV2-ABUT

for straight P.A.D. abutments



8926-SW

for Locator abutment, short



8927-SW

for Locator abutment, long



DSPDCLH-24

for angled screw hole, length 24 mm




DSPDCLH-32

for angled screw hole, length 32 mm

Parallel pins


ø 2.80 mm

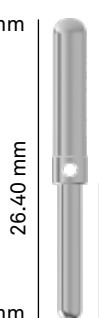

PP-2/28

Parallelism pin with one side of ø 2.00 mm and the other side of ø 2.80 mm


ø 2.80 mm

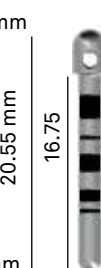

PPTS-2-28

Parallelism pin with depth notches, small version

ø 2.80 mm

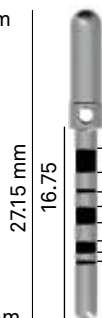

PPTL-2-28

Parallelism pin with depth notches, large version

ø 2.80 mm



PPS-2/3

Parallelism pin for short drills with notches at 5.00, 6.00, and 7.00 mm

Bone taps with hexagonal connector

	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
	SM-MSC-330	A-MSC-380	A-MSC-425	A-MSC-500
	-	SH-MS-380	SH-MS-425	SH-MS-500

Bone taps with right angle shank

	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
	SM-MS-330-CA	A-MS-380-CA	A-MS-425-CA	A-MS-500-CA

Extraction tools for intraoral removing of implants

	BC-EX230	short
	BL-EX230	long

Bone profiler for P.A.D. abutments

	transgingival height 3.00 mm, wide			
	ø 3.30 	ø 3.80 	ø 4.25 	ø 5.00 
	A-PAD-PS330-L	A-PAD-PS380-L	A-PAD-PS425-L	A-PAD-PS500-L
	transgingival height 5.00 mm, narrow			
	ø 3.30 	ø 3.80 	ø 4.25 	ø 5.00 
	A-PAD-PS330-S	A-PAD-PS380-S	A-PAD-PS425-S	A-PAD-PS500-S
	A-PAD-GUI-PS-230	profiler guide		

Conical osteotomes

				
__ ø 0.20	__ ø 0.90	__ ø 1.60	__ ø 2.00	__ ø 2.40
E-OS-020-PP	E-OS-090-PP	E-OS-160-PC	E-OS-200-PC	E-OS-240-PC
osteotome ø 0.20 flat tip	osteotome ø 0.90 flat tip	osteotome ø 1.60 concave tip	osteotome ø 2.00 concave tip	osteotome ø 2.40 concave tip

Conical osteotomes with flat tip

	h	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
	8.50-10.00	SH-OS-380-100-PP	SH-OS-425-100-PP	SH-OS-500-100-PP
	11.50	SH-OS-380-115-PP	SH-OS-425-115-PP	SH-OS-500-115-PP
	13.00	SH-OS-380-130-PP	SH-OS-425-130-PP	SH-OS-500-130-PP
	15.00	SH-OS-380-150-PP	SH-OS-425-150-PP	SH-OS-500-150-PP

Conical osteotomes with round tip

	h	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
	8.50-10.00	SH-OS-380-100-PR	SH-OS-425-100-PR	SH-OS-500-100-PR
	11.50	SH-OS-380-115-PR	SH-OS-425-115-PR	SH-OS-500-115-PR
	13.00	SH-OS-380-130-PR	SH-OS-425-130-PR	SH-OS-500-130-PR
	15.00	SH-OS-380-150-PR	SH-OS-425-150-PR	SH-OS-500-150-PR

Torque wrench with control lever



TWL



TWL-AVV-EX

Adaptor for torque wrench

Service mounter for impact procedures



MOU-EX230

mounter



CM2

mounter stop key

Depth gauge



PROF3

Spare O-ring



ORING180-088

Extractors for atraumatic removal of implants



EXTR-IMP-S

small



EXTR-IMP-M

medium



EXTR-IMP-L

large


**The extractor is not specific to the system connection; it must be selected according to the implant to be removed.*

Depth notches




The laser marking on the cutting edges and the use of stops allow for safe preparation and control of the instrument's insertion depth. The drills always create a hole longer than the desired implant. The oversizing (LS) is equal to the height of the tip of the drill being used.

LL: total length of the working part, including the tip.

LS: length of the overpreparation. This measurement must be calculated in addition to the length of the preparation hole.

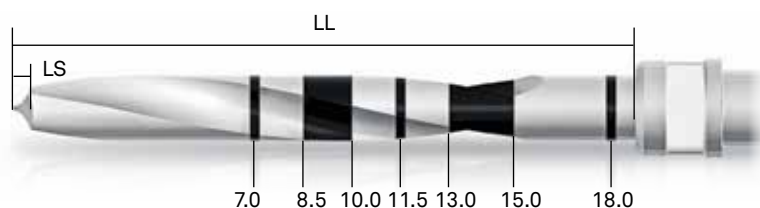

The working depth is also identified on bone taps and osteotomes through laser-etched notches, as follows:

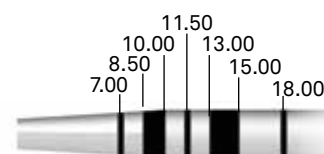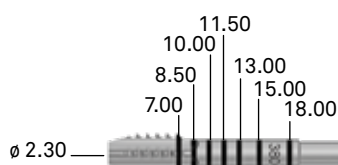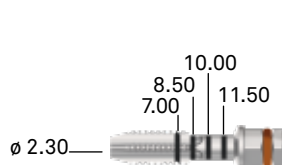

Prosthetic solutions of excellence

CONVERGENT PLAIN B-SPACE ABUTMENT FOR MULTIPLE SCREW RETAINED SOLUTIONS

A revolutionary study using synchrotron tomography has demonstrated that the synergy between converging morphology and surface, as opposed to divergent morphologies and non-UTM surfaces, enhances fibroblast adhesion and early activation, facilitating the formation of a biological seal and promoting the integration of peri-implant soft tissues.

Connective fibers appear more interconnected, dense, and stable.

UTM surface along their entire profile.

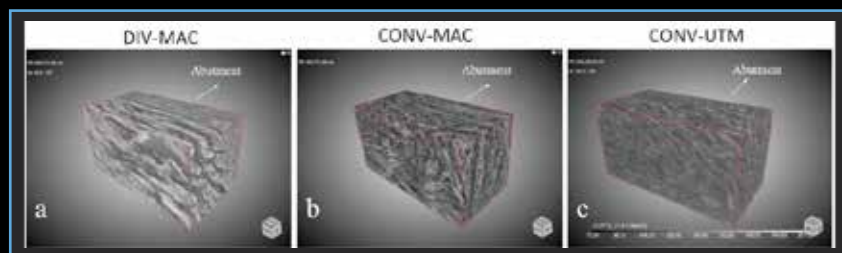


Image courtesy of Dr. L. Canullo, A. Giuliani, M. Furlani, M. Menini, A. Piattelli, G. Iezzi.



CONICO SYSTEM



Fixed prosthesis without screws and without cement

Easily removable and reinsertable by the clinician

Biological seal and maintenance of healthy and stable tissues over time

For both single and multiple crowns

Not removable by the patient

XA PROSTHESIS

Convergent morphology that provides more space for soft tissues and supports bony peaks.

They are designed for the One-Abutment-One-Time technique and platform switching solutions.

At the base of the posts for cemented solutions and abutments for screw-retained solutions, there are micro-rills that allow for the stabilization of the connective tissue.

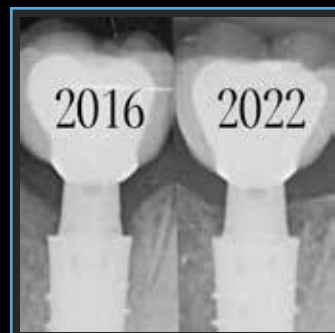


Image courtesy of Dr. Xavier Vela

SOLUTIONS FOR ANGLED SCREW HOLE

They are versatile prosthetic solutions that reposition the screw hole palatally to preserve aesthetics.

ORIENTABLE SCREW HOLE ABUTMENT



The castable cannula rotates on the spherical head of the abutment, allowing 360° of angulation. The cannula rotates on the spherical head of the abutment, enabling a 28° angle correction in all directions.



ANGLED SCREW HOLE INTERFASE



The hexagonal screwdriver with a special design allows for engaging the screw head even in the presence of a pronounced angle and a palatally repositioned screw hole.

Angled screw hole Interfases are an ideal solution for digital workflows.

PAD^e AESTHETIC SCREW-RETAINED DISPARALLEL PROSTHESIS



Easy-to-insert abutments

Maximum freedom and ergonomics for superstructures

Pink anodization for better mimetism

Their compact size makes them suitable for rehabilitation of partially or completely edentulous arches in limited horizontal and vertical spaces.

Specific sleeves are available for luting technique with a standard hole or an angled hole for the palatal repositioning of screw access.

Implant Platforms and Prosthetic Combinations

The connection platforms of Premium Power and Shelta Power implants are characterized by the presence of an internal hexagon and an external collar that is the same for all diameters, allowing for the selection of the desired abutment size. According to the chosen prosthetic protocol, it is possible to perform prosthetic Platform Switching with a smaller diameter abutment or use an abutment with a congruent diameter to the implant.

ø implant	3.30 ●	3.80 ●	4.25 ●	5.00 ●
ø compatible prosthetic component	3.30 ●	3.30* ●	3.30* ●	3.30* ●
	-	3.80 ●	3.80 ●	3.80 ●
	-	-	4.25 ●	4.25 ●
	-	-	-	5.00 ●



Premium Power ø 5.00 mm
with post ø 3.80 mm




Premium Power ø 5.00 mm
with post ø 4.25 mm



Premium Power ø 5.00 mm
with post ø 5.00 mm


Please note

*Prosthetic components of ø 3.30 mm result in Platform Switching with implants of ø 3.80, 4.25, and 5.00 mm, closing onto the collar instead of the platform (see side image). It is recommended to use these abutments for single crowns exclusively in the anterior regions (excluding premolars) and in the distal areas for supporting multiple prostheses.











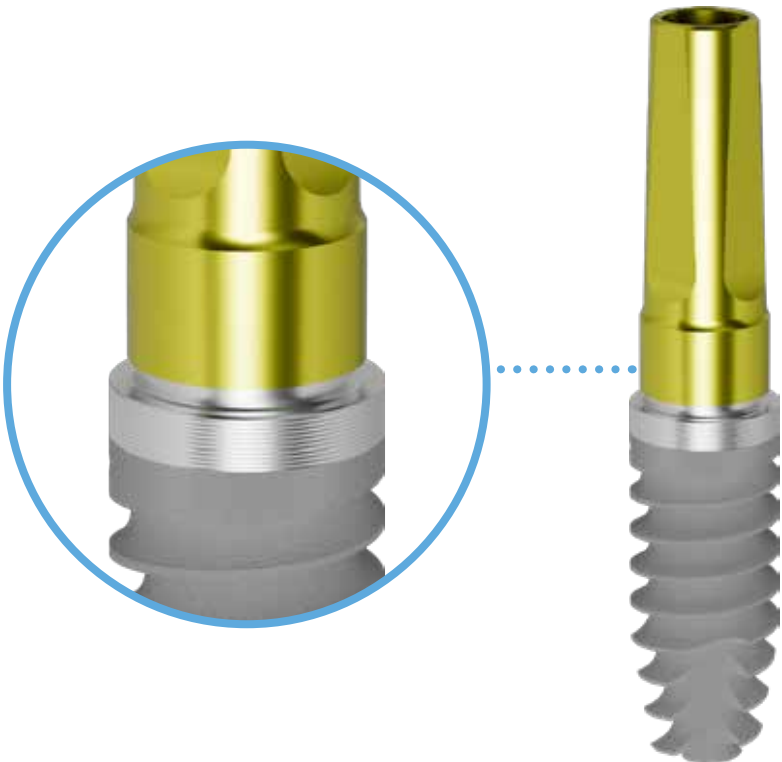


Healing abutments

With straight emergence

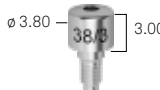h mm	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
2.00	 A-TMG-330-2	 A-TMG-380-2		
3.00	 A-TMG-330-3	 A-TMG-380-3		
5.00	 A-TMG-330-5	 A-TMG-380-5		
7.00	 A-TMG-330-7	-		

With anatomic emergence

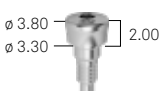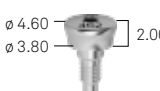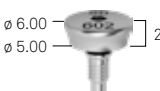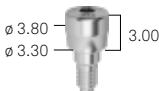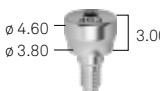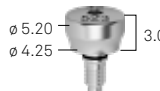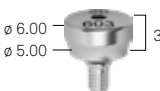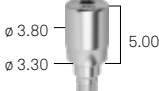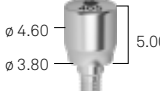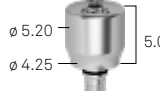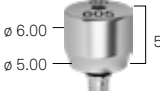h mm	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
2.00	 A-TMGR-330-2	 A-TMGR-380-2	 AS-TMGR-425-2	 AS-TMGR-500-2
3.00	 A-TMGR-330-3	 A-TMGR-380-3	 AS-TMGR-425-3	 AS-TMGR-500-3
5.00	 A-TMGR-330-5	 A-TMGR-380-5	 AS-TMGR-425-5	 AS-TMGR-500-5
7.00	-	 A-TMGR-380-7	 AS-TMGR-425-7	 AS-TMGR-500-7

With coronal enlargement

h mm $\varnothing$ 3.30 ● $\varnothing$ 3.80 ●




A-TMTCR-330-2

A-TMTCR-380-2



A-TMTCR-330-3

A-TMTCR-380-3

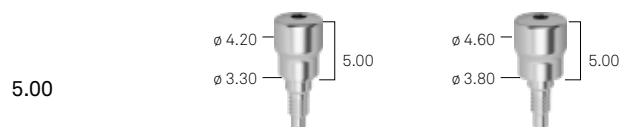

A-TMTCR-330-5

A-TMTCR-380-5

With wide profile

h mm $\varnothing$ 3.30 ● $\varnothing$ 3.80 ●




A-TMGB-330-2

A-TMGB-380-2



A-TMGB-330-3

A-TMGB-380-3

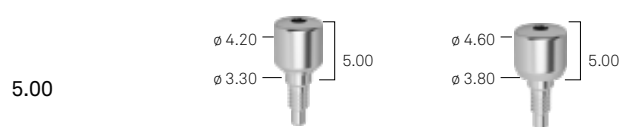

A-TMGB-330-4

A-TMGB-380-4

REGISTRA technique


Transfer screws

$\varnothing$ 3.30 ●


$\varnothing$ 3.80 ●

$\varnothing$ 4.25 ●


Anatomical transfer screw for REGISTRA technique, height 5.00 mm.

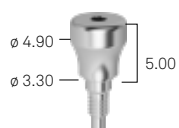

A-TO-330-5

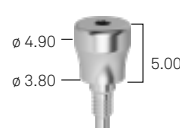


A-TO-380-5

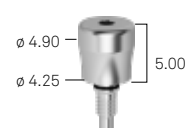

AS-TO-425-5

Anatomical transfer screw for REGISTRA technique, height 7.00 mm.

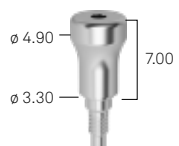

A-TO-330-7



A-TO-380-7




AS-TO-425-7


Anatomical transfer screw for REGISTRA technique, height 9.00 mm.



A-TO-330-9

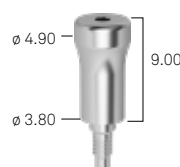

A-TO-380-9

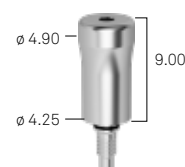

AS-TO-425-9

Transfer cap









Universal cap for REGISTRA transfer screws



CAP-TO





Impression and model phase

Pick-up transfers

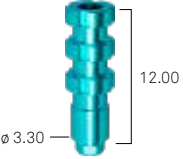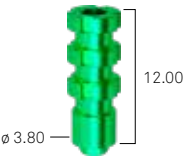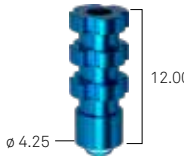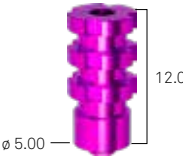
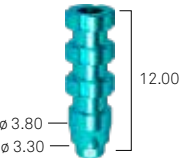
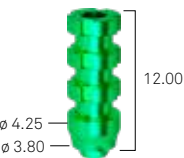
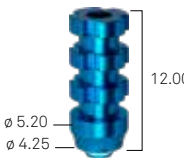
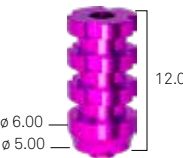
	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
straight emergence	 A-TRA-330	 A-TRA-380	 AS-TRA-425	 AS-TRA-500
anatomical emergence	 A-TRAR-330	 A-TRAR-380	 AS-TRAR-425	 AS-TRAR-500

VTRA2-180-15 screw included.

Pull-up transfers

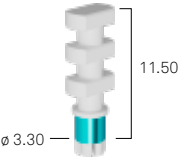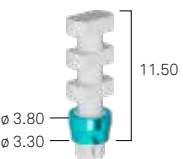
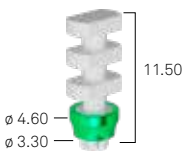
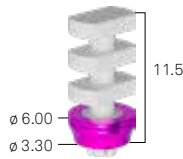
	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
straight emergence	 A-TRAP-330	-	-	-
anatomical emergence	 A-TRARP-330	 A-TRARP-380	 AS-TRARP-425	 AS-TRARP-500

Closed-tray transfers

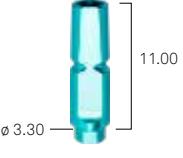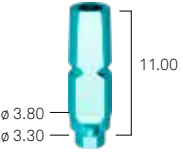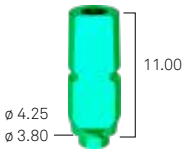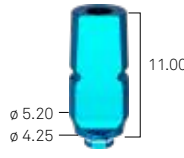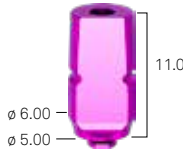	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
straight emergence	 A-TRAS-330	-	-	-
anatomical emergence	 A-TRAS-330	 A-TRAS-380	 AS-TRAS-425	 AS-TRAS-500


VTRA2-180-10 screw included.

Analogs

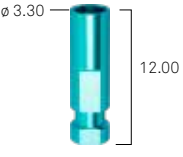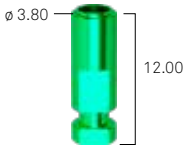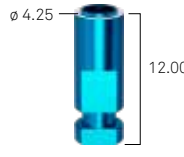	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
	 A-ANA-330	 A-ANA-380	 SH-ANA-425	 SH-ANA-500

Provisional posts


SIMPLE in PEEK with titanium base

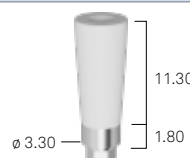


ø 3.30	ø 3.80	ø 4.25	ø 5.00	
A-MPSC-330	-	-	-	repositionable, straight emergence
A-MPSCR-330	A-MPSCR-380	AS-MPSCR-425	AS-MPSCR-500	repositionable, anatomical emergence

VM2-180 fixation screw included.


SIMPLE in titanium

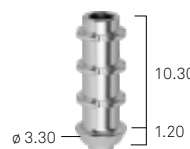

ø 3.30	ø 3.80	ø 4.25	ø 5.00	
A-MPSA-330-EX	A-MPSA-380-EX	-	-	repositionable, anatomical emergence
A-MPSA-330	A-MPSA-380	-	-	rotating, anatomical emergence

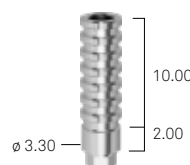

ø 3.30	ø 3.80	ø 4.25	ø 5.00	
A-MPSCI-330-EX	A-MPSCI-380-EX	AS-MPSCI-425-EX	AS-MPSCI-500-EX	repositionable, straight emergence
A-MPSCI-330	A-MPSCI-380	AS-MPSCI-425	AS-MPSCI-500	rotating, straight emergence

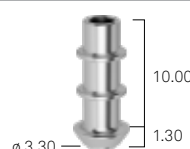

ø 3.30	ø 3.80	ø 4.25	ø 5.00	
A-MPS-330	A-MPS-380	-	-	rotating, wide emergence

VM2-180 fixation screw included.

In titanium with castable sleeve



ø 3.30	ø 3.80	ø 4.25	ø 5.00	
A-CTI-330	A-CTI-380	-	-	rotating
A-CCI-S	A-CCI-S	-	-	spare castable PMMA sleeve, fixation screw not included

Vite VM2-180 inclusa.


Posts for luting technique



unique platform

L-CT-340-ROT Gr. 5 titanium posts, h 14.00 mm, straight emergence

L-CCRCO-340-ROT Cobalt chrome titanium posts, h 14.00 mm, straight emergence

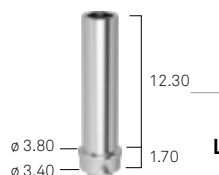







L-CTR-340-ROT Gr. 5 titanium posts, h 14.00 mm, anatomical emergence

L-CCRCOR-340-ROT Cobalt chrome posts, h 14.00 mm, anatomical emergence

L-VMS-180 screw included.

Premade posts

Straight

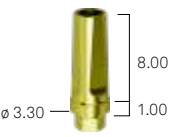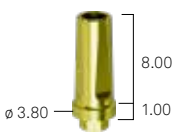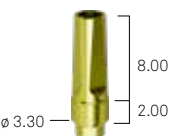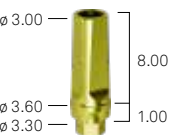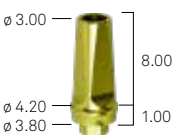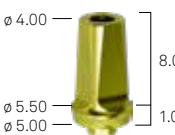straight emergence	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
transgingival h 1.00	 A-MD-330-1	 A-MD-380-1	 AS-MD-425-1	 AS-MD-500-1
transgingival h 2.00	 A-MD-330-2	 A-MD-380-2	 AS-MD-425-2	 AS-MD-500-2
transgingival h 4.00	 A-MD-330-4	 A-MD-380-4	 AS-MD-425-4	 AS-MD-500-4
anatomical emergence				
transgingival h 1.00	 A-MDR-330-1	 A-MDR-380-1	 AS-MDR-425-1	 AS-MDR-500-1
transgingival h 2.00	 A-MDR-330-2	 A-MDR-380-2	 AS-MDR-425-2	 AS-MDR-500-2
transgingival h 4.00	 A-MDR-330-4	 A-MDR-380-4	 AS-MDR-425-4	 AS-MDR-500-4

VM2-180 fixation screw included.

Angled

	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
angled at 15°, straight emergence transgingival h 1.75	 A-MA15-330	-	-	-
angled at 15°, anatomical emergence transgingival h 1.80	 A-MAR15-330	 A-MAR15-380	 AS-MAR15-425	 AS-MAR15-500
angled at 25°, anatomical emergence transgingival h 1.80	-	 A-MAR25-380	 AS-MAR25-425	-

VM2-180 fixation screw included.





Shoulderless premade straight posts

	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
repositionable h 15.00			-	-
	A-MD-330-15	A-MD-380-15		
rotating h 15.00			-	-
	A-MD-330-15-ROT	A-MD-380-15-ROT		

L-VMS-180 fixation screw included.

Prosthetic solutions for vertical technique

For vertical technique


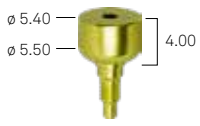
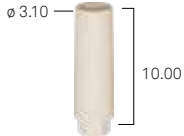
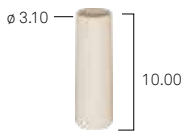
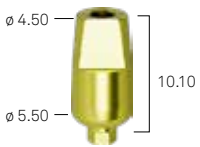
	unique platform
Gr. 5 titanium healing abutments	
	A-TMG-MEFL-330
REEF resin posts, repositionable*	
	A-PPF-330-EX
REEF resin posts, rotating*	
	A-PPF-330
Gr. 5 titanium preparable posts	
	A-MEFL-330
impression cap pack of 5 pcs.	
	CAP-MEFL-5

Resin and preparable posts: fixation screw VM2-180 included.

*Not available in USA






Preparable posts

Straight

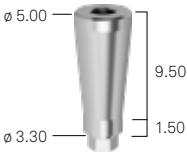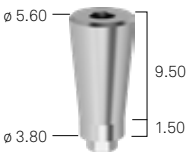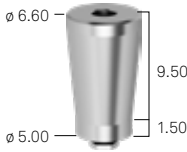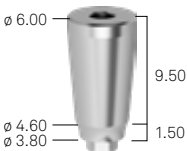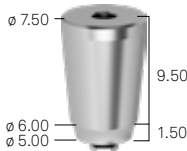	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
straight emergence	 A-MF-330	 A-MF-380	 AS-MF-425	 AS-MF-500
anatomical emergence	 A-MFR-330	 A-MFR-380	 AS-MFR-425	 AS-MFR-500





VM2-180 fixation screw included.

Precut

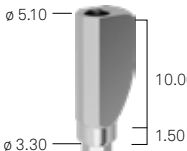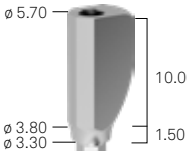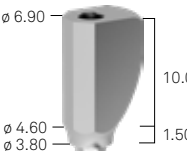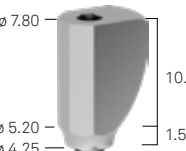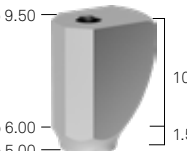	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
straight emergence	 A-MFP-330	-	-	-
anatomical emergence	 A-MFPR-330	 A-MFPR-380	 AS-MFPR-425	 AS-MFPR-500

VM2-180 fixation screw included.





SIMPLE

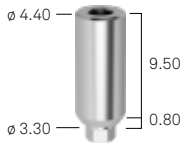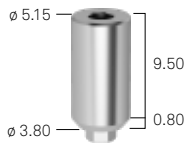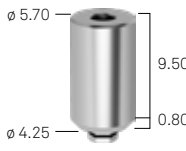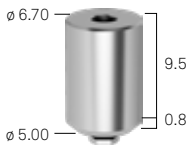	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
very wide emergence	 A-MFS-330	 A-MFS-380	 AS-MFS-425	 AS-MFS-500

VM2-180 fixation screw included.










XA prosthetic range - Cemented prosthesis

Posts for XA cemented restorations

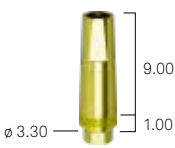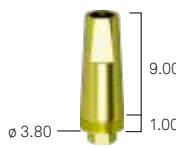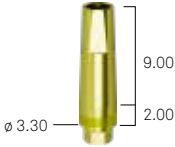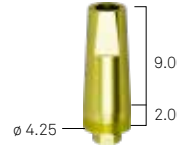transgingival h	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
1.00 mm	 SH-MD-F-330-1	 SH-MD-F-380-1	 SH-MD-F-425-1	-
2.00 mm	 SH-MD-F-330-2	 SH-MD-F-380-2	 SH-MD-F-425-2	-

L-VMS-180 fixation screw included.

Components for XA cemented restorations

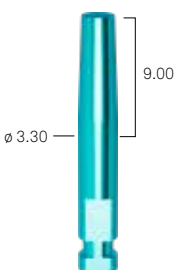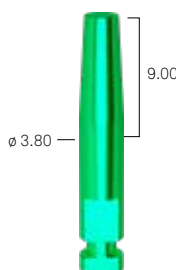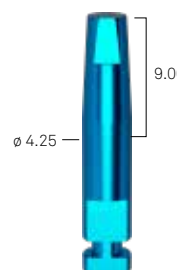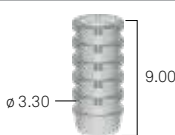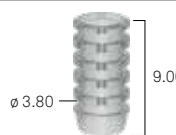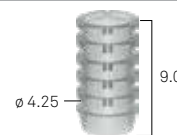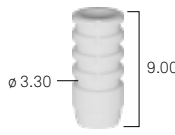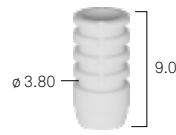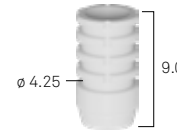	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
analog	 SH-ANA-MD-F-330	 SH-ANA-MD-F-380	 SH-ANA-MD-F-425	-
transfer	 SH-TRA-MD-F-330	 SH-TRA-MD-F-380	 SH-TRA-MD-F-425	-
PEEK caps	 SH-CP-MD-F-330	 SH-CP-MD-F-380	 SH-CP-MD-F-425	-

XA prosthetic range - Screw retained prosthesis

Posts for screw retained restorations

h $\varnothing$ 3.30 $\varnothing$ 3.80



A-ABU-F-TS-330-4 SH-ABU-F-TS-380-4



A-ABU-F-TS-330-5 SH-ABU-F-TS-380-5




A-ABU-F-TS-330-6 SH-ABU-F-TS-380-6

H 4.50 mm post: mm: SH-VABU180-200-4 screw included.
H 5.50 mm post: SH-VABU180-200-5 screw included.
H 6.50 mm post: SH-VABU180-200-6 screw included.

XA components for digital prosthesis


unique platform

scanbodies for laboratory
and intraoral use in titanium
for XA posts

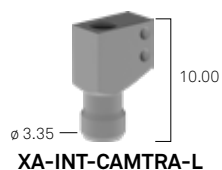


XA-INT-CAMTRA-L

digital analog for XA
abutment

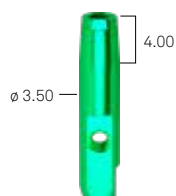

XA-ANABU-DIG

titanium premilled posts for
XA abutment
XA-CAMT-ABU-12-SW

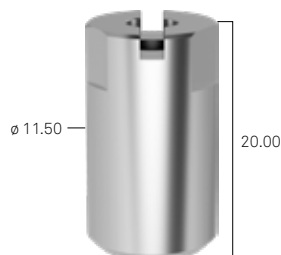

cobalt chrome premilled
posts for XA abutment
XA-CAMT-ABU-12-SW

T-Connect for multiple
crowns, XA abutment,
prof. 1.0



XA-TCIN-10-ROT

T-Connect for multiple
crowns, XA abutment



XA-TCIN-ROT

T-Connect for single crown
XA abutment




XA-TCIN

Laboratory components for XA screw retained restorations


unique platform

analog

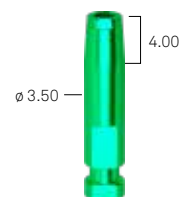


SH-ANABU-F-380

transfers

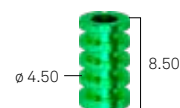

SH-TRABU-F-380

castable sleeves,
repositionable

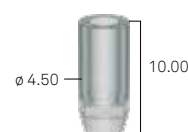

SH-CCABU-F-380

castable sleeves,
rotating




SH-CCABU-F-380-ROT

PEEK healing cap



SH-CG-ABU-F-380

Gr. 5 titanium sleeves,
repositionable

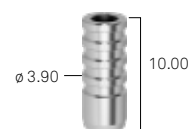

SH-CTABU-F-380

Gr. 5 titanium sleeves,
rotating



SH-CTABU-F-380-ROT

caps for luting technique



CAP-TS-DEF

Transfer: SH-VTRABU-F-200 screw included.
Castable and titanium sleeves: A-PLAIN-VP200 screw included.

Castable posts

With a metal base

	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
cobalt chrome base, repositionable	 A-UCRCO-330-EX	 A-UCRCO-380-EX	 AS-UCRCO-425-EX	 AS-UCRCO-500-EX
cobalt chrome base, rotating	 A-UCRCO-330	 A-UCRCO-380	 AS-UCRCO-425	 AS-UCRCO-500
titanium base, rotating <i>only for titanium fusions</i>	 A-UCTR-330-EX	 A-UCTR-380-EX	-	-
spare castable sleeve, fixation screw not included	 A-CCUCR-330	 A-CCUCR-380	 AS-CCUCR-425	 AS-CCUCR-500

VM2-180 fixation screw included.


Fully castable

	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
repositionable, straight emergence	 A-CC-330-EX	 A-CC-380-EX	-	-
repositionable, anatomical emergence	 A-CCR-330-EX	 A-CCR-380-EX	 AS-CCR-425-EX	 AS-CCR-500-EX
rotating, straight emergence	 A-CC-330	 A-CC-380	-	-
rotating, anatomical emergence	 A-CCR-330	 A-CCR-380	 AS-CCR-425	 AS-CCR-500

VM2-180 fixation screw included.

Solutions with angled screw hole

These solutions have a distinctive morphology, featuring a lateral window or an adjustable element that allows for the displacement of the hole for the passing screw to either the palatal or lingual position, thus avoiding unaesthetic vestibular holes.





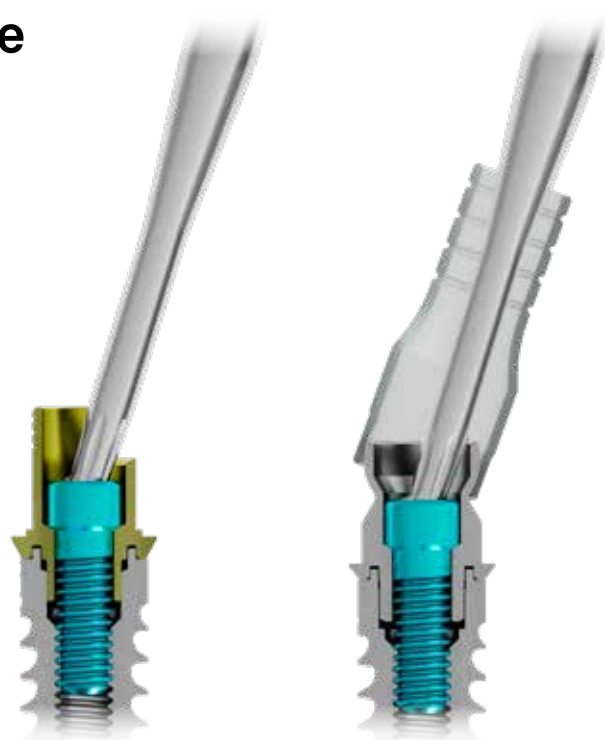


Angled screw hole Interfase

	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
rotating	 IND3PKR330/TIA	 IND3PKR380/TIA	-	-
repositionable	 IND3PKH330/TIA	 IND3PKH380/TIA	-	-

A-VMA-180 fixation screw not included.









Angled screw hole Abutment

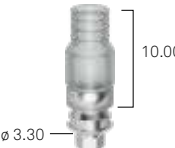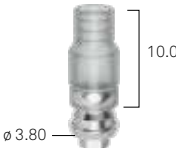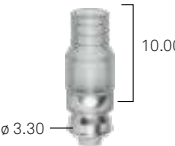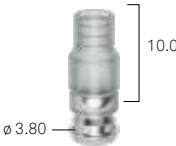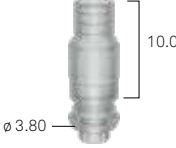	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
cobalt chrome base, repositionable	 PD3PKH330/CC	 PD3PKH380/CC	-	-
cobalt chrome base, rotating	 PD3PKR330/CC	 PD3PKR380/CC	-	-
fully castable, repositionable	 PD3PKH330/P	 PD3PKH380/P	-	-
fully castable, rotating	 PD3PKR330/P	 PD3PKR380/P	-	-

Vite A-VMA-180 fixation screw not included.



Prosthetic line PAD^e (Aesthetic screw-retained disparallel prosthesis)

Straight PADe abutments, for direct screwing on implants

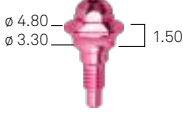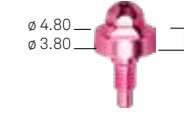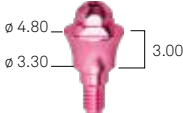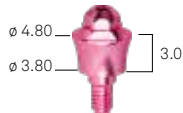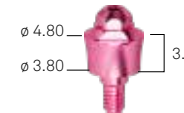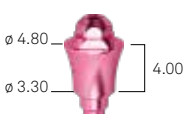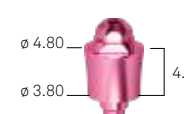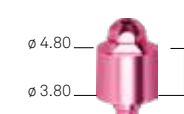transgingival h	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
1.50 mm	 A-PADE-AD330-15	 A-PADE-AD380-15	 AS-PADE-AD425-15	 AS-PADE-AD500-15
3.00 mm	 A-PADE-AD330-30	 A-PADE-AD380-15	 AS-PADE-AD425-30	 AS-PADE-AD500-30
4.00 mm	 A-PADE-AD330-40	 A-PADE-AD380-40	 AS-PADE-AD425-40	 AS-PADE-AD500-40

Packaged with dedicated carrier: for intraoral screwing, use the specific driver AVV2-ABUT.

Angled abutments

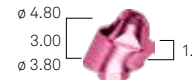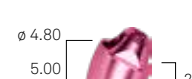transgingival h	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
Angled at 17° transgingival h 3.00	 A-PADE-AA330-173	 A-PADE-AA380-173	 AS-PADE-AA425-173	 AS-PADE-AA500-173
Angled at 17° transgingival h 5.00	 A-PADE-AA330-175	 A-PADE-AA380-175	 AS-PADE-AA425-175	 AS-PADE-AA500-175
Angled at 30° transgingival h 3.00	 A-PADE-AA330-303	 A-PADE-AA380-303	 AS-PADE-AA425-303	 AS-PADE-AA500-303
Angled at 30° transgingival h 5.00	 A-PADE-AA330-305	 A-PADE-AA380-305	 AS-PADE-AA425-305	 AS-PADE-AA500-305
Angolati a 45° transgingival h 3.00	 A-PADE-AA330-453	 A-PADE-AA380-453	-	-
Angolati a 45° transgingival h 5.00	 A-PADE-AA330-455	 A-PADE-AA380-455	-	-

PADR-VM-180 screw included.

PADR-VTRAL-140-M carrier included: for intraoral screwing, use the dedicated screwdriver from the HSM series.

P.A.D.r prosthetic components

Healing caps

Titanium healing cap



PADR-CG

PEEK healing cap



PADR-CGP

PADR-CG: PAD-VP-140 screw included.

PADR-CGP: PAD-VCGP-140 screw included.

Components for CAD-CAM prosthesis

Intraoral scanbody



PADR-INT-CAMTRA-L

Intraoral scanbody
for direct screw connection*



Y-AVV-SB-DIG-NX

Intraoral scanbody
for direct screw connection, long*



Y-AVV-SBL-DIG-NX

Digital analog for PADe abutment



PADR-ANA-DIG

Sleeve for luting technique non
engaging for abutments PADe with
cone h 4



PADR-TC-M-4

Sleeve for luting technique non
engaging for abutments PADe with
cone h 8 and reduction cut



PADR-TC-M-8

Sleeve for luting technique non
engaging for abutment PADe for
angled screw hole



PADR-TCA-M

Castable sleeve for PADR-TC-M and
PADR-TC-M-8, non engaging



PADR-TC-CC-M-8

Engaging castable sleeve for
PADR-TCA-M



PADR-TC-CC-S-8

*Scanbody and sleeve for luting technique: PAD-VP-140 screw included.
non engaging sleeve for luting technique for PADe abutment for angled
screw hole: screw PAD-VPA-140 not included.
Digital analog: VADIG-180 screw included.*

**Use the dedicated screwdriver YV-CV-MAN-M.*

Impression phase

Snap-on cap for impression



PADR-CAP-EX

Pick-up transfer, non engaging



PADR-TRA

Analog



PADR-ANA

Pick up transfer: PAD-VTRAL-140 screw included

P.A.D.r. components for prosthesis relining and luting technique

PEEK sleeve, non engaging



PADR-CP

Titanium sleeve, non engaging




PADR-CT

PMMA castable sleeve



PADR-CCEM

PMMA castable sleeve,
non engaging

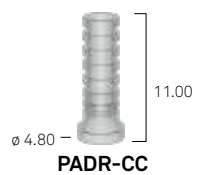

PADR-CC

Titanium sleeve without shoulder, for
aesthetic reconstructions



PADR-BAS

PEEK sleeve, Gr. 5 titanium and PMMA: PAD-VP-140 screw included.

Castable post with metal base

PMMA post with cobalt-chrome base,
rotating



PADR-UCRCO

PAD-VP-140 screw included.

PLAIN B-Space prosthetic range

Abutment

transgingival h	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
2.00 mm	 2.00 ø 3.30	 2.00 ø 3.30 ø 3.80	 2.00 ø 3.30 ø 4.25	 2.00 ø 3.30 ø 5.00
	A-PLAIN-ABU330-2	A-PLAIN-ABU-380C-2	AS-PLAIN-ABU-425C-2	AS-PLAIN-ABU-500C-2
3.00 mm	 3.00 ø 3.30	 3.00 ø 3.30 ø 3.80	 3.00 ø 3.30 ø 4.25	 3.00 ø 3.30 ø 5.00
	A-PLAIN-ABU330-3	A-PLAIN-ABU-380C-3	AS-PLAIN-ABU-425C-3	AS-PLAIN-ABU-500C-3
4.00 mm	 4.00 ø 3.30	 4.00 ø 3.30 ø 3.80	 4.00 ø 3.30 ø 4.25	 4.00 ø 3.30 ø 5.00
	A-PLAIN-ABU330-4	A-PLAIN-ABU-380C-4	AS-PLAIN-ABU-425C-4	AS-PLAIN-ABU-500C-4

Components for overstructure

unique platform	
healing cap	 5.00 ø 4.90 ø 3.30 A-PLAIN-CG330
castable sleeve	 10.00 ø 3.30 A-PLAIN-CC330
Gr. 5 titanium sleeve	 9.50 ø 3.30 1.30 A-PLAIN-CT330
Pick-up transfer	 12.00 ø 3.30 A-PLAIN-TRA-330
analog	 12.00 ø 3.30 A-PLAIN-ANA-330

Castable and titanium sleeves: A-PLAIN-VP200 fixation screw included.
Pick-up transfers: A-PLAIN-VTRA200 fixation screw included.

Components for digital prosthesis

unique platform	
Titanium scanbody for laboratory and intraoral use, for PLAIN B-Space abutment	 13.00 ø 3.30 A-INT-CAMTRA-PLAIN-L
Digital analog for PLAIN B-Space abutment	 12.00 ø 3.30 A-PLAIN-ANA-330-DIG
PLAIN B-Space sleeve for luting technique h. 6.00 mm	 5.50 ø 3.30 A-PLAIN-MDT-6
PLAIN B-Space sleeve for luting technique h. 10.00 mm	 10.00 ø 3.30 A-PLAIN-MDT-10

Scanbody: A-PLAIN-VP200 screw included.
Digital analog: VADIG-180 screw included.
PLAIN B-Space sleeves for luting technique: screw not included.
A-PLAIN-VP200 screw for straight hole or ABU-VMA-200 for angled screw hole available for separate ordering.




Conico conometric prosthesis

Conico straight posts

h	ø 3.30 ●	ø 3.80 ●
0.50	 A-MD-TS-330-05	-
1.00	 A-MD-TS-330-10	 A-MD-TS-380-10
2.00	 A-MD-TS-330-20	 A-MD-TS-380-20
3.00	 A-MD-TS-330-30	 A-MD-TS-380-30



VM2-180 fixation screw included.

Conico angled posts

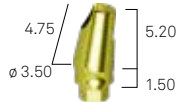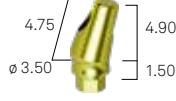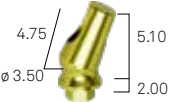
	ø 3.30 ●	ø 3.80 ●
angled at 7.5°	 A-MA07-TS-330-2	 A-MA07-TS-380-2
angled at 15°	 A-MA15-TS-330-2	 A-MA15-TS-380-2
angled at 22.5°	 A-MA20-TS-330-2	 A-MA20-TS-380-2

VM2-180 fixation screw included.

Conico components for digital prosthesis

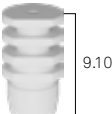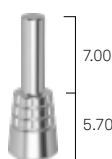
cap for partial or full-arch fixed prosthesis	 CAP2-TS-DEF
non-rotating cap for single element	 CAP2-TS-IND
cap for the removal	 CAP2-TS-REM
PEEK transfer cap for direct impression on Conico abutment	 CAP2-TS-IMP
Conico post analog	 ANA2-MD-TS
cap with pin, to be used for intraoral welding (welding bars available upon request)	 CAP2-TS-PIN

Accessories for Conico system

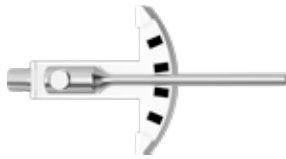parallelometer	 PAR-PP
parallelometer insert	 PAR-INS-AS
parallelometer driver	 CAP2-TS-PAR

Remember to always purchase the parallelometer with the insert relating to the implant connection on which it is to be used.

Conico components for digital prosthesis

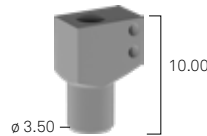unique platform	
Titanium scanbody for laboratory and intraoral use, for Conico abutments	 C-INT-CAMTRA-L
Digital analog for Conico abutment	 ANA2-MD-TS-DIG

Digital analog: VADIG_180 screw included.

Components for digital prosthesis

Scanbody for laboratory and intraoral use, titanium



Fixation screw VM2-180 included.

Digital analog



VADIG-180 screw included


Components for special technique

for P.A.D.r abutment

for PLAIN B-Space abutment

for XA posts

for Conico system

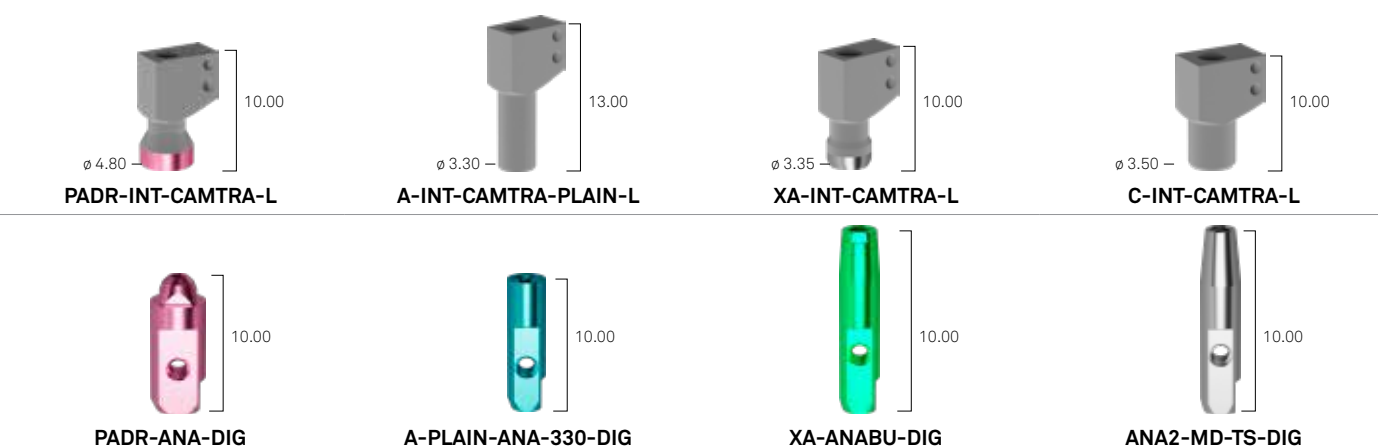

Scanbody for PLAIN abutment and scanbody for XA posts: A-PLAIN-VP200 screw included.

Scanbody for P.A.D.r abutment: PAD-VP-140 screw included.

Digital analog: VADIG-180 screw included.

Screws and screwdrivers for digital analog

unique platform

Screw for digital analog M 1.8



VADIG-180

Pin for digital analog



PIN-ANA-DIG

Carrier to insert the digital analog M1.8












INS-ANA-180

Carrier to insert the digital analog M 1.8, long




INS-ANA-180-L

T-Connect with with 4.00 mm cementing cone height

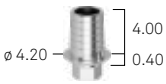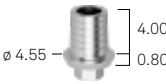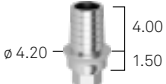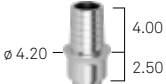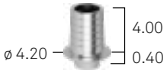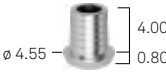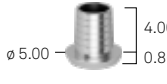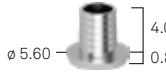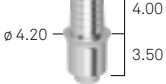	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
				
	A-BASTZR-S-330-4	A-BASTZR-S-380-4	AS-BASTZR-S-425-4	AS-BASTZR-S-500-4
repositionable			-	-
	A-TCR-S-330-4-1	A-TCR-S-380-4-1		
			-	-
	A-TCR-S-330-4-2	A-TCR-S-380-4-2		
			-	-
	A-TCR-S-330-4-3	A-TCR-S-380-4-3		
rotating				
	A-BASTZR-M-330-4	A-BASTZR-M-380-4	AS-BASTZR-M-425-4	AS-BASTZR-M-500-4
			-	-
	A-TCR-M-330-4-1	A-TCR-M-380-4-1		
			-	-
	A-TCR-M-330-4-2	A-TCR-M-380-4-2		
			-	-
	A-TCR-M-330-4-3	A-TCR-M-380-4-3		

VM2-180 fixation screw included.

Sleeves for wax-up modeling on T-Connect

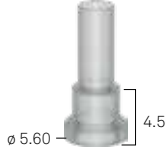	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
total h 12.00 mm				
	A-CCBAS-330-4	A-CCBAS-380-4	A-CCBAS-425-4	A-CCBAS-500-4

Angled screw hole Interface

	ø 3.30 ●	ø 3.80 ●
rotating		
	IND3PKR330/TIA	IND3PKR380/TIA
repositionable		
	IND3PKH330/TIA	IND3PKH380/TIA

A-VMA-180 fixation screw not included.

	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
	 A-BASTZR-S-330-6	 A-BASTZR-S-380-6	 AS-BASTZR-S-425-6	 AS-BASTZR-S-500-6
repositionable	 A-TCR-S-330-6-1	 A-TCR-S-380-6-1	-	-
	 A-TCR-S-330-6-2	 A-TCR-S-380-6-2	-	-
	 A-TCR-S-330-6-3	 A-TCR-S-380-6-3	-	-
	 A-BASTZR-M-330-6	 A-BASTZR-M-380-6	 AS-BASTZR-M-425-6	 AS-BASTZR-M-500-6
rotating	 A-TCR-M-330-6-1	 A-TCR-M-380-6-1	-	-
	 A-TCR-M-330-6-2	 A-TCR-M-380-6-2	-	-
	 A-TCR-M-330-6-3	 A-TCR-M-380-6-3	-	-

VM2-180 fixation screw included.

Sleeves for wax-up modeling on T-Connect

	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
total h 12.00 mm	 ø 4.20 6.60 A-CCBAS-330-6	 ø 4.55 4.50 A-CCBAS-380-6	 ø 5.00 4.50 A-CCBAS-425-6	 ø 5.60 4.50 A-CCBAS-500-6

Unica T-Connect

ø 3.30 ●


ø 3.80 ●

ø 4.25 ●


ø 5.00 ●




A-TC-UNI-330

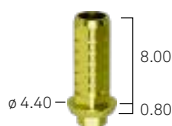


A-TC-UNI-380

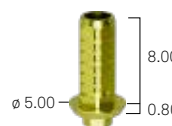


AS-TC-UNI-425

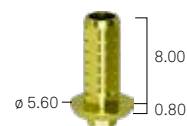


AS-TC-UNI-500

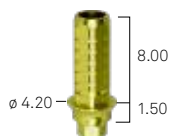


A-TC-UNI-330-1

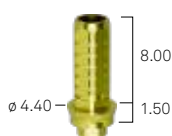

A-TC-UNI-380-1

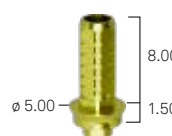


AS-TC-UNI-425-1

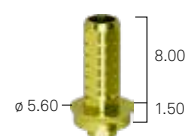


AS-TC-UNI-500-1


repositionable

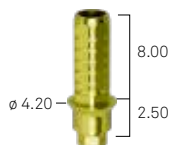

A-TC-UNI-330-2

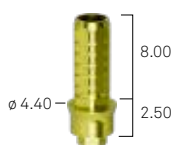


A-TC-UNI-380-2

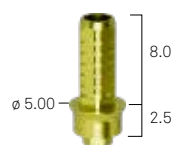


AS-TC-UNI-425-2




AS-TC-UNI-500-2

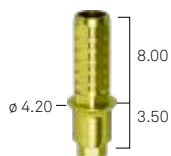


A-TC-UNI-330-3

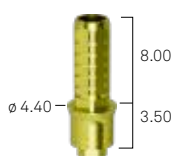

A-TC-UNI-380-3

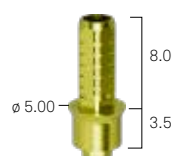

AS-TC-UNI-425-3

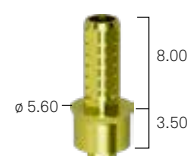

AS-TC-UNI-500-3



A-TC-UNI-330-ROT



A-TC-UNI-380-ROT



AS-TC-UNI-425-ROT



AS-TC-UNI-500-ROT



A-TC-UNI-330-1-ROT



A-TC-UNI-380-1-ROT




AS-TC-UNI-425-1-ROT



AS-TC-UNI-500-1-ROT

rotating

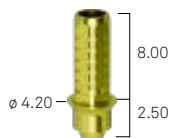

A-TC-UNI-330-2-ROT




A-TC-UNI-380-2-ROT



AS-TC-UNI-425-2-ROT




AS-TC-UNI-500-2-ROT

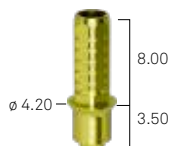

A-TC-UNI-330-3-ROT



A-TC-UNI-380-3-ROT

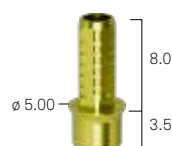

AS-TC-UNI-425-3-ROT



AS-TC-UNI-500-3-ROT

VM2-180 fixation screw included.

The angled hole screw must always be ordered separately A-VMA-180 because it is not included in the Interfase package. If you do not have the appropriate screwdriver DSPDCLH-24 or DSPDCLH-32, please order it separately.

Cutting guide and holder for Unica T-Connect



Cutting guide holder S
(small) for Unica T-connect
left side

TC-UNI-HOLD4-SX

Cutting guide holder S
(small) for Unica T-connect
right side

TC-UNI-HOLD4-DX


Cutting guide holder L (large)
for Unica T-connect right
side

TC-UNI-HOLD6-DX

Cutting guide holder L (large)
for Unica T-connect left side


TC-UNI-HOLD6-SX

4 mm cutting guide for
Unica T-Connect

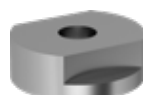

TC-UNI-DIMA-H4

6 mm cutting guide for
Unica T-Connect

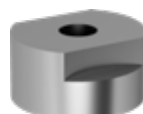

TC-UNI-DIMA-H6

Cutting guide for angled hole
for T-Connect Unica



TC-UNI-DIMA-ANG

Please note: STL file available for independently 3D printing of the holders.



See how easy it is to precisely
customize the T-Connect Unica
with the cutting guides and the
corresponding holder.

Unica T-Connect accessories

Intraoral testing O-ring for
T-Connect Unica




TC-UNI-ORING

Castable sleeve for
T-Connect, rotating



TC-UNI-CC

Small castable sleeve
for Unica T-Connect,
repositionable

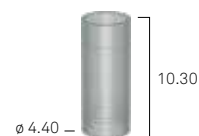

TC-UNI-CC-S

Large castable sleeve
for Unica T-Connect,
repositionable



TC-UNI-CC-L

O-ring inserter



C-UNI-ORING-INS

T-Connect for XA abutment

T-Connect for multiple
crowns, XA abutment,
prof. 1.0



XA-TCIN-10-ROT

T-Connect for multiple
crowns, XA abutment



XA-TCIN-ROT

T-Connect for single crown
XA abutment



XA-TCIN

Components for chairside digital prosthesis

Chairside T-Connect supports with 5.00 mm h cementation cone

transgingival h	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
				
	A-TCS-330S-C	A-TCS-380S-C	AS-TCS-425L-C	AS-TCS-500L-C
1.50 mm			-	-
	A-TCS-330S-C-15	A-TCS-380S-C-15		
2.50 mm			-	-
	A-TCS-330S-C-25	A-TCS-380S-C-25		

VM2-180 fixation screw included.

Premilled

for implants, in titanium

	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
repositionable				
titanium	A-CAMT-330-12-SW	A-CAMT-380-12-SW	AS-CAMT-425-12-SW	AS-CAMT-500-12-SW
cobalt-chrome	-	A-CAMC-380-12-SW	AS-CAMC-425-12-SW	AS-CAMC-500-12-SW

VM2-180 fixation screw included.

Platform level

	piattaforma unica	
repositionable		
titanium	A-CAMTABU330-CRE	A-CAMTABU380-CRE

VM2-180 fixation screw included.

for XA abutment, in titanium

	unique platform
	
titanium	XA-CAMT-ABU-12-SW
cobalt-chrome	XA-CAMC-ABU-12-SW

A-PLAIN-VP200 fixation screw included.

Solutions for overdenture

Locator abutment

transgingival h	ø 3.30 ●	ø 3.80 ●
1.00 mm	 ø 3.30 — 1.00 1670	 ø 3.80 — 1.00 1675
2.00 mm	 ø 3.30 — 2.00 1671	 ø 3.80 — 2.00 1676
3.00 mm	 ø 3.30 — 3.00 1672	 ø 3.80 — 3.00 1677
4.00 mm	-	 ø 3.80 — 4.00 1678
5.00 mm	 ø 3.30 — 5.00 1674	 ø 3.80 — 5.00 1679
6.00 mm	-	 ø 3.80 — 6.00 1680

Accessories for Locator abutment

	unique platform
4 aluminium analogs	 8530
4 aluminium transfers, 4 retainers cod. 8515 included	 8505
4 polyethylene parallelism pins	 8517
steel plate AISI 316 L	 9530
complete Locator Core Tool	 8393
steel push rod for inserting the retainers	 8397
steel driver for abutment screwing/unscrewing	 8390
retention jacket for the driver, to transfer the Locator abutments into the oral cavity	 8394

Retainers and spacer rings for Locator abutments

	unique platform
low retention capacity for impression taking, black, in polyethylene, 4 pcs packaging	 8515
retention 5 lb = 2268 g, in nylon, disparallelisms up to 10° per side, 4 pcs packaging	 8524
retention 3 lb = 1361 g, in nylon, disparallelisms up to 10° per side, 4 pcs packaging	 8527
retention 1,5 lb = 680 g, in nylon, disparallelisms up to 10° per side, 4 pcs packaging	 8529
retention 4 lb = 1814 g, in nylon, disparallelisms up to 20° per side, 4 pcs packaging	 8547
retention 2 lb = 907 g, in nylon, disparallelisms up to 20° per side, 4 pcs packaging	 8915
retention 1 lb = 450 g, in nylon, disparallelisms up to 20° per side, 4 pcs packaging	 8548
spacer rings in silicon rubber, 20 pcs packaging	 8514
kit: 2 Gr. 5 titanium caps, 2 spacer rings, 2 black retainers and 2 nylon retainers for each of the 3 different retention capacities	 8519-2
kit: 2 Gr. 5 titanium caps, 2 spacer rings, 2 black retainers and 2 nylon retainers for each of the 3 different retention capacities meant for disparallelisms	 8540-2
kit: 2 Gr. 5 steel caps, 2 spacer rings, 2 black retainers and 2 nylon retainers for each of the 3 different retention capacities	 8550-2

Screwdrivers for locator abutment

	unique platform
short	 8926-SW
long	 8927-SW

Locator R-Tx abutment*

transgingival h	ø 3.30 	ø 3.80 
1.00 mm	 ø 3.30 — 1.00 31110-01	 ø 3.80 — 1.00 31111-01
2.00 mm	 ø 3.30 — 2.00 31110-02	 ø 3.80 — 2.00 31111-02
3.00 mm	 ø 3.30 — 3.00 31110-03	 ø 3.80 — 3.00 31111-03
4.00 mm	 ø 3.80 — 4.00 31110-04	 ø 3.80 — 4.00 31111-04
5.00 mm	 ø 3.30 — 5.00 31110-05	 ø 3.80 — 5.00 31111-05
6.00 mm	 ø 3.80 — 6.00 31110-06	 ø 3.80 — 6.00 31111-06

* All-in-One packaging that includes an abutment, a spacer ring, a cap, and an insert for each retention (gray, blue, pink, white).

Caps and retainer for Locator R-Tx abutment

	unique platform
pack of 4 caps and 4 black retainers for impression taking	 30013
grey, without retention, pack of 4 pcs	 30001
blue, low retention, pack of 4 pcs	 30002
pink, medium retention, pack of 4 pcs	 30003
white, high retention, pack of 4 pcs	 30004
black, for the impression taking, pack of 4 pcs	 30012


Accessories for Locator R-Tx abutment

	unique platform
aluminium analogs for abutments ø 3.35 mm, pack of 4 pcs	 30014
aluminium analogs for abutments ø 4.00, pack of 4 pcs	 30015
aluminium analogs for abutments ø 5.00, pack of 4 pcs	 30016
transfers with low retention, pack of 4 pcs	 30017
spacer rings for relining phase, pack of 4 pcs	 30018
Locator core tool for Locator R-Tx	 30021
RA screwdriver for Locator R-Tx, hex 1.25 mm	 08008
digital screwdriver for Locator R-Tx, hex 1.25 mm	 08007


Ball attachments

ø 3.30 ●


ø 3.80 ●

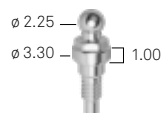


A-AS-330-1

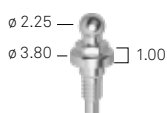


A-AS-380-1

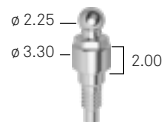


A-AS-330-2

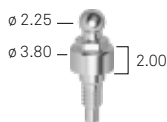

A-AS-380-2

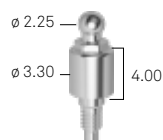

A-AS-330-4

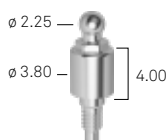

A-AS-380-4

Analog for ball attachments

unique platform



ANAS

Caps for ball attachments

unique platform

polyamide cap



040CRNSNDR8-SM*

steel container for polyamide cap external diameter ø 4.80 mm



041CAN-SM*

Gr. 5 titanium cap with cap in two parts, retention spring titanium and plastic, plastic mounting ring



CAP-TIT-1

spare plastic ring



AN-CAP-TIT-1

spare steel spring, average hardness



MOL1-CAP-TIT-1

spare steel spring, soft



MOL2-CAP-TIT-1

gold alloy cap, with 3 plastic rings and a tin spacer for the lab



CAP-1

* Pack of 6 pieces

Bars

castable bar, l 5.00 cm, h 3.00 mm, thickness 2.20 mm, ovoid-shaped profile with spacer



BARC-CAV-TIT

divisible bar attachment in titanium for oval bars, h 3.00 mm, thickness 2.20 mm



CAV-TIT

castable bar, l 5.00 cm, ø 2.20 mm



BARC

bar attachment in gold alloy for round bars ø 2.20 mm



CAV-375

Clinical and laboratory prosthetic screws

Standard fixation screws

	VM2-180	
	VM2-180-10	10 pcs packaging

Conical support fixation screws

	L-VMS-180	
	L-VMS-180-10	10 pcs packaging

Fixation screw for angled screw hole

	A-VMA-180
---	------------------

Screws for Pick-up transfer

	VTRA2-180-15	h 15.00
	VTRA2-180-20	h 20.00
	VTRA2-180-MAN	for manual screwing

Screws for closed-tray transfers

	VTRA2-180-10	h 10.00
--	---------------------	---------

Screws for XA technique

	SH-VABU180-200-4	for post h. 4.50 for screw-retained prosthesis
	SH-VABU180-200-5	for post h. 5.50 for screw-retained prosthesis
	SH-VABU180-200-6	for post h. 6.50 for screw-retained prosthesis
	A-PLAIN-VP200	for cemented prosthesis
	L-VMS-180	for cemented prosthesis
	SH-VTRABU-F-200	transfer screw for XA abutment

Recommended torque for tightening prosthetic components

8-10 Ncm	surgical cover screws, healing abutments, transfer screws
20-25 Ncm	prosthetic screws
25-30 Ncm	direct screwing on implants prosthetic components
8-10 Ncm	transfer screws

Warning

Fixation screws are always included in the packaging of each prosthetic component, unless otherwise indicated. Please use laboratory screws for laboratory phases and keep the screw included with the prosthetic component for the final fastening in the patient's mouth.

Screws for PAdE abutments

	PADR-VM-180	for PAdE abutment
	PADR-VM-180-10	for PAdE abutment 10 pcs packaging
	PAD-VTRAL-140	for Pick-up transfer, long
	PAD-VTRA-140	for Pick-up transfer, standard
	PAD-VTRAL-140-MAN	for Pick-up transfer, manual screwing
	PAD-VP-140	for prosthetic components on PAdE abutments
	PAD-VP-140-10	for prosthetic components on PAdE abutments 10 pcs packaging
	PAD-VPA-140	fixation screw for PADR-TCA-M to use with the dedicated screwdriver DSPDCLH-24 and DSPDCLH-32
	PADR-VTRAL-140-M	carrier for PAdE

Screws for CAD-CAM prosthesis

	A-PLAIN-VP200	for PLAIN abutment
	A-PLAIN-VP200-10	for PLAIN abutment, 10 pcs packaging
	A-PLAIN-VTRA200	for PLAIN abutment Pick-up transfers
	ABU-VMA-200	for angled screw hole

Screws for CAD-CAM prosthesis


	A-CAMTVABU180	for zirconium oxide single crowns and ECHO overstructures
	PAD-VCAM140	for zirconium oxide overstructure on abutment PAdE
	CAMPRON205-10	spare ring washers for the head of the screw
	VADIG-180	locking screw for digital analogs

Surgical sequences

Implant Site Preparation Sequences for Premium Power Implants

The illustrated sequence shown is an example of site preparation for a Premium Power implant with a length of 11.50 mm and a diameter of $\varnothing$ 5.00 mm.

The use of the drill stops corresponding to the desired preparation depth, included in the kit, is also recommended.

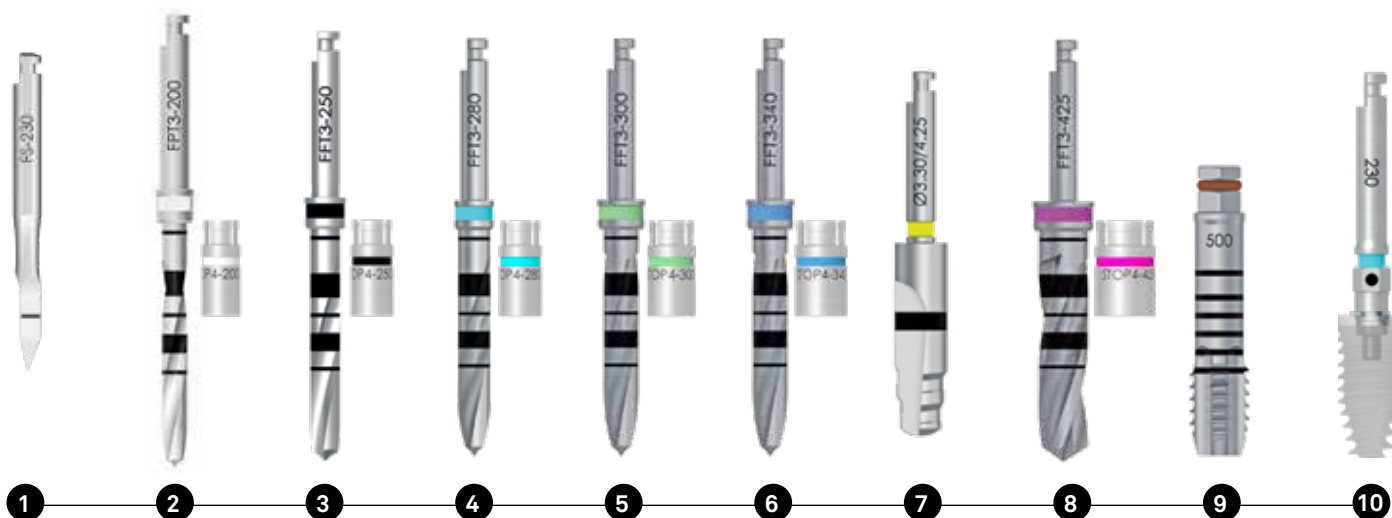


The instrument sequence shown below is applicable to all implants in the Premium Power range.

		$\varnothing$ 3.30 mm ●	$\varnothing$ 3.80 mm ●	$\varnothing$ 4.25 mm ●	$\varnothing$ 5.00 mm ●
1	Precision drill	✓	✓	✓	✓
2	Cylindrical pilot drill $\varnothing$ 2.00 mm	✓	✓	✓	✓
3	Cylindrical final drill $\varnothing$ 2.50 mm	✓	✓	✓	✓
4	Cylindrical final drill $\varnothing$ 2.80 mm	✓	✓	✓	✓
5	Cylindrical final drill $\varnothing$ 3.00 mm		✓	✓	✓
6	Cylindrical final drill $\varnothing$ 3.40 mm			✓	✓
7	Intermediate drill $\varnothing$ 3.30/4.25 mm				✓
8	Cylindrical final drill $\varnothing$ 4.25 mm				✓
9	Bone tap	SM-MS-330-CA	A-MS-380	A-MS-425	A-MS-500
10	Easy-Insert short driver	✓	✓	✓	✓

Implant Site Preparation Sequences for Premium Power Implants with Shorty drills

The illustrated sequence shown is an example of site preparation for a Premium Power implant with a length of 7 mm and a diameter of $\varnothing$ 5.00 mm.

The use of the drill stops corresponding to the desired preparation depth, included in the kit, is also recommended.

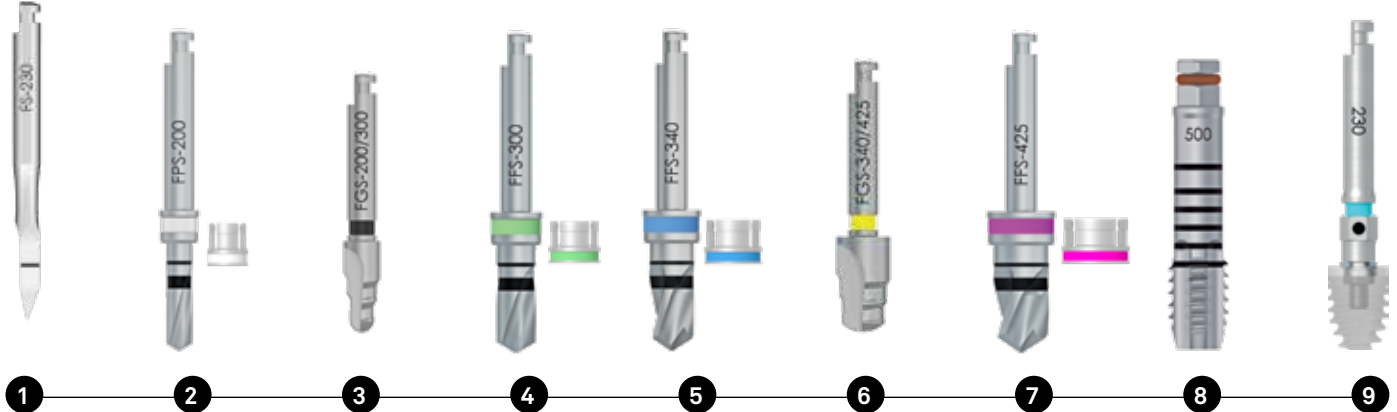

The instrument sequence shown below is applicable to all implants in the Premium Power range with 7 mm height.

	$\varnothing$ 4.25 mm ●	$\varnothing$ 5.00 mm ●
1 Precision drill	✓	✓
2 Cylindrical pilot drill $\varnothing$ 2.00 mm	✓	✓
3 Cylindrical final drill $\varnothing$ 2.50 mm	✓	✓
4 Cylindrical final drill $\varnothing$ 3.00 mm	✓	✓
5 Cylindrical final drill $\varnothing$ 3.40 mm	✓	✓
6 Intermediate drill $\varnothing$ 3.40/4.25 mm	✓	✓
7 Cylindrical final drill $\varnothing$ 4.25 mm	✓	✓
8 Bone tap	A-MS-425	A-MS-500
9 Easy-Insert short driver	✓	✓

Implant Site Preparation Sequences for Shelta Power Implants

The illustrated sequence shown is an example of site preparation for a Shelta Power implant with a length of 11.50 mm and a diameter of $\varnothing$ 5.00 mm.

The use of the drill stops corresponding to the desired preparation depth, included in the kit, is also recommended.



The instrument sequence shown below is applicable to all implants in the Shelta Power range.

	$\varnothing$ 3.80 mm ●	$\varnothing$ 4.25 mm ●	$\varnothing$ 5.00 mm ●
1 Precision drill	✓	✓	✓
2 Cylindrical pilot drill $\varnothing$ 2.00 mm	✓	✓	✓
3 Cylindrical final drill $\varnothing$ 2.50 mm	✓	✓	✓
4 Conical final drill $\varnothing$ 3.80 mm	✓	✓	✓
5 Conical final drill $\varnothing$ 4.25 mm		✓	✓
6 Conical final drill $\varnothing$ 5.00 mm			✓
7 Bone tap	SH-MS-380	SH-MS-425	SH-MS-500
8 Easy-Insert short driver	✓	✓	✓



rev. 06-26v.01



Sweden & Martina S.p.A.
Via Veneto, 10
35020 Due Carrare (PD), Italy
Tel. +39.049.9124300
info@sweden-martina.com
www.sweden-martina.com

Sweden & Martina Ltd
Unit 1b
Amberley Court, Whitworth Road
Crawley, West Sussex, RH11 7XL
Toll free 0800 1123575
info.uk@sweden-martina.com

Authorized Distributor
Sweden & Martina Inc.
78 John Miller Way
Unit 1021 - Kearny, New Jersey 07032
Toll free 1-844-8MARTINA (1-844-862-7846)
info.us@sweden-martina.com

www.sweden-martinainc.com

Sweden & Martina Mediterranea S.L.
España - info.es@sweden-martina.com
Sweden & Martina Lda
Portugal - info.pt@sweden-martina.com
Sweden Martina Romania srl - România
info.ro@sweden-martina.com
Sweden Martina France Sarl - France
info.fr@sweden-martina.com

SOME PRODUCTS MAY NOT BE REGULATORY/RELEASED FOR SALE IN ALL MARKETS.

The implants, standard prosthetic components and surgical instruments contained in this issue are Medical devices and are manufactured by Sweden & Martina S.p.A. except for:

- the Locator Abutment and related components which are medical devices manufactured by Zest Anchors, Inc., 2875 Loker Avenue East, Carlsbad, CA 92010, USA. Locator is a trademark by Zest Anchors, Inc.
- The Screwdrivers for angled screw hole are medical devices manufactured by Talladium España S.L., Avenida Blondel, 54 3º, 25002 Lleida, Spain.
- The TWL torque wrench with control lever and its adapter are medical devices manufactured by Elos Medtech Pinol A/S, Engvej 33, 3330 Gorlose, Denmark.
- Caps for ball attachments, medical devices manufactured by Rhein83 S.r.l., Via E. Zago 10/ABC, 40128 Bologna, Italy.

They conform to the ISO 9001 and ISO 13485 standards and are certified with the CE in compliance with Regulation (EU) Medical Devices n. 2017/745.
The Sweden & Martina plant manufactures Medical Devices in compliance with the CGMPs in force in the USA and in other countries worldwide.



All trademarks herein are the property of Sweden & Martina S.p.A. unless otherwise indicated.
This material is intended for laboratories and clinicians and is not intended for patient distribution.
This material is not to be redistributed, duplicated, or disclosed without the express written consent of Sweden & Martina S.p.A.
For additional product information, including indications, contraindications, warnings, precautions, and potential adverse effects, see Sweden & Martina S.p.A. website.

The contents are updated at the time of publication. Check with the company for any subsequent updates.