

PRAMAPOWER

International Product Overview





Sweden & Martina



Dental Arena


PRAMA POWER

Table of contents

4	Prama Slim Power implants
6	Prama Power implants
14	Prama Power Nasal and Pterygoid implants
16	UTM and ZirTi: two surfaces in synergy
19	The advantages of the unique Collex connection
20	Guidelines for the placement of Prama Power using the intramucosal protocol
24	Prama Power the implant for a complete digital workflow
26	Ergonomics of surgical procedures
28	Prama Power Surgical Kit
30	Surgical Kits for Guided Surgery
34	Components for Guided Surgery
36	Optional integration sets
38	Optional surgical instruments
42	Excellence in prosthetic solutions
43	Solutions for angled screw hole
44	Impression and model phase
45	Healing abutments
46	Prama IN healing abutments
47	Provisional posts
48	Premade posts
49	Preparable posts
49	Castable posts
50	Solutions with angled screw hole
52	Prosthetic line PAD [®] (Aesthetic screw-retained disparallel prosthesis)
54	PLAIN B-Space prosthetic range
55	Conico conometric prosthesis
56	Components for digital prosthesis
62	Premilled
63	Overdenture solutions
66	Implant site preparation sequences
70	Clinical and laboratory prosthetic screws
71	Bibliography


Prama Slim Power

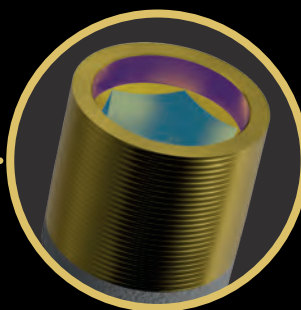SLIM $\varnothing$ 3.30




**COLLEX
CONNECTION**

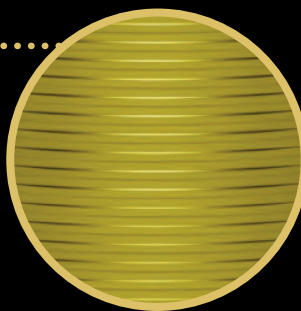For a simplified
prosthetic approach



UTM NECK

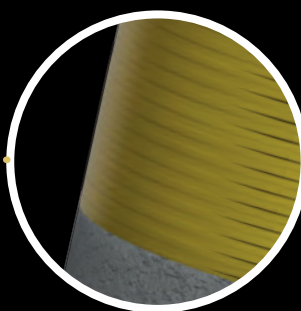The surface designed for
the interface between soft
and hard tissues




CORONAL MORPHOLOGY

More stable in thin cortical



**REVERSE BUTTRESS
THREAD**

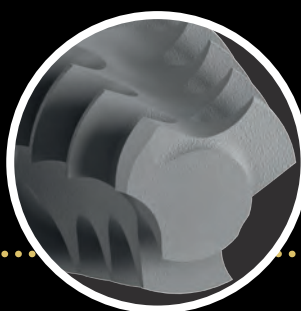Sharp and penetrating
even in dense bone



APEX

Tapered with a deeper spiral
throughout its length, ensuring
stability, especially in
post-extraction cases

Prama Slim Power implants*

ø 3.30 mm




	SHORT NECK h 1.80 mm	REGULAR NECK h 2.80 mm
h fixture		
6.00 mm	-	 LA-ZT-330-060
8.50 mm	 LAS-ZT-330-085	 LA-ZT-330-085
10.00 mm	 LAS-ZT-330-100	 LA-ZT-330-100
11.50 mm	 LAS-ZT-330-115	 LA-ZT-330-115
13.00 mm	 LAS-ZT-330-130	 LA-ZT-330-130
15.00 mm	-	 LA-ZT-330-150
surgical fixation screw (included)	 L-VT-340	


* Sweden & Martina dental implants are implantable medical devices indicated for the rehabilitation of patients with partial or complete edentulism. They are intended for surgical placement in the mandibular or maxillary bone. Sweden & Martina implant systems are designed for both one-stage and two-stage surgical procedures. Due to their diameter, implants with a diameter of 3.30 mm or less are recommended in cases of reduced bone width. They may be used for the restoration of single crowns replacing incisors and canines. Due to their reduced dimensions, they are not suitable for single-tooth restorations in the premolar or molar regions. In posterior regions, these implants may only be used for bar-retained or bridge-supported restorations in combination with larger-diameter implants.

Prama Power implants

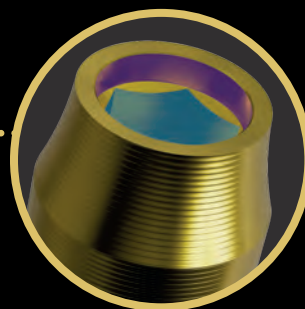
PRAMA POWER



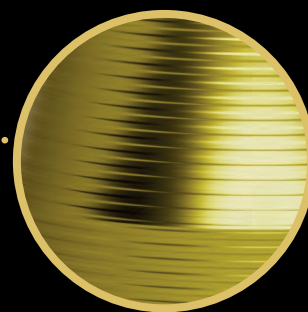
COLLEX CONNECTION
Common to all diameters for a simplified prosthetic approach



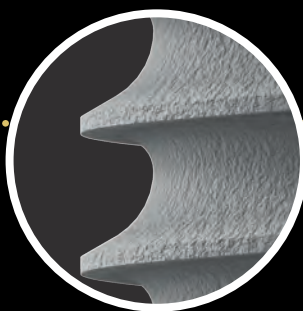UTM NECK
The surface designed for the interface between soft and hard tissues



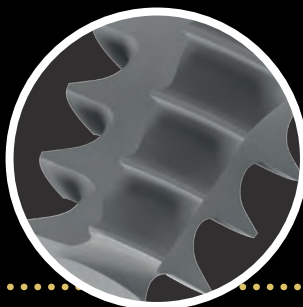REVERSE BUTRESS THREAD
Sharp and penetrating even in dense bone




APICAL NOTCHES
Designed to facilitate insertion and clot decomposition




CONICAL CORE AND CONICAL SPIRAL
Provides excellent sensitivity during insertion

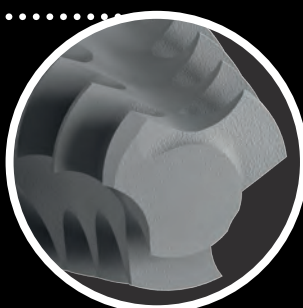APEX
Tapered with a deeper spiral throughout its length, ensuring stability, especially in post-extraction cases



PRAMA RF POWER - PRAMA RF SL POWER

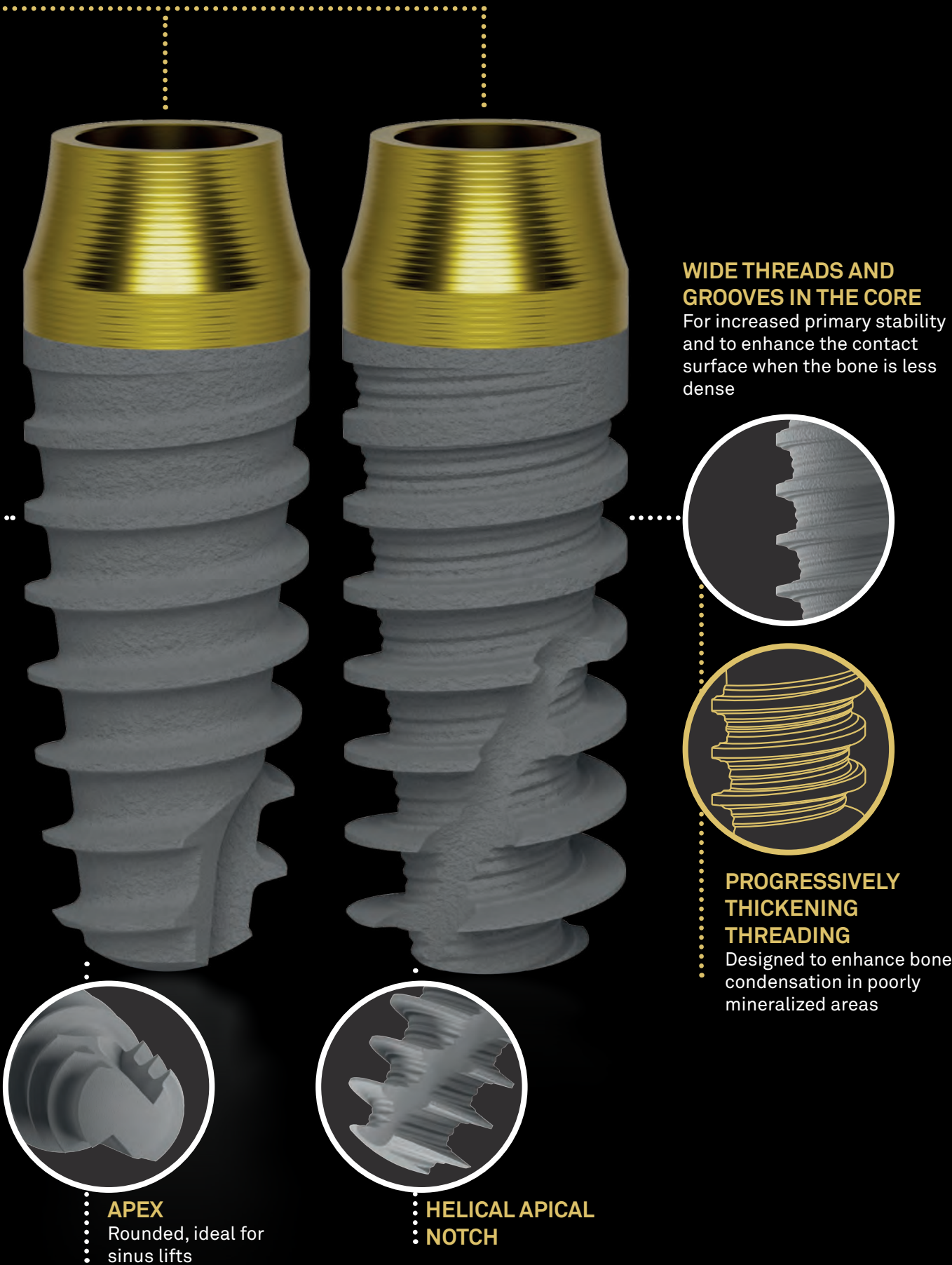

Short Neck h 1.80 mm



PRAMA POWER			
h fixture	ø 3.80 ●	ø 4.25 ●	
8.50 mm	 LAS-ZT-380-085	 LAS-ZT-425-085	
10.00 mm	 LAS-ZT-380-100	 LAS-ZT-425-100	
11.50 mm	 LAS-ZT-380-115	 LAS-ZT-425-115	
13.00 mm	 LAS-ZT-380-130	 LAS-ZT-425-130	
surgical cover screw (included)	 L-VT-340		



PRAMA RF POWER

ø 3.80 ●

ø 4.25 ●



LSS-ZT-380-085



LSS-ZT-425-085



LSS-ZT-380-100



LSS-ZT-425-100



LSS-ZT-380-115



LSS-ZT-425-115

-

-



L-VT-340

Regular Neck h 2.80 mm



PRAMA POWER				
h fixture	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●	
6.00 mm	 LA-ZT-380-060	 LA-ZT-425-060	 LA-ZT-500-060	
8.50 mm	 LA-ZT-380-085	 LA-ZT-425-085	 LA-ZT-500-085	
10.00 mm	 LA-ZT-380-100	 LA-ZT-425-100	 LA-ZT-500-100	
11.50 mm	 LA-ZT-380-115	 LA-ZT-425-115	 LA-ZT-500-115	
13.00 mm	 LA-ZT-380-130	 LA-ZT-425-130	 LA-ZT-500-130	
15.00 mm	 LA-ZT-380-150	 LA-ZT-425-150	-	
surgical cover screw (included)		 L-VT-340		



PRAMA RF POWER			PRAMA RF SL POWER		
ø 3.80 ●	ø 4.25 ●	ø 5.00 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
-	 LS-ZT-425-060	 LS-ZT-500-060	-	-	-
 LS-ZT-380-085	 LS-ZT-425-085	 LS-ZT-500-085	 LS-ZT-380SL-085	 LS-ZT-425SL-085	 LS-ZT-500SL-085
 LS-ZT-380-100	 LS-ZT-425-100	 LS-ZT-500-100	 LS-ZT-380SL-100	 LS-ZT-425SL-100	 LS-ZT-500SL-100
 LS-ZT-380-115	 LS-ZT-425-115	 LS-ZT-500-115	 LS-ZT-380SL-115	 LS-ZT-425SL-115	 LS-ZT-500SL-115
 LS-ZT-380-130	 LS-ZT-425-130	 LS-ZT-500-130	 LS-ZT-380SL-130	 LS-ZT-425SL-130	 LS-ZT-500SL-130
 LS-ZT-380-150	 LS-ZT-425-150	-	 LS-ZT-380SL-150	 LS-ZT-425SL-150	-



L-VT-340

Long Neck h 3.80 mm



PRAMA POWER			
h fixture	ø 3.80 ●	ø 4.25 ●	
8.50 mm	 LAL-ZT-380-085	 LAL-ZT-425-085	
10.00 mm	 LAL-ZT-380-100	 LAL-ZT-425-100	
11.50 mm	 LAL-ZT-380-115	 LAL-ZT-425-115	
13.00 mm	 LAL-ZT-380-130	 LAL-ZT-425-130	
15.00 mm	-	-	
surgical cover screw (included)	 L-VT-340		




PRAMA RF POWER		PRAMA RF SL POWER	
ø 3.80 ●	ø 4.25 ●	ø 3.80 ●	ø 4.25 ●
 LSL-ZT-380-085	 LSL-ZT-425-085	 LSL-ZT-380SL-085	 LSL-ZT-425SL-085
 LSL-ZT-380-100	 LSL-ZT-425-100	 LSL-ZT-380SL-100	 LSL-ZT-425SL-100
 LSL-ZT-380-115	 LSL-ZT-425-115	 LSL-ZT-380SL-115	 LSL-ZT-425SL-115
 LSL-ZT-380-130	 LSL-ZT-425-130	 LSL-ZT-380SL-130	 LSL-ZT-425SL-130
-	-	 LSL-ZT-380SL-150	 LSL-ZT-425SL-150
 L-VT-340			


PRAMA POWER: NASAL AND PTERYGOID

In cases of severe atrophy, there are alternatives to bone regeneration that can enhance treatment predictability and allow for immediate loading.

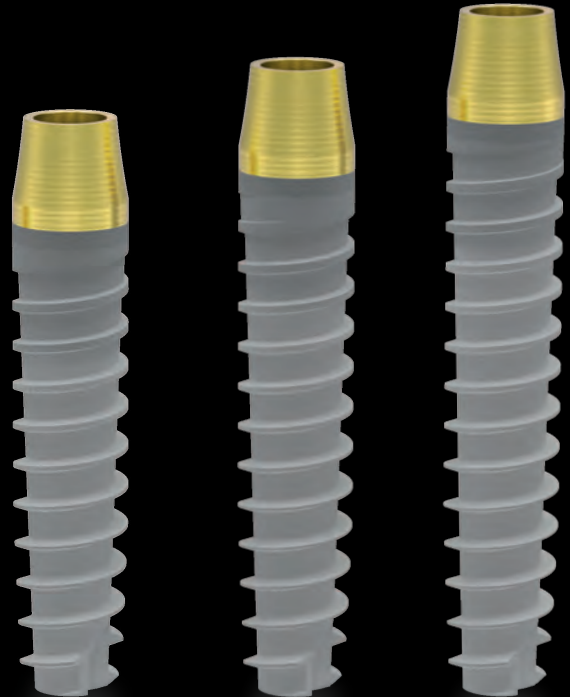The expansion of the Prama Power range enables the use of nasal and pterygoid positioning to increase primary stability and the biomechanical resistance of the entire rehabilitation.



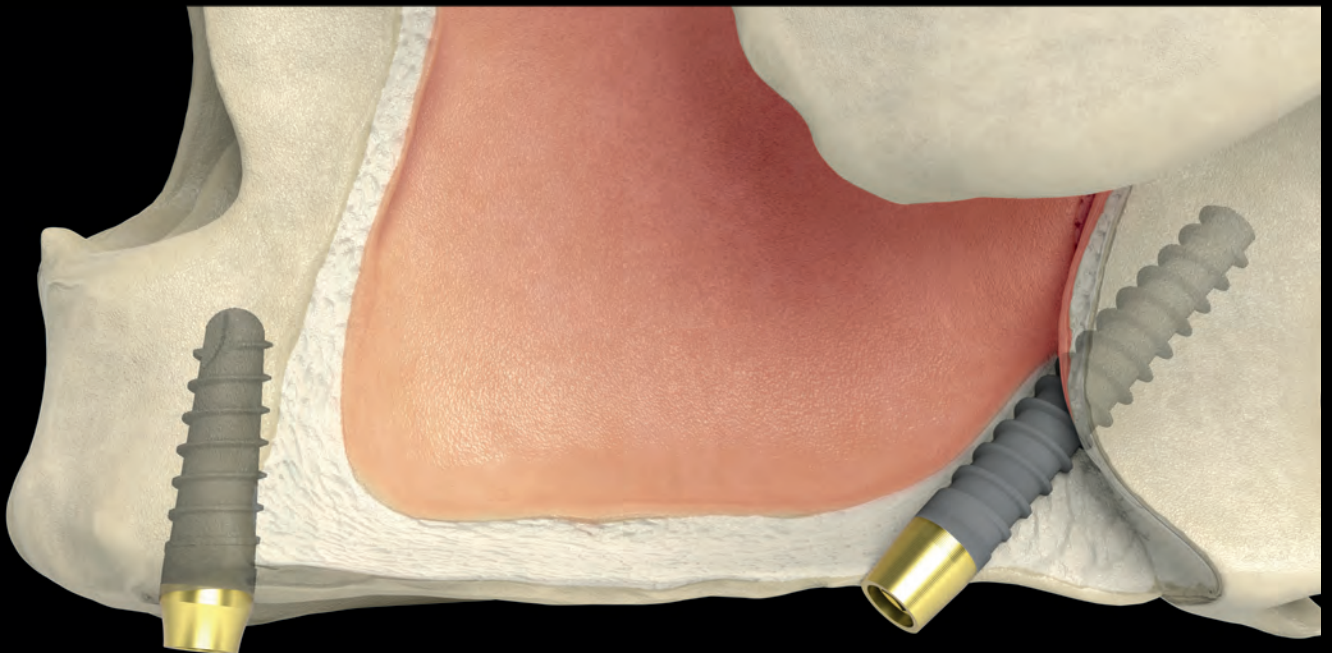These techniques, accessible to many intermediate to advanced implantologists, prove to be simple when combined with thorough pre-surgical planning that identifies the bone to engage, the implant inclination, and anatomical limitations.




Prama Power: Nasal and Pterygoid



SHORT NECK 1.80 mm		REGULAR NECK 2.80 mm		LONG NECK 3.80 mm		
h fixture	ø 3.80 ● ø 4.25 ●	ø 3.80 ● ø 4.25 ●	ø 3.80 ● ø 4.25 ●	ø 3.80 ● ø 4.25 ●	ø 3.80 ● ø 4.25 ●	
18.00 mm	 LSS-ZT-380SL-180	 LSS-ZT-425SL-180	 LS-ZT-380SL-180	 LS-ZT-425SL-180	 LSL-ZT-380SL-180	 LSL-ZT-425SL-180
20.00 mm	 LSS-ZT-380SL-200	 LSS-ZT-425SL-200	-	-	 LSL-ZT-380SL-200	 LSL-ZT-425SL-200
22.00 mm	 LSS-ZT-380SL-220	 LSS-ZT-425SL-220	-	-	 LSL-ZT-380SL-220	 LSL-ZT-425SL-220
surgical cover screw (included)	 L-VT-340					

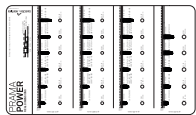Prama Power, Prama RF Power, Prama Power Nasal and Pterygoid X-ray templates



L-L100	real size
L-L120	20% larger size
L-L130	30% larger size


Drills: To insert Prama Power Nasal and Pterygoid implants, dedicated drills have been realized, with the relative stops, that allow to prepare all the new implant lengths.

Drills			Depth stops
	FP-200-L-SW	Pilot drill	STOP-2025-180-SW STOP-2025-190-SW STOP-2025-210-SW
	FK250-L-SW	Final conical drill	STOP-2025-220-SW STOP-2025-240-SW STOP-2025-260-SW
	FK380-L-SW	Final conical drill	STOP-3842-180-SW STOP-3842-190-SW STOP-3842-210-SW
	FK425-L-SW	Final conical drill	STOP-3842-220-SW STOP-3842-240-SW STOP-3842-260-SW


An optional tray is available to house the drills and stops for Prama Power Nasal and Pterygoid implants (CODE EXTR-26-TRAY-INT).
Not all products may be released for sale in all markets.

UTM and ZirTi: the synergy between two surfaces

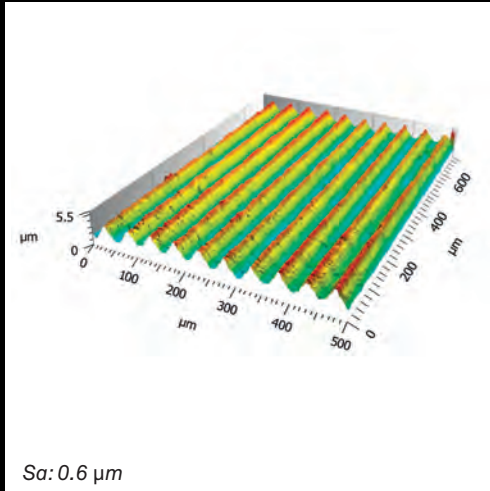UTM surface: the substrate designed for soft tissues



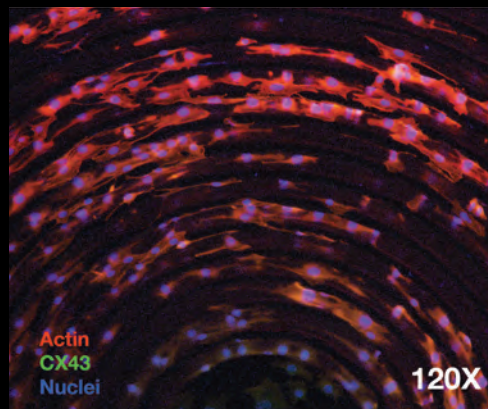It's not a machined surface, but a regular micro-threading that provides the correct sizing to enhance cellular organization.




Arrangement of fibroblasts on UTM surface after 72 hours, observed with fluorescence microscopy - in vitro experiment.

The micro-thread on the implant neck provides a guide for the unidirectional movement of cells, with **the biological advantage of faster, organized, and oriented deposition of collagen fibers.**

The resulting clinical benefit is a quicker healing process long-term maintenance of **healthy and stable tissues** over time.

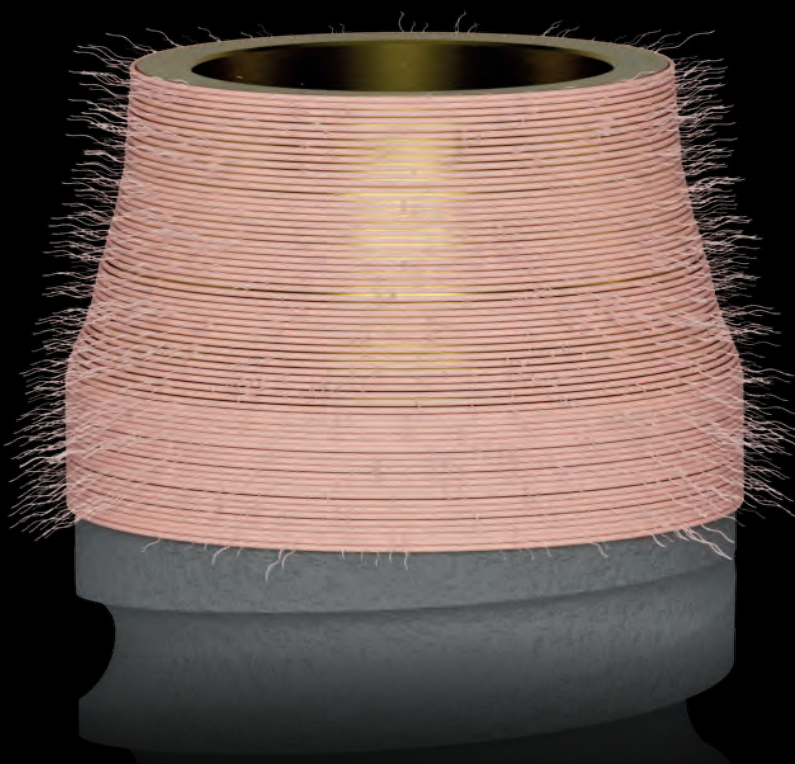Image courtesy of Prof. Macaluso, University of Parma, Department of Medicine and Surgery - Dentistry Center



Convergence and UTM

Greater fiber density

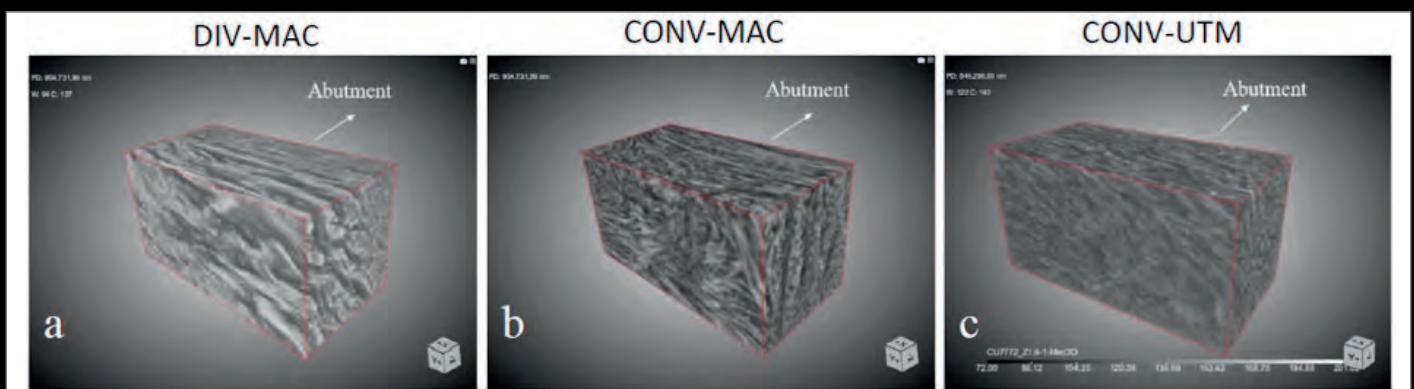
CLINICAL
ORAL IMPLANTS
RESEARCH

Influence of abutment macro- and micro-geometry on morphologic and morphometric features of peri-implant connective tissue

L. Canullo, A. Giuliani, M. Furlani, M. Menini, A. Piattelli, G. Iezzi

A revolutionary human histological study, published in Clinical Oral Implants Research, using **synchrotron tomography**, has demonstrated that the **synergy between the convergent profile of the transmucosal component and the UTM surface improves the adhesion and early activation of fibroblasts**, promoting the formation of a biological seal and enhancing the integration of peri-implant soft tissues. **Synchrotron tomography**, along with the significant difference in connective tissue height (vertical thickness) between convergent and divergent groups, also highlighted the significantly positive effect of the **3D micro-ridged UTM surface** (CONV-UTM) compared to the machined convergent surface (CONV-MAC) and even more so compared to the divergent surface (DIV-MAC) in terms of the density of interwoven collagen bundles.


Polarized light microscopy of the longitudinal section of the peri-abutment soft tissue around UTM convergent abutment. The light blue (red arrows) indicated the collagen fibers parallel to the vertical profile of the implant-abutment unit, whereas the dark blue color showed the collagen fibers that run perpendicularly or circularly to the implant-abutment unit.



The study demonstrated the **positive effect** of the **convergent** macrogeometry compared to traditional divergent one in terms of collagen fiber density. The convergent geometry of the implant-prosthetic portion in contact with the soft tissues plays a key role in the quality of peri-implant healing.

UTM and ZirTi: the synergy between two surfaces


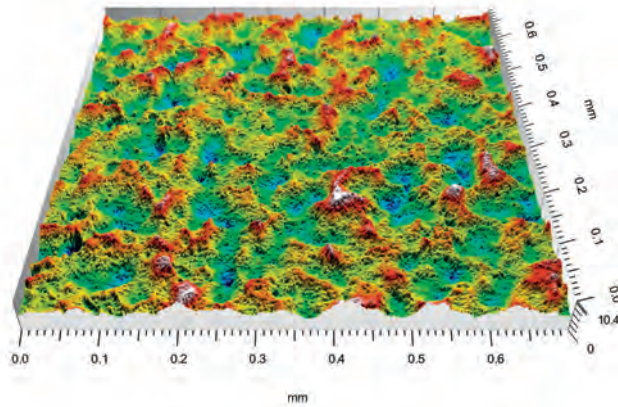
ZirTi surface: ideal for osseointegration



Sa: 1.4-1.7 μm




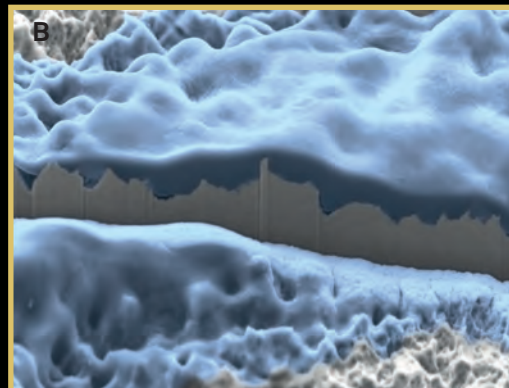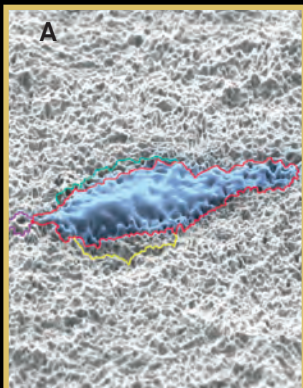The ZirTi surface is obtained with a sandblasting process with zirconium oxide, followed by an acid etching, techniques that contribute to a topography to which the osteoblasts have shown to intimately adapt and to activate rapid proliferation.



Enlargements of primary osteoblasts on the ZirTi surface: in images A and B, the **intimate adaptation of a cell on the peaks** that characterize the ZirTi can be observed, thanks to a section obtained with FIB (Focused Ion Beam).

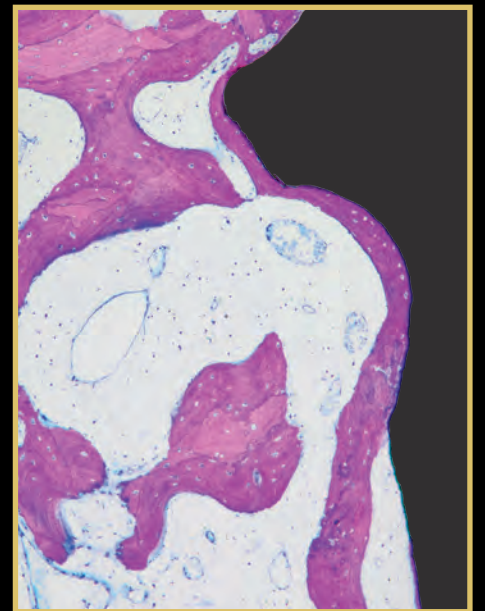Image courtesy of Prof. Macaluso, University of Parma, Department of Medicine and Surgery - Dentistry Center




Histology in humans at 3 months: even in bone D4, excellent mineralization is observed in contact with the ZirTi surface.

Image courtesy of prof. Marco Gargari and dr. Antonio Rocci.




**Bibliography on
ZirTi surface.**

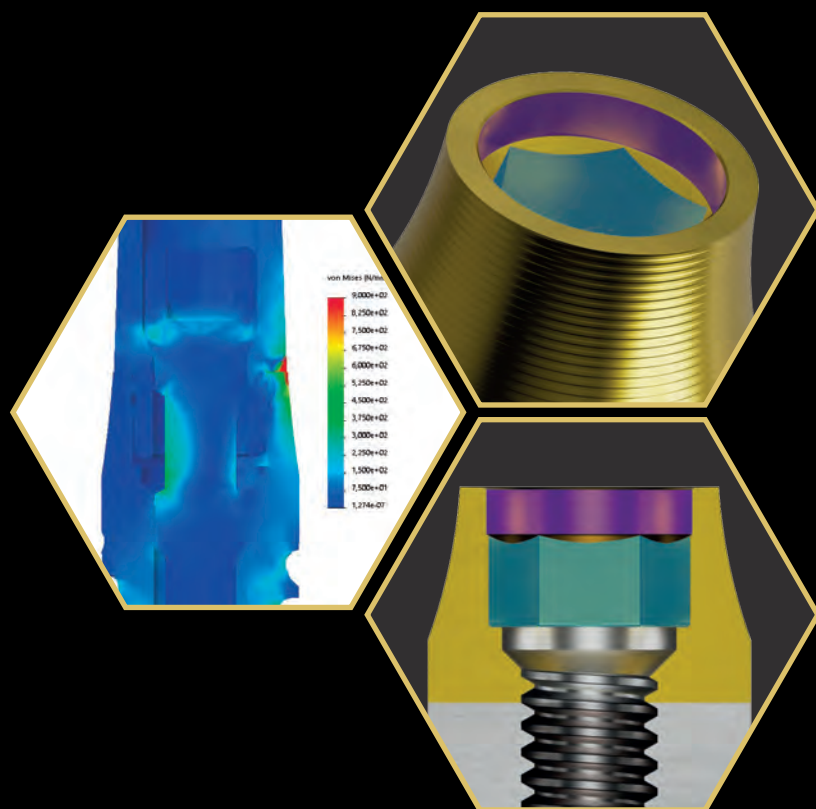The advantages of COLLEX CONNECTION



COLLEX CONNECTION:

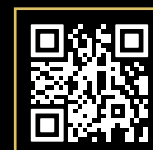
several studies have documented its robustness and prosthetic stability compared to other connections without a collar.



Bibliography on
Collex Connection

INTERNAL HEXAGON connection with a collar of the same size for all diameters.

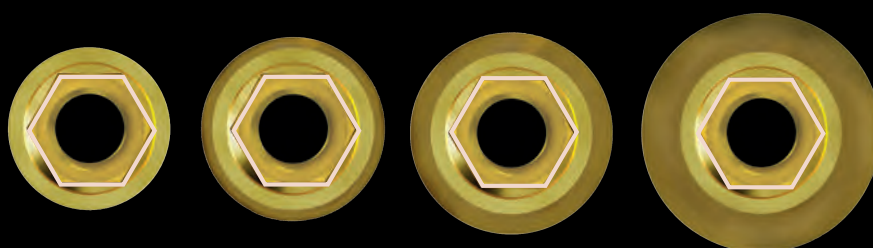
FOUR DIFFERENT IMPLANT DIAMETERS, only one prosthetic connection



Horizontal and vertical platform switching TISSUES STABILITY

For many years now, scientific literature has been focused on the concept of **platform switching**, a prosthetic rehabilitation technique that involves using abutments with a smaller diameter than the implant platform. This technique aims to improve the biomechanical distribution of prosthetic load and, more importantly, to distance the prosthetic connection from the cervical bone, thus moving the critical point of bacterial infiltration away from the bone.

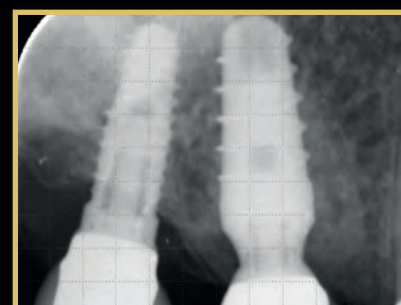It is from this experience that the idea of **maximizing the platform switching with the Prama Power neck** arises, utilizing both its horizontal and vertical components.



Courtesy of Dentisti Vignato

Guidelines for the insertion of PRAMA with **Intramucosal protocol**

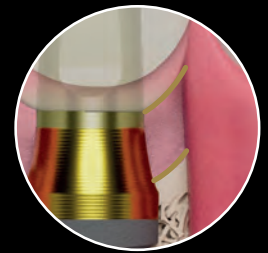
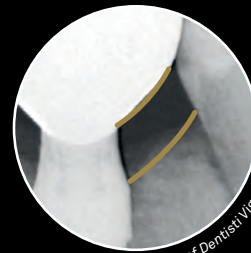
PRE-OP SITUATION

The insertion of Prama implant should be planned starting from the prosthesis.

The desired profile of the parabola is to be taken as a reference, which must be harmonious with the festooning of the adjacent elements or of the contralateral element.

FINAL SITUATION

Literature and clinical experience show that after the healing a bone remodelling occurs, that is parallel to the emergence profile of the prosthetic crown, approximately 3 mm from the latter. Therefore, the positioning of the implant should be such as to maximize the biological space available for the connective tissue. The use of a vertical post will allow the crown to close at the desired level, with its zenith positioned at least 1 mm below the gum margin. Bone remodeling and regrowth of the soft tissues will take place in an optimal manner and maximizing the aesthetic and functional results of the rehabilitation.



Courtesy of Dentisti Vignaro



SURGICAL PHASE

In light of what we have just seen, the implant should be placed sufficiently deep. Therefore, the advice is to choose the implant at least one length shorter than the depth of the prepared site. For example, if the available bone allows the insertion of a 10 mm implant, a 10 mm long hole will be prepared, but a 8.5 mm implant will be placed, and the insertion depth will be modulated up to the level in such a way as to take advantage of the intramucosal potential of the Prama neck.

Do you know why the synergy between UTM micromorphology e convergent macromorphology of the Prama neck has a significant impact on the quality of peri-implant soft tissues?

Do you know why the synergy between UTM micromorphology and the convergent macromorphology of the Prama neck has such a significant impact on the quality of peri-implant soft tissues?

Read the complete article published on the Clinical Oral Implants Research



What if the patient has a
VERY THIN GINGIVAL BIOTYPE?



PRE-OP SITUATION

In case of particularly thin gingival biotype (approximately 1 mm of gum), the protocol for intramucosal insertion and for subsequent rehabilitation of a Prama implant does not change. Just further deepen the level of insertion of the implant into the bone, procedure possible thanks to the UTM surface of the Prama neck that works in a excellent way also in contact with hard tissues.



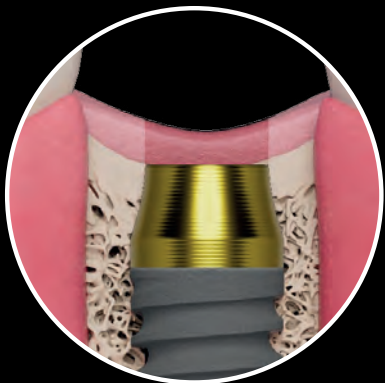
FINAL SITUATION

Planning always should start from the prosthesis: once designed the position of the prosthetic crown, whose zenith will always go assumed approximately 1 mm below the gingival margin, the insertion of the implant follows the protocol described in the previous pages.

SURGICAL PHASE



In case of thin biotype, since a deeper placement is needed, it is recommended to choose an shorter implant length. In the example on the left, the drill is used at the depth of 11.5 mm for the insertion of a 8.5 mm long implant.



PRAMA POWER

THE IMPLANT FOR A COMPLETE DIGITAL WORKFLOW

Digital Atelier

A team for all services of:

- Design;
- Custom made prosthesis production;
- Guided surgery;
- 3D printing;
- Assistance.



Courtesy of dr. Davide Di Paola

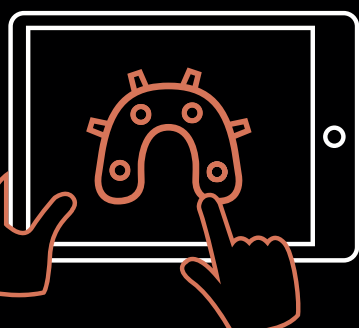


DATA ACQUISITION



3D PRINTING
AND MILLING





GUIDED SURGERY



ECHOPLAN PRO A

ECHOPLAN PRO S

Guided Implantology

With the use of the Archiplan software it is possible to create surgical templates. Guided surgery is made particularly simple thanks to the availability of specific dedicated instruments: surgical kits complete with drills, mucotomes, levelers, drivers, screwdrivers and everything needed to perform guided surgery with Prama Power. The drills with integrated stop allow you to save time and obtain predictable results. For details on the contents of the guided surgery kits with which to insert Prama Power, see page 26.



CAD-CAM

Custom Made Prosthesis

The Digital Atelier can:

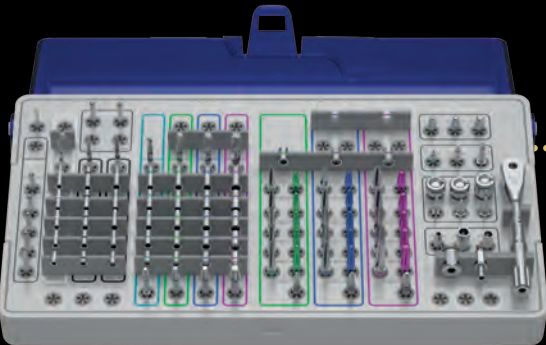
- scan physical models received from customers;
- acquire virtual models already scanned intraorally or in the lab by customers and sent in the form of files, coming from all open scanner systems;
- virtually model any prosthetic product on customer prescription;
- produce custom-made products with CAD CAM technology.



Download the Exocad guides here

Ergonomics of Surgical Procedures

The **SURGICAL KITS** are comprehensive, **compact, ergonomic**, and contain all the necessary tools for the surgery of **Prama Power implants**.



The color codes, also present on the surgical trays, guide the correct sequence of instrument usage.



Both the drills and the stops are characterized by a color code that guides operators in tool selection.



The **SURGICAL TRAYS** are made of **Radel**, a polymer that maintains its characteristics even after numerous cycles of washing and autoclave sterilization.

CYLINDRICAL or CONICAL DRILLS, based on different preparations, are supplied complete with depth stops.

The drills are made of **stainless steel** for surgical use and maintain their characteristics unchanged over time and through multiple uses.



The drills are characterized by the presence of clearly visible depth notches.

The drills are **highly sharp, centered, and perform well**. The morphology of the conical drills allows for the collection of valuable bone for regeneration needs.





The kits include a convenient **RATCHET**, which can be used both as a fixed wrench and as a dynamometric wrench for controlling the implant insertion torque and tightening of prosthetic screws.

- The torque can be adjusted from 10 to 70 Ncm;
- The torque adjustment is highly accurate because each individual ratchet is individually calibrated.

Included is a wide range of **PROSTHETIC DRIVERS** for digital use, for the dynamometric ratchet and for micromotor.

The range is suitable for different oral access needs, from extra-short tools to longer instruments for Toronto-type prosthetic rehabilitations.

The slightly tapered coupling between the driver and the screw allows for adequate retention, ensuring secure transport of the screw in the oral cavity.



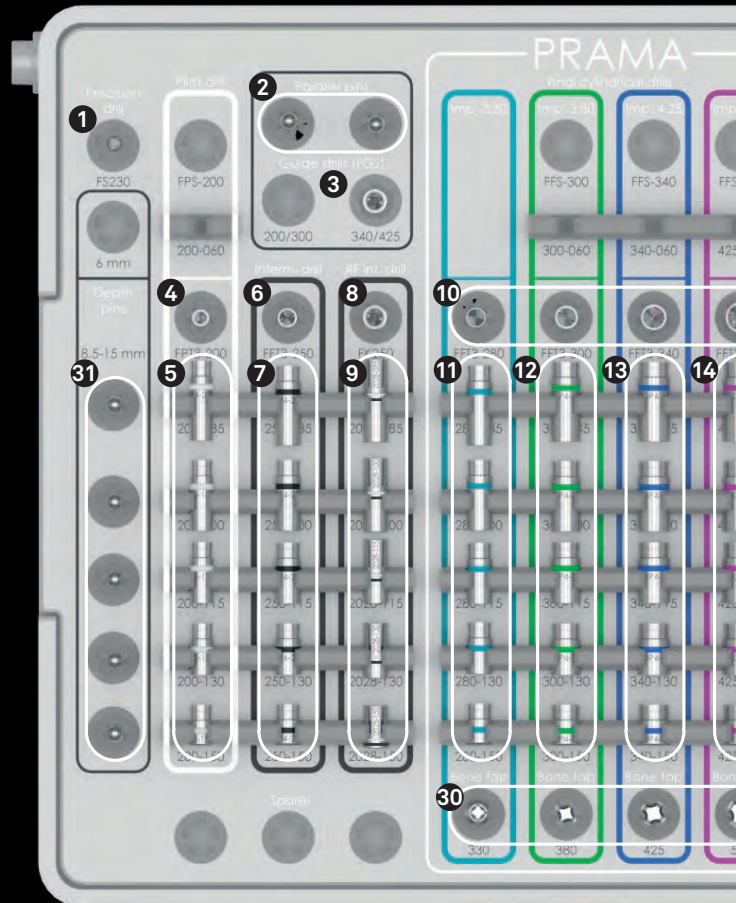
The surgical procedure involves a countless approach that ensures visibility and ergonomics with the use of the patented **EASY INSERT** drivers:

- it can be easily removed from implant wells even in the presence of high insertion torque;
- the dodecagonal design prevents deformation of the connection edges, ensuring the maintenance of precision and prosthetic stability;
- available in the kit in various configurations, both for use with a micromotor and the ratchet;
- the instruments have a hexagon on the shaft, aligned with the connection hexagon, to facilitate rotational orientation of the fixture.

Prama Power surgical kit


Complete surgical kit containing all the necessary tools for the insertion of Prama Power implants (code ZPRAMA-INT).
All instruments can be purchased individually as spare parts.

- 1  **Precision drill**
FS-230
- 2  **Parallelism pins**
PPTL-2-28
- 3  **Intermediate drill**
• FGS-340/425
- 4  **Pilot drill***
• FPT3-200-LXS
- 5  **Stops for pilot drills**
• STOP4-200-085
• STOP4-200-100
• STOP4-200-115
• STOP4-200-130
• STOP4-200-150
- 6  **Intermediate cylindrical drill**
ø 2.50 mm*
• FFT3-250-LXS
- 7  **Stop for intermediate cylindrical drill**
ø 2.50 mm
• STOP4-250-085
• STOP4-250-100
• STOP4-250-115
• STOP4-250-130
• STOP4-250-150
- 8  **Intermediate conical drill**
• SE-FK250
- 9  **Stops for Intermediate conical drill**
CSR-STOP-2028-085
CSR-STOP-2028-100
CSR-STOP-2028-115
CSR-STOP-2028-130
CSR-STOP-2028-150
- 10  **Final cylindrical drill***
• FFT3-280-LXS for implants ø 3.30 mm
• FFT3-300-LXS for implants ø 3.80 mm
• FFT3-340-LXS for implants ø 4.25 mm
• FFT3-425-LXS for implants ø 5.00 mm
- 11  **Stops for conical drills**
FFT3-280-LXS
• STOP4-280-085
• STOP4-280-100
• STOP4-280-115
• STOP4-280-130
• STOP4-280-150
- 12  **Stops for cylindrical drills**
FFT3-300-LXS
• STOP4-300-085
• STOP4-300-100
• STOP4-300-115
• STOP4-300-130
• STOP4-300-150
- 13  **Stops for cylindrical drills**
FFT3-340-LXS
• STOP4-340-085
• STOP4-340-100
• STOP4-340-115
• STOP4-340-130
• STOP4-340-150
- 14  **Stops for cylindrical drills**
FFT3-425-LXS
• STOP4-425-085
• STOP4-425-100
• STOP4-425-115
• STOP4-425-130
• STOP4-425-150
- 15  **Conical drills**
• SH-FK380-085
• SH-FK380-100
• SH-FK380-115
• SH-FK380-130
• SH-FK380-150
- 16  **Replies**
• SH-380-085-RP
• SH-380-100-RP
• SH-380-115-RP
• SH-380-130-RP
• SH-380-150-RP
- 17  **Conical drills**
• SH-FK425-085
• SH-FK425-100
• SH-FK425-115
• SH-FK425-130
• SH-FK425-150
- 18  **Replies**
• SH-425-085-RP
• SH-425-100-RP
• SH-425-115-RP
• SH-425-130-RP
• SH-425-150-RP
- 19  **Conical drills**
• SH-FK500-085
• SH-FK500-100
• SH-FK500-115
• SH-FK500-130
• SH-FK500-150
- 20  **Replies**
• SH-500-085-RP
• SH-500-100-RP
• SH-500-115-RP
• SH-500-130-RP
• SH-500-150-RP



L-TRAYL-INT instrument carrying case in Radel.

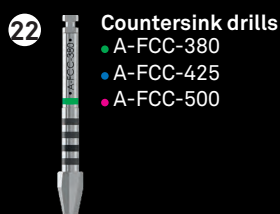- 21** Stops for conical drills
- SH-STOP4-FK380
 - SH-STOP4-FK425
 - SH-STOP4-FK500



- 23** Easy Insert Drivers for implant insertion
- EASYC4-EX230-CA



EASYL4-EX230-CA




EASY4-EX230-EX



- 25** Torque control ratchet
- CRI5-KIT



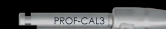- 26** Extensions and adapter
- PROF-CAL3
- Extension for surgical drills




- BPM-15
- Extensions for instruments with hexagonal connector



- B-AVV-CA3
- Mechanical adapter for instruments with hexagonal connector



- 27** Manual adapter for instruments with hexagonal connector
- AVV3-MAN-DG



- 28** Adapter for instruments with hexagonal connector
- AVV-CA-DG-EX



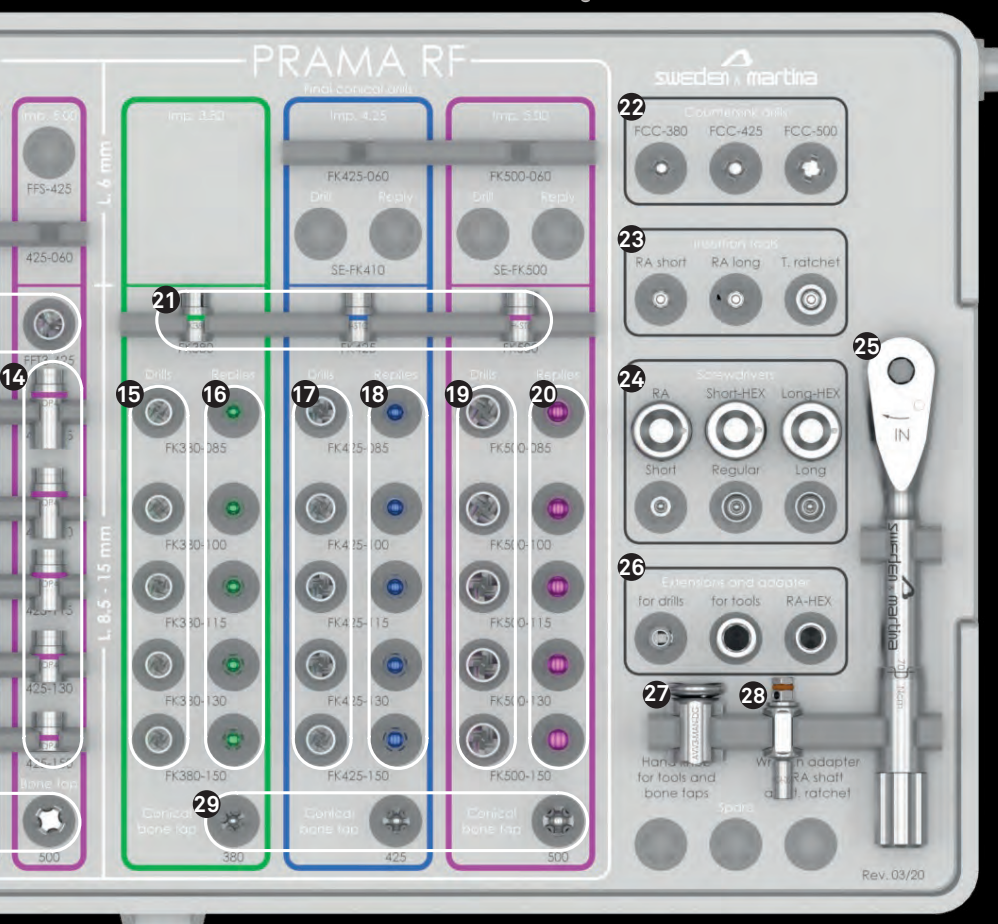
- 31** Pins for neck positioning
- L-PP-085
 - L-PP-100
 - L-PP-115
 - L-PP-130
 - L-PP-150



- 24** Digital screwdriver
- HSMXS-20-DG HSM-20-DG HSML-20-DG



- Prosthetic screwdriver
- HSM-20-CA HSM-20-EX HSML-20-EX



- 29** Conical bone taps
- SH-MS-380-CA
 - SH-MS-425-CA
 - SH-MS-500-CA

- 30** Cylindrical bone taps
- SM-MS-330
 - A-MS-380
 - A-MS-425
 - A-MS-500

X-ray templates included in the kit, they can be reordered separately with the codes L-L100 (real size), L-L120 (20% larger size), L-L130 (30% larger size).

1 Mucotomes
GS-MUC-415
GS-MUC-550

2 Bone profilers
GS-LC-415
GS-LC-550

3 Initial drills and progressive drill
GS-F200-415
GS-F200-550
GS-F2030-550

4 Drivers
HSM-20-EX HSML-20-EX HSM-20-DG HSML-20-DG HSM-20-CA

Optional drills for the preparation of

18

**Drills for the preparation of
h 18.00 mm implants
ø 4.25 mm**

GS-F200-180-550

GS-F340-180-550

17

**Drills for the preparation of
h 18.00 mm implants
ø 3.80 mm**

GS-F200-180-415
GS-F300-180-415

16
Pins and drill for pins
GS-PIN (4 pcs) GS-FPN-148

15 Adapter for instruments with hexagonal connector
AVV3-MAN-DG

14

**Torque control
ratchet
CR15-KIT**

13

Drivers for mounters
EASY4-EX250-EX



EASYC4-EX250-CA

12 **Countersink drills**
 GS-FCS-A380 for $\varnothing$ 3.80 mm implant
 GS-FCS-A425 for $\varnothing$ 4.25 mm implant
 GS-FCS-A500 for $\varnothing$ 5.00 mm implant

Optional drills for preparing 6.00 mm high implants can be ordered also cumulatively as a supplementary set with the code GS-PROA-INTEGRA-060.

5



**Optional drills for h 7.00 mm
ø 3.80 mm implants**
GS-F200-060-415
GS-F250-060-415
GS-F280-060-415
GS-F300-060-415

6



**Optional drills for h 7.00 mm
ø 4.25 and 5.00 mm implants**
GS-F200-060-550
GS-F300-060-550
GS-F340-060-550
GS-F425-060-550

7



**Intermediate drills
for ø 4.15 mm sleeve**
GS-F200-085-415
GS-F200-100-415
GS-F200-115-415
GS-F200-130-415
GS-F200-150-415

GS-F250-085-415
GS-F250-100-415
GS-F250-115-415
GS-F250-130-415
GS-F250-150-415

8



**Final drills for
ø 4.15 mm sleeve**
GS-F280-085-415
GS-F280-100-415
GS-F280-115-415
GS-F280-130-415
GS-F280-150-415

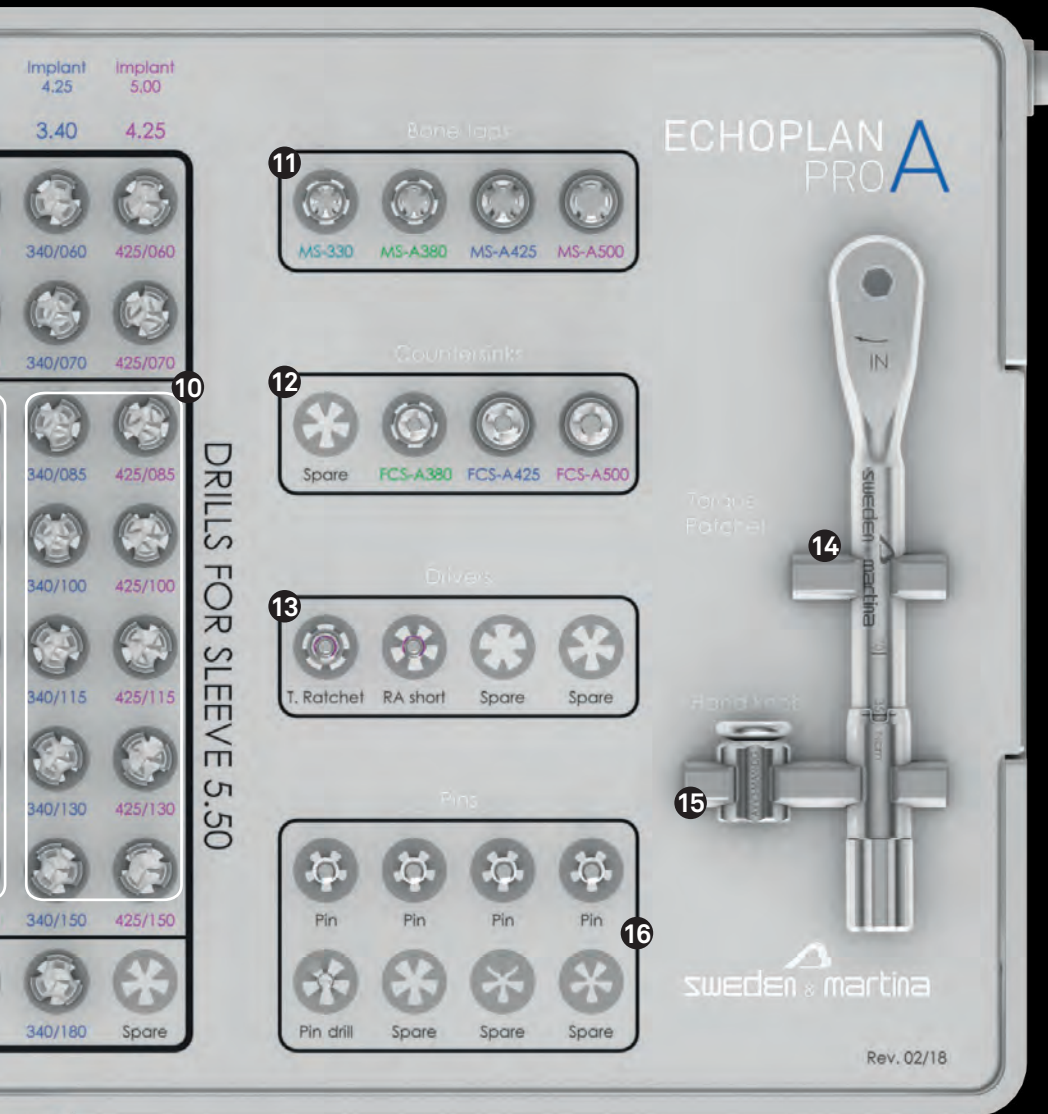
GS-F300-085-415
GS-F300-100-415
GS-F300-115-415
GS-F300-130-415
GS-F300-150-415

9



**Intermediate drills for
ø 5.50 mm sleeve**
GS-F200-085-550
GS-F200-100-550
GS-F200-115-550
GS-F200-130-550
GS-F200-150-550

GS-F300-085-550
GS-F300-100-550
GS-F300-115-550
GS-F300-130-550
GS-F300-150-550



11



Bone taps
GS-MS-330 per impianto ø 3.30 mm
GS-MS-A380 per impianto ø 3.80 mm
GS-MS-A425 per impianto ø 4.25 mm
GS-MS-A500 per impianto ø 5.00 mm

10



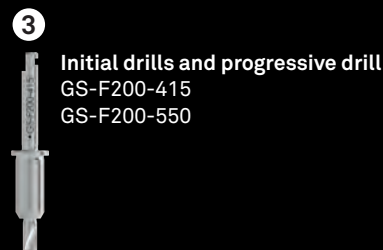
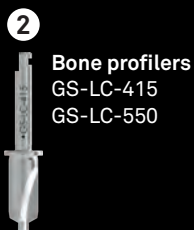
**Final drills for ø 4.25 mm implants
with ø 5.50 mm sleeve**
GS-F340-085-550
GS-F340-100-550
GS-F340-115-550
GS-F340-130-550
GS-F340-150-550

**Final drills for ø 5.00 mm implants
with ø 5.50 mm sleeve**
GS-F425-085-550
GS-F425-100-550
GS-F425-115-550
GS-F425-130-550
GS-F425-150-550

1 Mucotomes
GS-MUC-415
GS-MUC-550

2 Bone profilers
GS-LC-415
GS-LC-550

3 Initial drills and progressive drill
GS-F200-415
GS-F200-550



18

GS-F200-060-415
GS-FK410U-060-415
GS-FK410-060-415

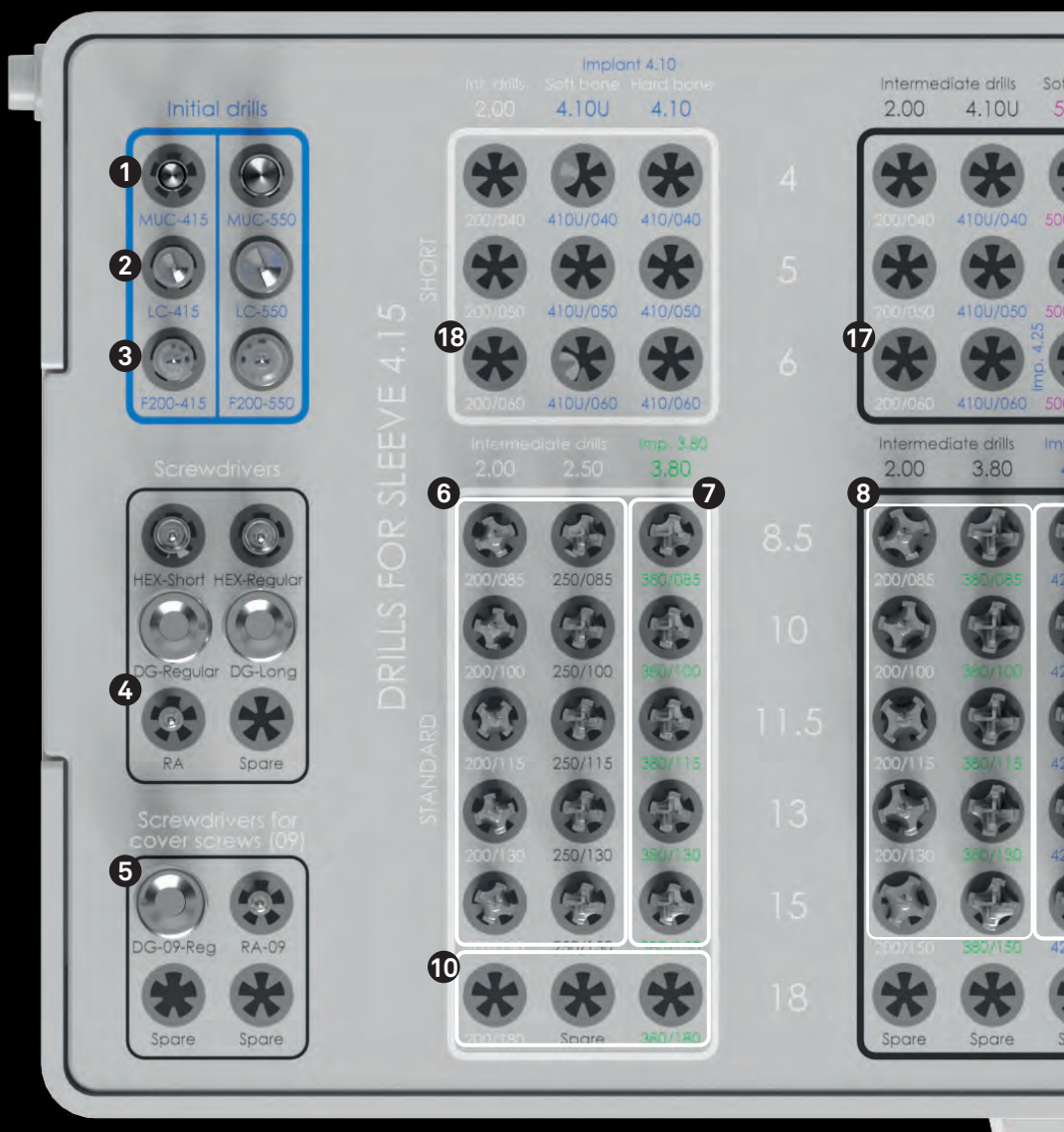
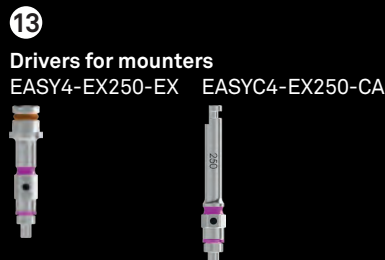


17

- GS-F200-060-550
- GS-FK410U-060-550
- GS-FK500U-060-550
- GS-FK500-060-550



16 Adapter for instruments with hexagonal connector
AVV3-MAN-DG



Guided surgery accessories

implant ø		3.30 ●		3.80 ● 4.25 ● - 5.00 ●
		 4.00  GS-B415-EX-6* GS-B415-SET-100*	 4.00  GS-B415-EX-6* GS-B415-SET-100*	 4.00  GS-B550-EX-6* GS-B550-SET-100*
mounter	SHORT NECK	-	 7.20 GS-MOU2-L415-072	 7.20 GS-MOU2-L500-072
	REGULAR NECK	 9.00 GS-MOU-A330	 6.20 GS-MOU-L415	 6.20  9.00 GS-MOU-L550 GS-MOU-A380SP**
	LONG NECK	-	 5.20 GS-MOU2-L415-052	 5.20 GS-MOU2-L500-052
fixation screw included		 GS-VTMOU-180**	 GS-VTMOU-180**	 GS-VTMOU-180**
easy mounter		 9.00 GS-EASY-AS-415-EX	 9.00 GS-EASY-L-415-EX	 9.00 GS-EASY-L-550-EX
		 9.00 GS-EASY-AS-415-CA	 9.00 GS-EASY-L-415-CA	 9.00 GS-EASY-L-550-CA
easy mounter with depth indication		-	 9.00 GS-EASYPE-L-415-EX	 9.00 GS-EASYPE-L-550-EX
		-	 9.00 GS-EASYPE-L-415-CA	 9.00 GS-EASYPE-L-550-CA

*The engaging guide sleeves are available in packs of 6 pieces, and in sets of 100 pieces, not sold individually.
 **It is recommended, with this mounter, to adjust the implant site preparation to ensure the torque value does not exceed 50 Ncm.
 ***The fixation screw is included with the mounter but can also be ordered individually as a spare.

Bone taps without stop

implant ø	3.30 ●	3.80 ●	4.25 ●	5.00 ●
Prama Power				
	GS-MSPE-A330	GS-MSPE-A380	GS-MSPE-A425	GS-MSPE-A500
Prama RF Power, Prama RF SL Power				
		GS-MSPE-S380	GS-MSPE-S425	GS-MSPE-S500

Countersink drills without stop

implant ø	3.30 ●	3.80 ●	4.25 ●	5.00 ●
Prama Power				
	-	GS-FCSPE-380	GS-FCSPE-425	GS-FCSPE-500

Sleeves for pins

6 item pack



Guidelines for Prama Power
insertion with guided surgery



Optional integration sets

L-INTEGRA-060 set


The L-INTEGRA-060 drill set includes the components necessary for the insertion of the Prama Power and Prama RF Power h 6.00 mm implants, the housings of which are provided inside the Prama Power kit.

Initial drills Prama Power h 6.00 mm*	
pilot drill	Shorty intermediate drill
	
FPS-200	FGS-200/300
	-
STOPS4-200-060	-

* Instruments included in ZSHORTY-INT

Final drills Prama Power h 6.00 mm*		
ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
		
FFS-300	FFS-340	FFS-425
		
STOPS4-300-060	STOPS4-340-060	STOPS4-425-060

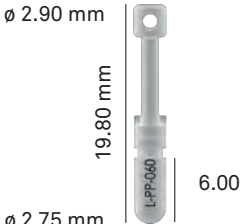* Instruments included in ZSHORTY-INT

Parallelism posts

L-PP-060

Final drills Prama RF Power h 6.00 mm		
ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
-		
-	SE-FK410	SE-FK500
-		
-	SE-STOP-FK410-060	SE-STOP-FK500-060

Reply for Prama RF Power h 6.00 mm	
	
SE-410-060-RP	SE-500-060-RP

Set L-INTEGRA set

L-INTEGRA set includes all the necessary components for the insertion of Prama Slim Power ø 3.30 mm, the countersink drills and the posts for neck positioning. These components are included in the ultimate version of the Prama Power surgical kit, but not in the previous version of the surgical kit, thus they are still available for those who have the previous version.



Posts for neck positioning		
	h	
	8.5	L-PP-085
	10.00	L-PP-100
	11.50	L-PP-115
	13.00	L-PP-130
	15.00	L-PP-150

Instruments for Prama Slim Power implant			
final drill	stop		
		h	
		8.50	STOP4-280-085
		10.00	STOP4-280-100
		11.50	STOP4-280-115
		13.00	STOP4-280-130
		15.00	STOP4-280-150
FFT3-280-LXS			


Optional surgical instruments (not included in Prama Power surgical kit)

Long precision drill	Intermediate drills	
		
FSL-230	FG-200/280XS	FG-330/425XS
	intermediate drill 200/280	intermediate drill 330/425

Drills for distal sectors				
ø 2.00	ø 2.80	ø 2.90	ø 3.00	ø 3.20
				
FPT5-200-LXS	FFT5-280-LXS	FFT5-290-LXS	FFT5-300-LXS	FFT5-320-LXS
ø 3.30	ø 3.40	ø 3.60	ø 4.25	ø 4.45
				
FFT5-330-LXS	FFT5-340-LXS	FFT5-360-LXS	FFT5-425-LXS	FFT5-445-LXS

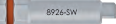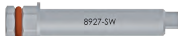
Screwdrivers not included in the kit		
	HSMXL-20-EX	with hexagonal connector, extra long
	HSMXS-20-CA	with right angle shank, extra short
	L-HSM-EX	Full Head, with hexagonal connector, short
	L-HSML-EX	Full Head, with hexagonal connector, long
	L-HSMXL-EX	Full Head, with hexagonal connector, extra long
	L-HSM-CA	Full Head, with right angle shank
	AVV-CAP-TIT-1	instrument for maintaining the titanium cap for ball attachments
	BASCC-EX	for ball attachments, with hexagonal connector
	AVV2-ABUT	for straight P.A.D. abutments
	8926-SW	for Locator abutment, short
	8927-SW	for Locator abutment, long
	DSPDCLH-24	for angled screw hole, length 24 mm
	DSPDCLH-32	for angled screw hole, length 32 mm

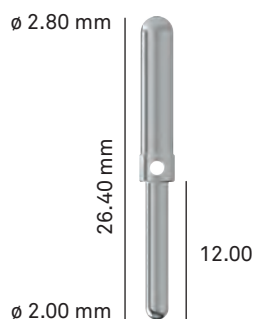Parallel pins




PP-2/28

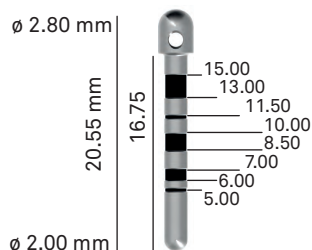Parallelism pin with one side of $\varnothing$ 2.00 mm and the other side of $\varnothing$ 2.80 mm



PPTS-2-28

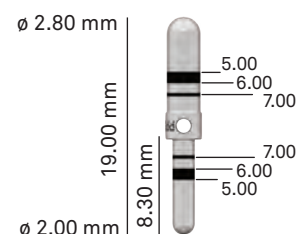Parallelism pin with depth notches, small version



PPS-2/3

Parallelism pin for short drills with notches at 5.00, 6.00, and 7.00 mm

Bone taps

$\varnothing$ 3.30

$\varnothing$ 3.80

$\varnothing$ 4.25

$\varnothing$ 5.00

Bone taps with hexagonal connector



SM-MSC-330

A-MSC-380

A-MSC-425

A-MSC-500



-

SH-MS-380

SH-MS-425

SH-MS-500

Bone taps with right angle shank



SM-MS-330-CA

A-MS-380-CA

A-MS-425-CA

A-MS-500-CA

Extraction tools for intraoral removing of implants



BC-EX230

short



BL-EX230

long

Extractors for atraumatic removal of implants



EXTR-IMP-S

small



EXTR-IMP-M

medium



EXTR-IMP-L

large

Bone profiler for P.A.D. abutments

unique platform



A-PAD-PS330-L

transgingival height 3.00, wide

A-PAD-PS330-S

transgingival height 5.00, narrow



A-PAD-GUI-PS-230

guide for profiler

Osteotomes

				
— $\varnothing$ 0.20	— $\varnothing$ 0.90	— $\varnothing$ 1.60	— $\varnothing$ 2.00	— $\varnothing$ 2.40
E-OS-020-PP	E-OS-090-PP	E-OS-160-PC	E-OS-200-PC	E-OS-240-PC
osteotome $\varnothing$ 0.20 flat tip	osteotome $\varnothing$ 0.90 flat tip	osteotome $\varnothing$ 1.60 concave tip	osteotome $\varnothing$ 2.00 concave tip	osteotome $\varnothing$ 2.40 concave tip

Osteotomes with flat tip

	h	$\varnothing$ 3.80 ●	$\varnothing$ 4.25 ●	$\varnothing$ 5.00 ●
	8.50-10.00	SH-OS-380-100-PP	SH-OS-425-100-PP	SH-OS-500-100-PP
	11.50	SH-OS-380-115-PP	SH-OS-425-115-PP	SH-OS-500-115-PP
	13.00	SH-OS-380-130-PP	SH-OS-425-130-PP	SH-OS-500-130-PP
	15.00	SH-OS-380-150-PP	SH-OS-425-150-PP	SH-OS-500-150-PP

Osteotomes with round tip

	h	$\varnothing$ 3.80 ●	$\varnothing$ 4.25 ●	$\varnothing$ 5.00 ●
	8.50-10.00	SH-OS-380-100-PR	SH-OS-425-100-PR	SH-OS-500-100-PR
	11.50	SH-OS-380-115-PR	SH-OS-425-115-PR	SH-OS-500-115-PR
	13.00	SH-OS-380-130-PR	SH-OS-425-130-PR	SH-OS-500-130-PR
	15.00	SH-OS-380-150-PR	SH-OS-425-150-PR	SH-OS-500-150-PR

Torque wrench with control lever

	TWL	
	TWL-AVV-EX	Adaptor for torque wrench

Depth gauge

	PROF3
--	--------------

Spare O-ring

	ORING180-088
---	---------------------

Service mounter for impact procedures

	MOU-EX230	mounter
	CM2	mounter stop key


Soft tissue depth gauge

	MP-LMRC-180	h 10 mm
--	--------------------	---------


Depth notches

The laser marking on the cutting edges and the use of stops allow for safe preparation and control of the instrument's insertion depth. The drills always create a hole longer than the desired implant. The oversizing (LS) is equal to the height of the tip of the drill being used.

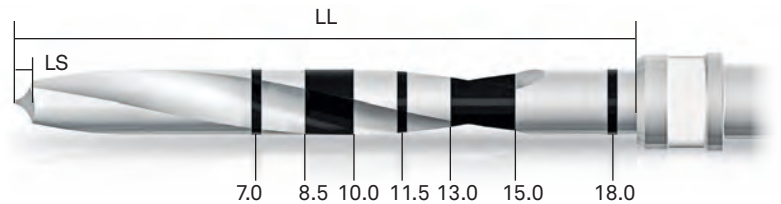LL: total length of the working part, including the tip.
LS: length of the overpreparation. This measurement must be calculated in addition to the length of the preparation hole.



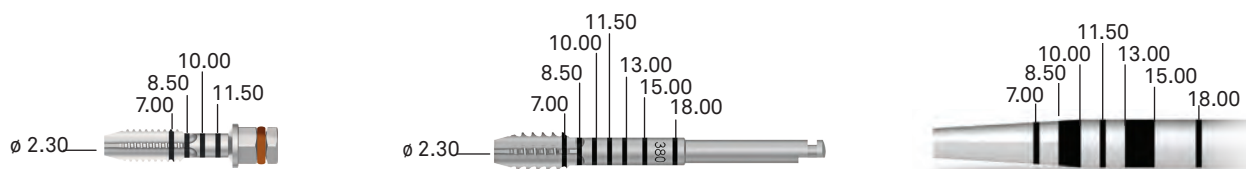The working depth is also identified on bone taps and osteotomes through laser-etched notches, as follows:




The working depth is also identified on bone taps without a stop and Countersink drills through laser markings, as follows:



Prosthetic solutions of excellence

Marginless solutions

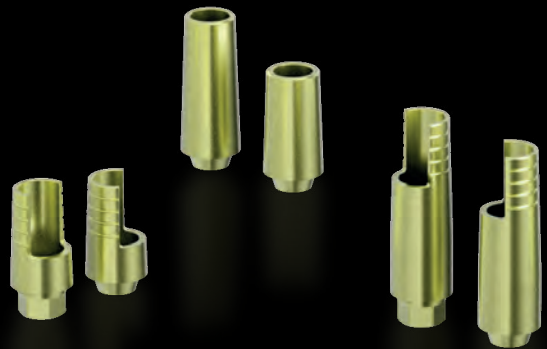Marginless and with margin options that leave to the clinician the decision on where to close the prosthetic crown margin (L-MD or L-MDT). These solutions promote soft tissues thickening because they leave more space for connective tissue.



Prama IN solutions

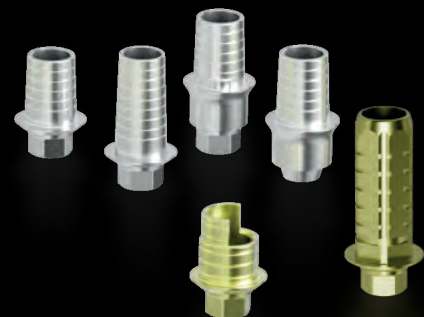
When single rehabilitation concerns a last molar element, an option is to use a “Prama IN” prosthetic solution which closes around the implant neck, allowing greater robustness to the rehabilitation, also in light of the important typical chewing loads of that area.

The “Prama IN” prosthetic solutions may also represent a valid option in cases where, due to factors related to the patient’s immune system response, bone and tissue remodeling was not the desired one.




T-CONNECT

Standard or cuttable T-Connects for straight or angled screw hole.



CONICO system

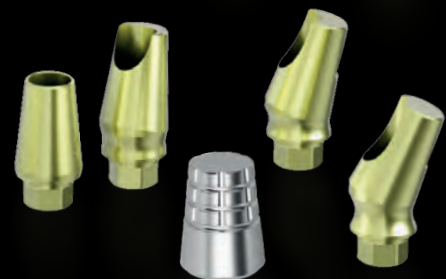- Conometric prosthesis without screws or cement;
- easy to remove: at any time the clinician can remove the prosthesis to perform a follow up and/or the hygiene;
- the conometric technique allows restorations from single elements to entire arches;
- biological seal and healthy and stable tissues over time.



PAD^e abutments for multiple prosthetic

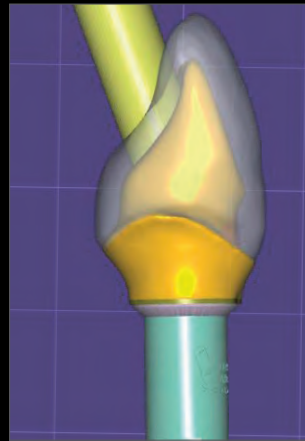
- PAD^e: Easy to insert;
- maximum freedom and ergonomics of the superstructures;
- pink anodization for better mimetism;
- their limited size makes them suitable for rehabilitations of partially or totally edentulous arches in reduced horizontal and vertical spaces.



Solutions for angled screw hole

In the case of implants in the aesthetic area, which need an angulation between the implant and the post axis, it is recommended to place the implants deeper in order to be able to use bases or posts for luting with angled screw hole, so that the screw has a palatal access, consistent with the aesthetic needs.

Solutions for luting are available with or without support base for the crown, in order to leave the clinician the option to choose where to close the crown.







Courtesy of dr. Paolo Nardinocchi



L-MDT for angled screw hole


Impression and model phase

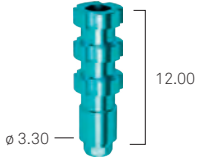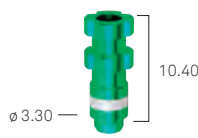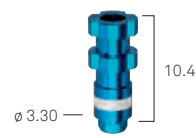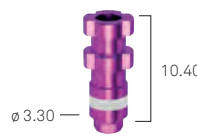Pick-up transfers








straight emergence			
ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
			
A-TRA-330	L-TRA-380	L-TRA-425	L-TRA-500

VTRA2-180-15 screw included.



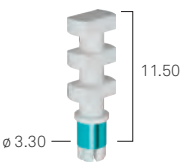
Pull-up transfers


straight emergence
unique platform

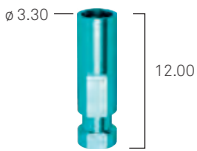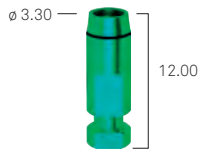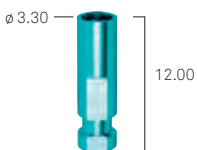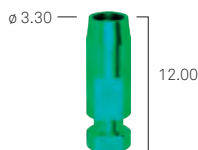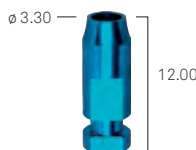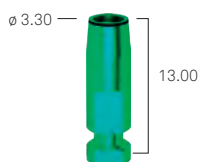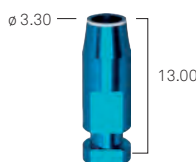A-TRAP-330

Analogs




ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
Short Neck			
			-
A-ANA-330	LS-ANA-380	LS-ANA-425	
Regular Neck			
			
A-ANA-330	L-ANA-380	L-ANA-425	L-ANA-500
Long Neck			
-			
	LL-ANA-380	LL-ANA-425	LL-ANA-500

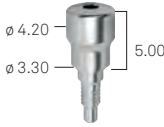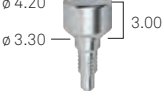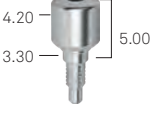Healing abutments

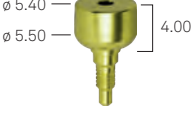Healing abutments					
transgingival h	2.00 mm	3.00 mm	4.00 mm	5.00 mm	7.00 mm
with straight emergence, unique platform			-		
	A-TMG-330-2	A-TMG-330-3	-	A-TMG-330-5	A-TMG-330-7
with straight emergence with coronal enlargement, unique platform			-		-
	A-TMTCR-330-2	A-TMTCR-330-3	-	A-TMTCR-330-5	-
with wide profile, unique platform					
	A-TMGB-330-2	A-TMGB-330-3	A-TMGB-330-4		

Healing abutment	
unique platform	
titanium Gr. 5	
A-TMG-MEFL-330	

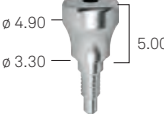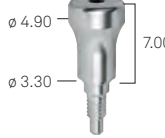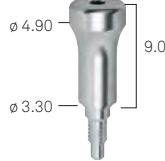REGISTRA technique

Anatomical transfer screw for REGISTRA technique			
h	5.00 mm	7.00 mm	9.00 mm
anatomical transfer screw for REGISTRA technique, unique platform			
	A-TO-330-5	A-TO-330-7	AS-TO-330-9


Transfer cap	
Universal cap for REGISTRA transfer screws	
	
CAP-TO	

Prama IN healing abutments (closing around the implant neck)



in Titanium					
transgingival h	closing around the neck h	ø 3.80 	ø 4.25 	ø 5.00 	
2.50	0.50*	5.00  2.50 L-TMG-380-25-05	5.40  2.50 L-TMG-425-25-05	6.10  2.50 L-TMG-500-25-05	
3.50	0.50*	5.00  3.50 L-TMG-380-35-05	5.40  3.50 L-TMG-425-35-05	6.10  3.50 L-TMG-500-35-05	
2.50	1.50**	5.00  2.50 L-TMG-380-25-15	5.40  2.50 L-TMG-425-25-15	6.10  2.50 L-TMG-500-25-15	
3.50	1.50**	5.00  3.50 L-TMG-380-35-15	5.40  3.50 L-TMG-425-35-15	6.10  3.50 L-TMG-500-35-15	
in REEF resin					
transgingival h	closing around the neck h	ø 3.80 	ø 4.25 	ø 5.00 	
6.00	0.50*	5.00  6.00 L-TMGPF-380-05	5.40  6.00 L-TMGPF-425-05	6.10  6.00 L-TMGPF-500-05	
6.00	1.50**	5.00  6.00 L-TMGPF-380-15	5.40  6.00 L-TMGPF-425-15	6.10  6.00 L-TMGPF-500-15	

Prama IN healing abutments in REEF resin: VM2-180 fixation screw included
*Only for Regular and Long neck. | *Only for Regular neck.

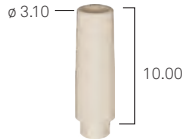Provisional posts

in REEF resin		
unique platform		
	A-PPF-330-EX	in REEF resin, repositionable
	A-PPF-330	in REEF resin, rotating

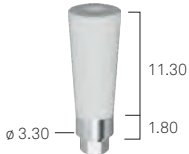
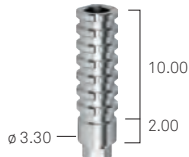
VM2-180 fixation screw included.

Simple in PEEK and Titanium Gr. 5		
unique platform		
	A-MPSC-330	Simple in PEEK, base in Titanium Gr. 5, repositionable
	A-MPSCI-330-EX	Simple in Titanium Gr. 5, repositionable
	A-MPSCI-330	Simple in Titanium Gr. 5, rotating

VM2-180 fixation screw included.



In Titanium with castable sleeve		
unique platform		
	A-CTI-330	rotating
	A-CCI-S	spare castable PMMA sleeve, fixation screw not included

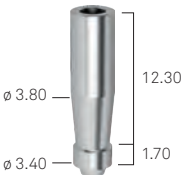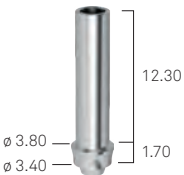VM2-180 fixation screw included.

Posts for luting technique		
unique platform		
	L-CT-340-ROT	Gr. 5 titanium posts, h 14.00 mm, straight emergence
	L-CCRCO-340-ROT	Cobalt chrome titanium posts, h 14.00 mm, straight emergence
	L-CTR-340-ROT	Gr. 5 titanium posts, h 14.00 mm, anatomical emergence
	L-CCRCOR-340-ROT	Cobalt chrome posts, h 14.00 mm, anatomical emergence

L-VMS-180 screw included.

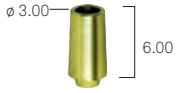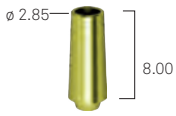Premade posts

Straight				
unique platform		post h		
		6.00	L-MD-340-6	repositionable
		8.00	L-MD-340-8	repositionable
		6.00	L-MD-340-6-ROT	rotating
		8.00	L-MD-340-8-ROT	rotating





L-VMS-180 screw included.

Straight, marginless	
unique platform	post h
	15.00 A-MD-330-15 repositionable

L-VMS-180 screw included.


Straight, marginless	
unique platform	post h
	15.00 A-MD-330-15-ROT rotating

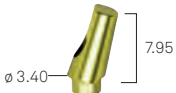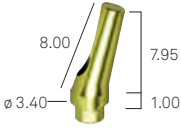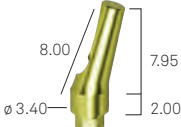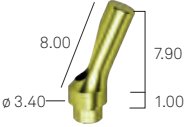L-VMS-180 screw included.

Angled				
unique platform		transgingival h		
		0.00	L-MA15-340*	
		1.00	L-MA15-340-1**	angled at 15°, repositionable
		2.00	L-MA15-340-2**	
		1.00	L-MA25-340-1**	angled at 25°, repositionable
		2.00	L-MA25-340-2**	





**L-VMS-180 fixation screw included, for *L-MA15-340 post: L-VM-180 screw included. The screw L-VM-180 has to be used with the screwdriver L-HSM-CA.

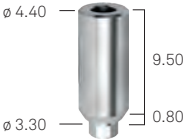Preparable posts

Preparable posts	
unique platform	
straight emergence	 A-MF-330
wide emergence	 A-MFS-330

VM2-180 fixation screw included.

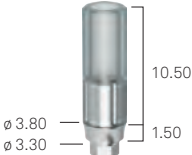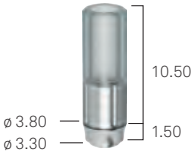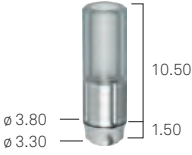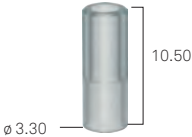Castable posts

With metal base	
unique platform	
cobalt chrome base, repositionable	 A-UCRCO-330-EX
cobalt chrome base, rotating	 A-UCRCO-330
titanium base, rotating <i>only for titanium fusions</i>	 A-UCTR-330-EX
spare castable sleeve, fixation screw not included	 A-CCUCR-330

VM2-180 fixation screw included.

Fully castable	
unique platform	
repositionable, straight emergence	 A-CC-330-EX
rotating, straight emergence	 A-CC-330



VM2-180 fixation screw included.

Solutions with angled screw hole







These solutions have a distinctive morphology, featuring a lateral window or an adjustable element that allows for the displacement of the hole for the passing screw to either the palatal or lingual position, thus avoiding unaesthetic vestibular holes.

Angled screw hole Interfase			
unique platform			
cobalt chrome base, repositionable		10.00 fully castable, repositionable	
	ø 3.30 PD3PKH330/CC		ø 3.30 PD3PKH330/P
cobalt chrome base, rotating		10.00 fully castable, rotating	
	ø 3.30 PD3PKR330/CC		ø 3.30 PD3PKR330/P

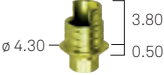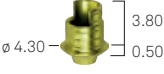

Angled screw hole Interfase	
unique platform	
repositionable	
	ø 4.30 3.80 0.50 IND3PKH330/TIA
rotating	
	ø 4.30 3.80 0.50 IND3PKR330/TIA

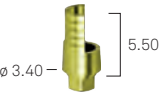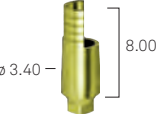
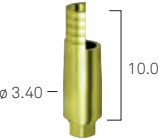
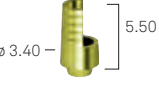
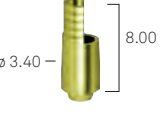
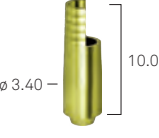
A-VMA-180 fixation screw not included.

Angled screw hole marginless Interfase			
unique platform			
cone height	6.00 mm	8.00 mm	10.00 mm
repositionable			
	ø 3.40 5.50 L-MDT-340-6	ø 3.40 8.00 L-MDT-340-8	ø 3.40 10.00 L-MDT-340-10
rotating			
	ø 3.40 5.50 L-MDT-340-6-ROT	ø 3.40 8.00 L-MDT-340-8-ROT	ø 3.40 10.00 L-MDT-340-10-ROT

For these components use the standard screw VM2-180 or the angled hole screw A-VMA-180. Screws are not included.




ATTENTION: The angled hole screw must always be ordered separately because it is not included in the Interfase package. If you do not have the appropriate screwdriver DSPDCLH-24 or DSPDCLH-32, please order it separately.

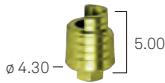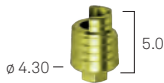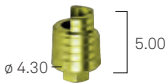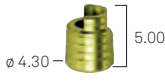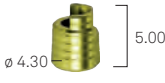Prama IN T-Connect closing around the implant neck at 0.50 mm

cone height		ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
repositionable	5.00 mm	 L-TCIN-380-05	 L-TCIN-425-05	 L-TCIN-500-05
rotating	5.00 mm	 L-TCIN-380-05-ROT	 L-TCIN-425-05-ROT	 L-TCIN-500-05-ROT




Only for Regular and Long Neck | For these components use the standard VM2-180 fixation screw or the angled screw hole one A-VMA-180. Screws are not included.

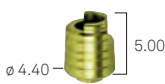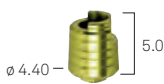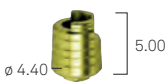Prama IN T-Connect closing around the implant neck at 1.00 mm, Regular Neck

cone height		ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
repositionable	5.00 mm	 L-TCIN-380-10	 L-TCIN-425-10	 L-TCIN-500-10




Only for Regular Neck | For these components use the standard VM2-180 fixation screw or the angled screw hole one A-VMA-180. Screws are not included.

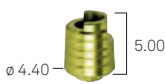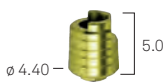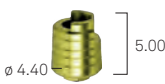Prama IN T-Connect closing around the implant neck at 1.00 mm, Long Neck

cone height		ø 3.80	ø 4.25	ø 5.00
repositionable	5.00 mm	 LL-TCIN-380-10	 LL-TCIN-425-10	 LL-TCIN-500-10




Only for Long Neck | For these components use the standard VM2-180 fixation screw or the angled screw hole one A-VMA-180. Screws are not included.

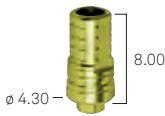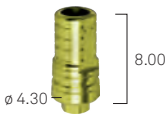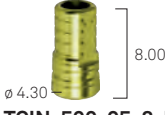Cuttable Prama IN T-Connect T-CUT closing around the implant neck at 0.50 mm

cone height		ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
8.00 mm		 L-TCIN-380-05-8	 L-TCIN-425-05-8	 L-TCIN-500-05-8
8.00 mm		 L-TCIN-380-05-8-ROT	 L-TCIN-425-05-8-ROT	 L-TCIN-500-05-8-ROT



Only for Regular Neck | For these components use the standard VM2-180 fixation screw or the angled screw hole one A-VMA-180. Screws are not included.

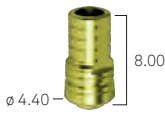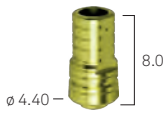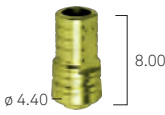Cuttable Prama IN T-Connect T-CUT closing around the implant neck at 1.00 mm, Regular Neck

cone height		ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
8.00 mm		 L-TCIN-380-10-8	 L-TCIN-425-10-8	 L-TCIN-500-10-8

Only for Regular Neck | For these components use the standard VM2-180 fixation screw or the angled screw hole one A-VMA-180. Screws are not included.

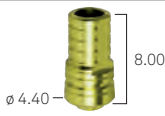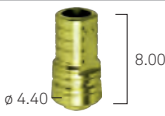Cuttable Prama IN T-Connect T-CUT closing around the implant neck at 1.00 mm, Long Neck




cone height		ø 3.80	ø 4.25	ø 5.00
8.00 mm		 LL-TCIN-380-10-8	 LL-TCIN-425-10-8	 LL-TCIN-500-10-8

Only for Long Neck | For these components use the standard VM2-180 fixation screw or the angled screw hole one A-VMA-180. Screws are not included.

ATTENTION: The angled hole screw must always be ordered separately because it is not included in the Interfase package. If you do not have the appropriate screwdriver DSPDCLH-24 or DSPDCLH-32, please order it separately.

Prosthetic line PAD^e (Aesthetic screw-retained disparallel prosthesis)

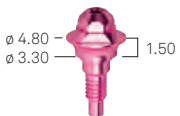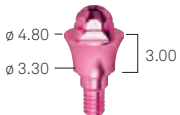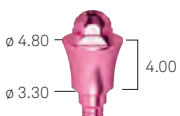Straight PADe abutments, for direct screwing on implants with Collex One connection

transgingival h	unique platform
1.50 mm	 A-PADE-AD330-15
3.00 mm	 A-PADE-AD330-30
4.00 mm	 A-PADE-AD330-40

Carrier PADR-VTRAL-140-M included: for intraoral screwing use the AVV2-ABUT screwdriver.

Angled abutments

	unique platform
Angled at 17° transgingival h 3.00 mm	 A-PADE-AA330-173
Angled at 17° transgingival h 5.00 mm	 A-PADE-AA330-175
Angled at 30° transgingival h 3.00 mm	 A-PADE-AA330-303
Angled at 30° transgingival h 5.00 mm	 A-PADE-AA330-305
Angled at 45° transgingival h 3.00 mm	 A-PADE-AA330-453
Angled at 45° transgingival h 5.00 mm	 A-PADE-AA330-455




Screw PADR-VM-180 included.
Carrier PADR-VTRAL-140-M included: for intraoral screwing use the appropriate screwdriver from the HSM series.

PADe prosthetic components

Healing caps	
Titanium healing cap	 PADR-CG
PEEK healing cap	 PADR-CGP

PADR-CG: PAD-VP-140 screw included. | PADR-CGP: PAD-VCGP-140 screw included.

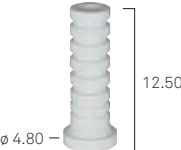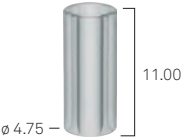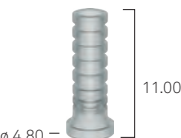Components for prosthesis relining and luting technique

PEEK sleeve, non engaging	 PADR-CP
Titanium sleeve, non engaging	 PADR-CT
PMMA castable sleeve	 PADR-CCEM
PMMA castable sleeve, non engaging	 PADR-CC
Titanium sleeve without shoulder, for aesthetic reconstructions	 PADR-BAS

Sleeves in PEEK, Titanium Gr. 5 and PMMA: screw PAD-VP-140 included.

Castable post with metal base

PMMA post with cobalt-chrome base, rotating	 PADR-UCRCO
---	--

PAD-VP-140 screw included.

Impression phase

Snap-on cap



PADR-CAP-EX

Transfer pick-up, rotating



PADR-TRA

Analog



PADR-ANA

Pick up transfer: PAD-VTRAL-140 screw included.

Components for CAD-CAM prosthesis

Intraoral scanbody



PADR-INT-CAMTRA-L

Intraoral scanbody for direct screw connection**



Y-AVV-SB-DIG-NX

Intraoral scanbody for direct screw connection, long*



Y-AVV-SBL-DIG-NX

Digital analog for abutment



PADR-ANA-DIG

Sleeve for luting technique non engaging for abutments with cone h 4



PADR-TC-M-4

Sleeve for luting technique non engaging for abutments with cone h 8 and reduction cut



PADR-TC-M-8

Sleeve for luting technique non engaging for abutment for angled screw hole



PADR-TCA-M*

Sleeve for luting technique repositionable for abutments with cone h 4



PADR-TC-S-4

Sleeve for luting technique repositionable for abutments with cone h 8 and reduction cut



PADR-TC-S-8

Sleeve for luting technique repositionable for abutment for angled screw hole



PADR-TCA-S*

Castable sleeve for PADR-TC-M and PADR-TC-M-8, non engaging



PADR-TC-CC-M-8

Engaging castable sleeve for PADR-TCA-M



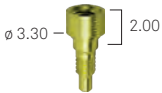
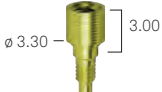
PADR-TC-CC-S-8




Scanbody and sleeve for luting technique: PAD-VP-140 screw included.
Digital analog: VADIG-180 screw included.

*Non engaging sleeve for luting technique for abutment for angled screw hole: screw PAD-VPA-140 not included.




**Use the dedicated screwdriver YV-CV-MAN-M.

PLAIN B-Space prosthetic range

Abutment	
transgingival h	unique platform
2.00 mm	
	A-PLAIN-ABU330-2
3.00 mm	
	A-PLAIN-ABU330-3
4.00 mm	
	A-PLAIN-ABU330-4

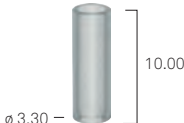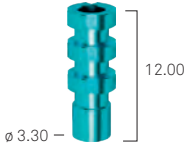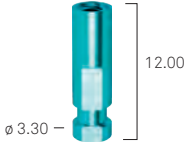
Components for overstructure	
	unique platform
healing cap	
	A-PLAIN-CG330
castable sleeve	
	A-PLAIN-CC330
Gr. 5 titanium sleeve	
	A-PLAIN-CT330
Pick-up transfer	
	A-PLAIN-TRA-330
analog	
	A-PLAIN-ANA-330

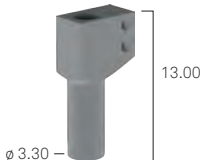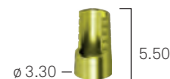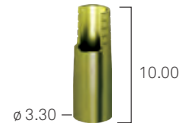Castable and titanium sleeves: A-PLAIN-VP200 fixation screw included.
 Pick-up transfers: A-PLAIN-VTRA200 fixation screw included.


Components for digital prosthesis	
	unique platform
Titanium scanbody for laboratory and intraoral use, for PLAIN B-Space abutment	
	A-INT-CAMTRA-PLAIN-L
Digital analog for PLAIN B-Space abutment	
	A-PLAIN-ANA-330-DIG
PLAIN B-Space sleeve for luting technique h. 6.00 mm	
	A-PLAIN-MDT-6
PLAIN B-Space sleeve for luting technique h. 10.00 mm	
	A-PLAIN-MDT-10




Scanbody: A-PLAIN-VP200 screw included.
 Digital analog: VADIG-180 screw included.
 PLAIN B-Space sleeves for luting technique: screw not included.
 A-PLAIN-VP200 screw for straight hole or ABU-VMA-200 for angled screw hole available for separate ordering.

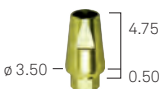Conico conometric prosthesis

Straight	
h	unique platform
0.50	 A-MD-TS-330-05
1.00	 A-MD-TS-330-10
2.00	 A-MD-TS-330-20
3.00	 A-MD-TS-330-30

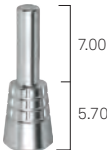VM2-180 fixation screw included.

Conico components	
cap for partial or full-arch fixed prosthesis	 CAP2-TS-DEF
non-rotating cap for single element	 CAP2-TS-IND
cap for the removal	 CAP2-TS-REM
PEEK transfer cap for direct impression on Conico abutment	 CAP2-TS-IMP
Conico post analog	 ANA2-MD-TS
cap with pin, to be used for intraoral wel (welding bars available upon request)	 CAP2-TS-PIN

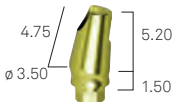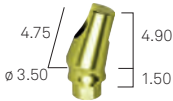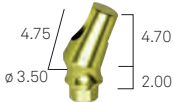Angled	
	unique platform
angled at 7.5°	 A-MA07-TS-330-2
angled at 15°	 A-MA15-TS-330-2
angled at 22.5°	 A-MA20-TS-330-2

VM2-180 fixation screw included.

Accessories for Conico system	
parallelometer	 PAR-PP
parallelometer insert	 PAR-INS-AS
parallelometer driver	 CAP2-TS-PAR


Remember to always purchase the parallelometer with the insert relating to the implant connection on which it is to be used.

Components for Conico digital prosthesis	
	unique platform
Titanium scanbody for laboratory and intraoral use, for Conico abutment	 C-INT-CAMTRA-L
Digital analog for Conico abutment	 ANA2-MD-TS-DIG

VADIG-180 VM2-180 screw included.

Components for digital prosthesis

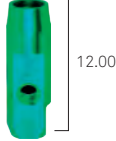
Components for impression and model phase - platform level

	ø 3.30 ●	ø 3.80 ●	ø 4.25 ●	ø 5.00 ●
scanbody*	 L-INT-CAMTRA-L	 L-INT-CAMTRA-L	 L-INT-CAMTRA-L	 L-INT-CAMTRA-L
digital analog short neck**	 A-ANA-330-DIG	 LS-ANA-380-DIG	 LS-ANA-425-DIG	-
digital analog regular neck**	 A-ANA-330-DIG	 L-ANA-380-DIG	 L-ANA-425-DIG	 L-ANA-500-DIG
digital analog long neck**	-	 LL-ANA-380-DIG	 LL-ANA-425-DIG	 LL-ANA-500-DIG

* VM2-18 screw included. | ** VADIG-180 screw included.

Components for impression and model phase - abutment level

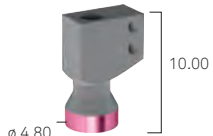unique platform

for P.A.D.r abutments	 10.00 ∅ 4.80 PADR-INT-CAMTRA-L	 10.00 PADR-ANA-DIG
for PLAIN B-Space abutments	 13.00 ∅ 3.30 A-INT-CAMTRA-PLAIN-L	 10.00 A-PLAIN-ANA-330-DIG
for Conico abutments	 10.00 ∅ 3.50 C-INT-CAMTRA-L	 10.00 ANA2-MD-TS-DIG



Screws and screwdrivers for digital analogs

unique platform

					
VADIG-180 screw for digital analog M1.8	PIN-ANA-DIG pin for digital analog	INS-ANA-140 carrier for digital analog M 1.4 P.A.D.r	INS-ANA-180 carrier for digital analog M1.8	INS-ANA-180-L carrier for digital analog M1.8, long	INS-ANA-200 carrier for digital analog M 2.0 Conico

Straight premade posts





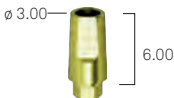
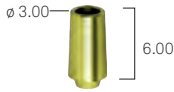
unique platform

post h	6.00	8.00
repositionable	 L-MD-340-6	 L-MD-340-8
rotating	 L-MD-340-6-ROT	 L-MD-340-8-ROT

L-VMS-180 screw included.

Angled screw hole marginless Interfase

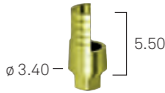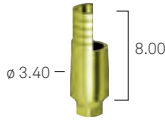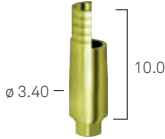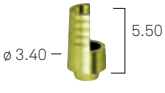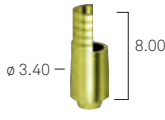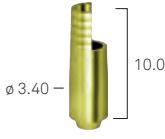unique platform









cone h	6.00 mm	8.00 mm	10.00 mm
repositionable	 L-MDT-340-6	 L-MDT-340-8	 L-MDT-340-10
rotating	 L-MDT-340-6-ROT	 L-MDT-340-8-ROT	 L-MDT-340-10-ROT

For these components use the standard screw VM2-180 or the angled hole screw A-VMA-180. Screws are not included.

T-Connect with cementing cone

unique platform

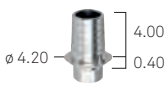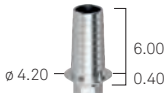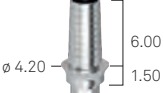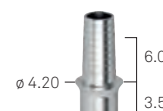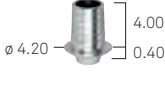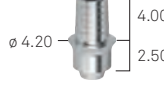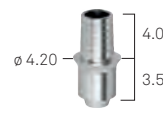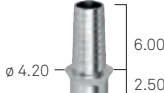repositionable

cone h 4.00 mm	 A-BASTZR-S-330-4	 A-TCR-S-330-4-1	 A-TCR-S-330-4-2	 A-TCR-S-330-4-3
cone h 6.00 mm	 A-BASTZR-S-330-6	 A-TCR-S-330-6-1	 A-TCR-S-330-6-2	 A-TCR-S-330-6-3
rotating				
cone h 4.00 mm	 A-BASTZR-M-330-4	 A-TCR-M-330-4-1	 A-TCR-M-330-4-2	 A-TCR-M-330-4-3
cone h 6.00 mm	 A-BASTZR-M-330-6	 A-TCR-M-330-6-1	 A-TCR-M-330-6-2	 A-TCR-M-330-6-3

VM2-180 screw included.

Chairside T-Connect with cementing cone

unique platform cone h 5.00

transgingival h



A-TCS-330S-C

1.50 mm



A-TCS-330S-C-15

2.50 mm



A-TCS-330S-C-25

VM2-180 screw included.

Sleeves for wax-up modeling on T-Connect

unique platform

total h 12.00 mm
cementing cone
h 4.00 mm



A-CCBAS-330-4

total h 12.00 mm
cementing cone
h 6.00 mm

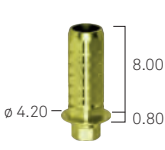


A-CCBAS-330-6


T-Connect Unica

unique platform


repositionable




A-TC-UNI-330

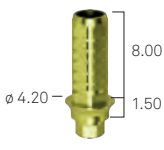

A-TC-UNI-330-1

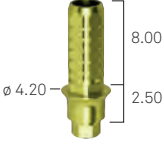


A-TC-UNI-330-2

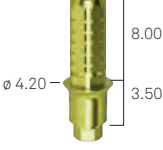


A-TC-UNI-330-3

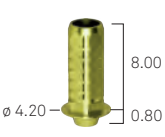
rotating




A-TC-UNI-330-ROT

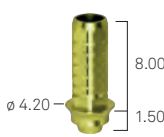

A-TC-UNI-330-1-ROT

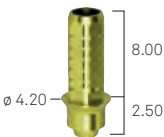

A-TC-UNI-330-2-ROT

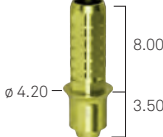

A-TC-UNI-330-3-ROT

VM2-180 screw included.

ATTENTION: The angled hole screw must always be ordered separately because it is not included in the Interfase package. If you do not have the appropriate screwdriver DSPDCLH-24 or DSPDCLH-32, please order it separately.

Cutting guide and holder for T-Connect Unica



Cutting guide holder S
(small) for T-Connect Unica
left side

TC-UNI-HOLD4-SX

Cutting guide holder S
(small) for T-Connect Unica
right side

TC-UNI-HOLD4-DX

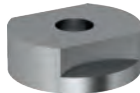
Cutting guide holder L (large)
for T-Connect Unica
right side

TC-UNI-HOLD6-DX

Cutting guide holder L (large)
for T-Connect Unica left side


TC-UNI-HOLD6-SX

4 mm cutting guide for
T-Connect Unica



TC-UNI-DIMA-H4

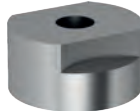6 mm cutting guide for
T-Connect Unica



TC-UNI-DIMA-H6

Cutting guide for angled hole
for T-Connect Unica



TC-UNI-DIMA-ANG

Please note: STL file available for independently 3D printing of the holders.



See how easy it is to precisely
customize the T-Connect Unica
with the cutting guides and the
corresponding holder.

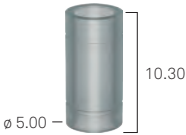
T-Connect Unica accessories

Intraoral testing O-ring for
T-Connect Unica




TC-UNI-ORING

Castable sleeve for
T-Connect, rotating




TC-UNI-CC

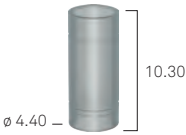Small castable sleeve
for T-Connect Unica,
repositionable



TC-UNI-CC-S

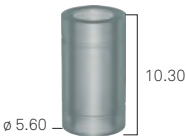Large castable sleeve
for T-Connect Unica,
repositionable



TC-UNI-CC-L

O-Ring inserter



TC-UNI-ORING-INS

Cutting guide and holder for T-Connect Prama IN



Cutting guide holder S (small) for T-Connect Prama IN left side
TCIN-HOLD4-SX

Cutting guide holder S (small) for T-Connect Prama IN right side
TCIN-HOLD4-DX

Cutting guide holder L (large) for T-Connect Prama IN right side
TCIN-HOLD6-DX

Cutting guide holder L (large) for T-Connect Prama IN left side
TCIN-HOLD6-SX

4 mm cutting guide for T-Connect Unica
TCIN-DIMA-H5

6 mm cutting guide for T-Connect Unica
TCIN-DIMA-H6

Cutting guide for angled hole for T-Connect Unica
TCIN-DIMA-ANG

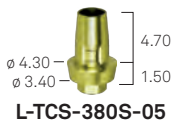
Please note: STL file available for independently 3D printing of the holders.

Chairside Prama IN T-Connect closing around the neck at 0.50 mm*

ø 3.80 ●

ø 4.25 ●

ø 5.00 ●



L-TCS-380S-05



L-TCS-425L-05



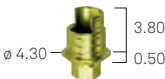
L-TCS-500L-05

*Only for Regular and Long Neck. L-VMS screw included.

Angled screw hole Interfase


unique platform

repositionable



IND3PKH330/TIA

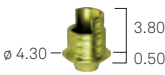rotating



IND3PKR330/TIA

A-VMA-180 screw not included.

T-Connect Prama IN accessories

T-Connect Prama IN
O-ring



TCIN-ORING

O-Ring inserter



TC-UNI-ORING-INS

unique platform

repositionable



in Titanium

A-CAMT-330-12-SW



in Cobalt Chrome

A-CAMC-330-12-SW

VM2-180 screw included.

Ø 3.80 ●

repositionable



in Titanium

L-CAMT-38005-12-SW



L-CAMT-42505-12-SW



L-CAMT-50005-12-SW

VM2-180 screw included. | *Only for Regular and Long Neck.

unique platform

repositionable






in Titanium

A-CAMTABU330-CRE

VM2-180 screw included.

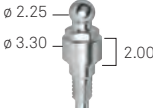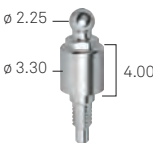Solutions for overdenture

Ball attachments	
unique platform	
 A-AS-330-1	
 A-AS-330-2	
 A-AS-330-4	
Analog for ball attachments	
unique platform	
 ANAS	
Caps for ball attachments	
unique platform	
polyamide cap	 040CRNSNDR8-SM*
steel container for polyamide cap external diameter 4.80 mm	 041CAN-SM*
Gr. 5 titanium cap with cap in two parts, retention spring titanium and plastic, plastic mounting ring	 CAP-TIT-1
spare plastic ring	 AN-CAP-TIT-1
spare steel spring, average hardness	 MOL1-CAP-TIT-1
spare steel spring, soft	 MOL2-CAP-TIT-1
gold alloy cap, with 3 plastic rings and a tin spacer for the lab	 CAP-1
* Pack of 6 pieces.	

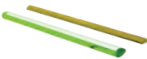
Overdenture on bars	
castable bar, l 5.00 cm, h 3.00 mm, thickness 2.20 mm, ovoid-shaped profile with spacer	 BARC-CAV-TIT
divisible bar attachment in titanium for oval bars, h 3.00 mm, thickness 2.20 mm	 CAV-TIT
castable bar, l 5.00 cm, 2.20 mm	 BARC
bar attachment in gold alloy for round bars 2.20 mm	 CAV-375

Locator abutment

transgingival h

unique platform

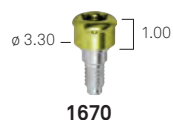1.00 mm



1670

2.00 mm



1671

3.00 mm



1672

4.00 mm

-

5.00 mm



1674

6.00 mm

-

Accessories for Locator abutment

unique platform

4 aluminium analogs



8530

4 aluminium transfers,
4 retainers cod. 8515 included



8505

4 polyethylene parallelism pins



8517

steel plate AISI 316 L



9530

complete Locator Core Tool



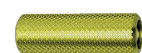
8393

steel push rod for
inserting the retainers



8397

steel driver for abutment
screwing/unscrewing



8390

retention jacket for the driver, to
transfer the Locator abutments
into the oral cavity



8394

Retainers and spacer rings for Locator abutments

unique platform

low retention capacity for impression taking,
black, in polyethylene,
4 pcs packaging



8515

retention 5 lb = 2268 g, in nylon,
disparallelisms up to 10° per side,
4 pcs packaging



8524

retention 3 lb = 1361 g, in nylon,
disparallelisms up to 10°
per side, 4 pcs packaging



8527

retention 1,5 lb = 680 g, in nylon,
disparallelisms up to 10°
per side, 4 pcs packaging



8529

retention 4 lb = 1814 g, in nylon,
disparallelisms up to 20°
per side, 4 pcs packaging



8547

retention 2 lb = 907 g, in nylon,
disparallelisms up to 20°
per side, 4 pcs packaging



8915

retention 1 lb = 450 g, in nylon,
disparallelisms up to 20°
per side, 4 pcs packaging



8548

spacer rings in silicon rubber,
20 pcs packaging



8514

kit: 2 Gr. 5 titanium caps, 2 spacer rings,
2 black retainers and 2 nylon retainers for
each of the 3 different retention capacities



8519-2

kit: 2 Gr. 5 titanium caps, 2 spacer rings,
2 black retainers and 2 nylon retainers for
each of the 3 different retention capacities
meant for disparallelisms



8540-2

kit: 2 Gr. 5 steel caps, 2 spacer rings,
2 black retainers and 2 nylon retainers for
each of the 3 different retention capacities



8550-2

Locator R-Tx abutment*	
transgingival h	unique platform
1.00 mm	 31110-01
2.00 mm	 31110-02
3.00 mm	 31110-03
4.00 mm	 31110-04
5.00 mm	 31110-05
6.00 mm	 31110-06

* All-in-One packaging that includes an abutment, a spacer ring, a cap, and an insert for each retention (gray, blue, pink, white).


Caps and retainer for Locator R-Tx abutment	
	unique platform
pack of 4 caps and 4 black retainers for impression taking	 30013
grey, without retention, pack of 4 pcs	 30001
blue, low retention, pack of 4 pcs	 30002
pink, medium retention, pack of 4 pcs	 30003
white, high retention, pack of 4 pcs	 30004
black, for the impression taking, pack of 4 pcs	 30012

Accessories for Locator R-Tx abutment	
	unique platform
aluminium analogs for abutments $\varnothing$ 3.35 mm pack of 4 pcs	 30014
aluminium analogs for abutments $\varnothing$ 4.00, pack of 4 pcs	 30015
aluminium analogs for abutments $\varnothing$ 5.00, pack of 4 pcs	 30016
transfers with low retention, pack of 4 pcs	 30017
spacer rings for relining phase, pack of 4 pcs	 30018
Locator core tool for Locator R-Tx	 30021
RA screwdriver for Locator R-Tx, hex 1.25 mm	 08008
digital screwdriver for Locator R-Tx, hex 1.25 mm	 08007

Implant site preparation sequences

Implant site preparation sequences for Prama Power implants with Shorty drills

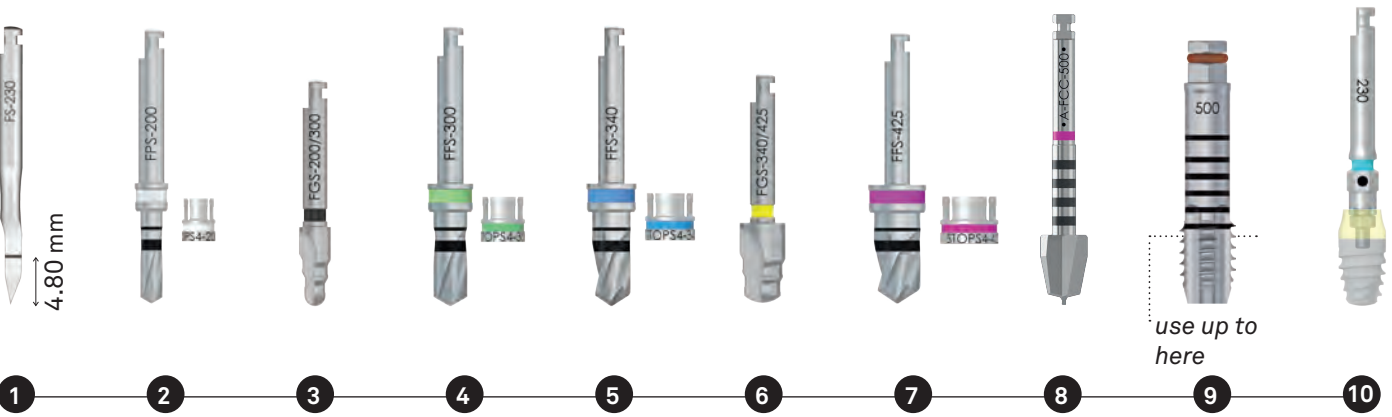The graphic sequence shown illustrates an example of preparation for a Prama Power implant with a height of 6.00 mm and a diameter of 5.00 mm. The use of the stops corresponding to the desired preparation depth, included in the kit, is also recommended.




The instrument sequence shown below is applicable to all implants in the Prama Power range.

	ø 3.30 mm	ø 3.80 mm	ø 4.25 mm	ø 5.00 mm
1 Precision drill		✓	✓	✓
2 Cylindrical pilot drill ø 2.00 mm	✓	✓	✓	✓
3 Intermediate drill ø 2.00/3.00 mm	✓	✓	✓	✓
4 Final cylindrical drill ø 3.00 mm		✓	✓	✓
5 Final cylindrical drill ø 3.40 mm			✓	✓
6 Intermediate drill ø 3.40/4.25 mm				✓
7 Cylindrical final drill ø 4.25 mm				✓
8 Countersink ø 5.00 mm		✓	✓	✓
9 Bone tap	A-MS-330	A-MS-380	A-MS-425	A-MS-500
10 Short Easy Insert driver	✓	✓	✓	✓

The use of the highlighted instruments is left to the clinician's judgment and experience, depending on the characteristics of the recipient bone.

Implant site preparation sequences for Prama Power implants

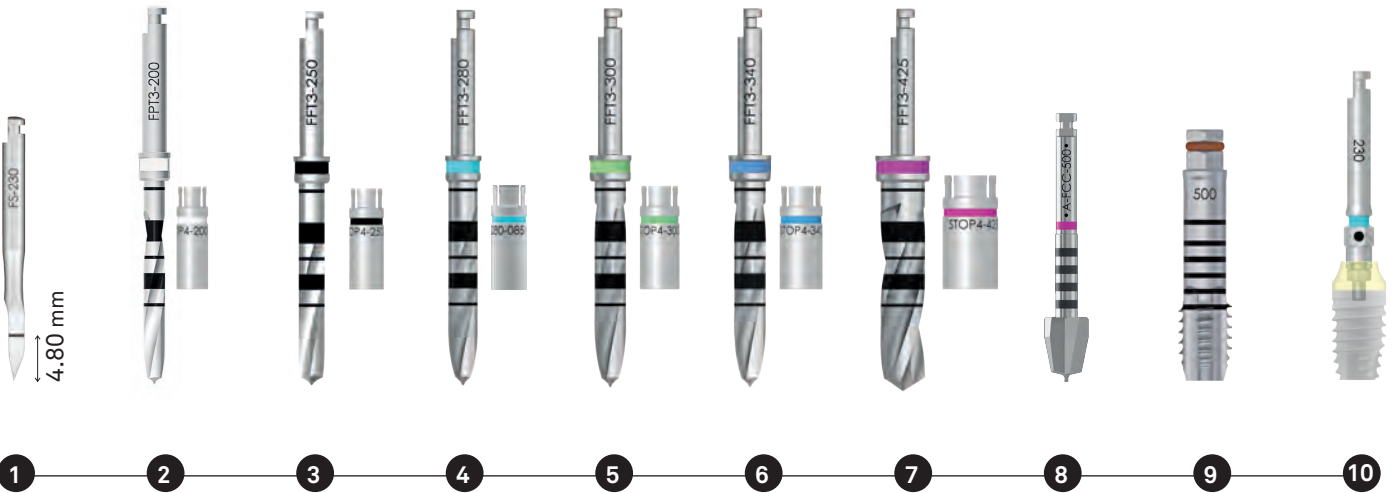The graphic sequence shown illustrates an example of site preparation for a Prama Power implant with a length of 8.50 mm and a diameter of 5.00 mm. The use of the corresponding depth stops, included in the kit, is also recommended to achieve the desired preparation depth.



The instrument sequence shown below applies to all implants in the Prama Power range.

		ø 3.30 mm 	ø 3.80 mm 	ø 4.25 mm 	ø 5.00 mm 	
1	Precision drill	✓	✓	✓	✓	
2	Cylindrical pilot drill ø 2.00 mm	✓	✓	✓	✓	
3	Cylindrical final drill ø 2.50 mm	✓	✓	✓	✓	
4	Cylindrical final drill ø 2.80 mm	✓	✓	✓	✓	
5	Cylindrical final drill ø 3.00 mm		✓	✓	✓	
6	Cylindrical final drill ø 3.40 mm			✓	✓	
7	Cylindrical final drill ø 4.25 mm				✓	
8	Countersink ø 5.00 mm					✓
9	Bone Tap	A-MS-330	A-MS-380	A-MS-425	A-MS-500	
10	Easy Insert Short Driver	✓	✓	✓	✓	

The use of the highlighted instruments is left to the clinician's discretion, depending on the characteristics of the recipient bone.

Implant site preparation sequences for Prama Power RF implants with Shorty drills

The graphic sequence shown illustrates an example of site preparation for a Prama Power implant with a length of 6.00 mm and a diameter of 5.00 mm. The use of the corresponding depth stops, included in the kit, is also recommended to achieve the desired preparation depth.



The instrument sequence shown below applies to all 6 mm-length implants in the Prama Power RF range.

	ø 4.25 mm ●	ø 5.00 mm ●
1 Precision drill	✓	✓
2 Cylindrical pilot drill ø 2.00 mm	✓	✓ *
3 Conical final drill ø 4.10 mm	✓	✓ *
4 Conical final drill ø 5.00 mm		✓ *
5 Countersink ø 5.00 mm	✓	✓
6 Easy Insert Short driver	✓	✓

*Instruments included in the L-INTEGRA-060 set, also available for individual purchase.

Implant site preparation sequences for Prama Power RF and Prama Power RF SL implants

The graphic sequence refers to the $\varnothing$ 5.00 mm implant with a length of 8.50 mm.
For longer implant lengths, the corresponding final drills should be used, identified by the matching suffix; for example, SH-FK425-100 for the 10 mm length, SH-FK425-115 for the 11.50 mm length, and so on.
The use of the corresponding depth stops, included in the kit, is also recommended to achieve the desired preparation depth.



The instrument sequence shown below applies to all implants in the Prama Power RF and Prama Power RF SL ranges.

		$\varnothing$ 3.80 mm ●	$\varnothing$ 4.25 mm ●	$\varnothing$ 5.00 mm ●
1	Precision drill	✓	✓	✓
2	Cylindrical pilot drill $\varnothing$ 2.00 mm	✓	✓	✓
3	Conical final drill $\varnothing$ 2.50 mm	✓	✓	✓
4	Conical final drill $\varnothing$ 3.80 mm	✓	✓	✓
5	Conical final drill $\varnothing$ 4.25 mm	✓	✓	✓
6	Conical final drill $\varnothing$ 5.00 mm		✓	✓
7	Final drill $\varnothing$ 5.00 mm			✓
8	Bone tap	SH-MS-380-CA	SH-MS-425-CA	SH-MS-500-CA
9	Easy-Insert driver	✓	✓	✓

The use of the highlighted instruments is left to the clinician's discretion, depending on the characteristics of the recipient bone.

Clinical and laboratory prosthetic screws

Standard fixation screw		
	VM2-180	
	VM2-180-10	10 pcs packaging
Conical support fixation screw		
	L-VMS-180	
	L-VMS-180-10	10 pcs packaging
Full Head Fixation screw		
	L-VM-180	
	L-VM-180-10	confezione da 10 unità
Fixation screw for angled screw hole		
	A-VMA-180	
Screws for Pick-up transfer		
	VTRA2-180-15	h 15.00
	VTRA2-180-20	h 20.00
	VTRA2-180-MAN	for manual screwing
Screw for closed-tray transfers		
	VTRA2-180-10	h 10.00
Screws for PLAIN abutment		
	A-PLAIN-VP200	for PLAIN abutment
	A-PLAIN-VP200-10	for PLAIN abutment, 10 pcs packaging
	A-PLAIN-VTRA200	for PLAIN abutment Pick-up transfers
	ABU-VMA-200	for angled screw hole

Screws for PADe abutments		
	PADR-VM-180	for PADe abutment
	PADR-VM-180-10	for PADe abutment 10 pcs packaging
	PAD-VTRAL-140	for Pick-up transfer, long
	PAD-VTRA-140	for Pick-up transfer, standard
	PAD-VTRAL-140-MAN	for Pick-up transfer, manual screwing
	PAD-VP-140	for prosthetic components on PADe abutments
	PAD-VP-140-10	for prosthetic components on PADe abutments 10 pcs packaging
	PAD-VPA-140	fixation screw for PADR-TCA-M to use with the dedicated screwdriver DSPDCLH-24 and DSPDCLH-32
	PADR-VTRAL-140-M	carrier for PADe

Screws for CAD-CAM prosthesis		
	A-CAMTVABU180	for zirconium oxide single crowns and ECHO overstructures
	PAD-VCAM140	for zirconium oxide overstructure on abutment PADe
	CAMPRON205-10	spare ring washers for the head of the screw
	VADIG-180	locking screw for digital analogs

Recommended torque for tightening prosthetic components	
8-10 Ncm	surgical cover screws, healing abutments, transfer screws
20-25 Ncm	prosthetic screws
25-30 Ncm	direct screwing on implants prosthetic components
8-10 Ncm	transfer screws

Warning
 Fixation screws are always included in the packaging of each prosthetic component, unless otherwise indicated. Please use laboratory screws for laboratory phases and keep the screw included with the prosthetic component for the final fastening in the patient's mouth.

Bibliography

- a. Canullo L., Tallarico M., Pradíes G., Marinotti F., Loi I., Cocchetto R.,
Soft and hard tissue response to an implant with a convergent collar in the esthetic area: preliminary report at 18 months,
Int J Esthet Dent 2017; 12:2-19:
- b. Canullo L, Menini M, Covani U, Pesce P,
Clinical outcomes of using a prosthetic protocol to rehabilitate tissue-level implants with a convergent collar in the esthetic zone: a 3-year prospective study,
J Prosthet Dent. 2019 Jun 18. pii: S0022-3913(19)30077-0. doi: 10.1016/j.prosdent.2018.12.022
- c. Valente N.A., Wu M., Toti P., Derchi G., Barone A,
Concave/convergent versus parallel/divergent implant transmucosal profiles impact on hard and soft peri-implant tissues: a systematic review with meta-analysis,
Int J Prosthodont. Sep/Oct 2020;33(5):553-564. doi: 10.11607/ijp.6726
- d. Cabanes Gumbau G, Pascual-Moscardò A, Penarrocha-Oltra D, Garcia-Mira B, Aizcorbe-Vicente J, Penarrocha-Diago M,
Volumetric variation of peri-implant soft tissues in convergent collar implants and crowns using the biologically oriented preparation technique,
Med Oral Patol Oral Cir Bucal. 2019 Sep 1;24(5):e643-e651. doi: 10.4317/medoral.22946
- e. Agustín-Panadero R, Martínez-Martínez N, L. Fernandez-Estevan, J. Faus-Lopez, F. Sola-Ruiz,
Influence of transmucosal Area Morphology on Peri-implant bone loss in tissue-level implants,
The International Journal of Oral & Maxillofacial Implants, 2019; 34: 947-952. doi: 10.11607/jomi.7329
- f. Marchio V, Derchi G, Cinquini C, Miceli M, Gabriele M, Alfonsi F, Barone A,
Tissue level implants in healthy versus medically compromised patients: a cohort comparative study,
Minerva Stomatologica 2020 mese;69(0):000–000, Doi: 10.23736/S0026-4970.20.04359-9
- g. Prati C, Zamparini F, Canullo L, Pirani C, Botticelli D, Gandolfi MG,
Factors affecting soft and hard tissues around two-piece transmucosal implants: a 3 year prospective cohort study,
J Oral Maxillofac Implants 2020;35, doi: 10.11607/jomi.7778
- h. Agustín-Panadero, R.; Bermúdez-Mulet, I.; Fernández-Estevan, L.; Fernanda Solá-Ruiz, M.; Marco-Pitarch, R.; García-Selva, M.; Zubizarreta-Macho, Á.; León-Martínez, R.,
Peri-Implant Behavior of Tissue Level Dental Implants with a Convergent Neck.,
Int. J. Environ. Res. Public Health 2021, 18, 5232, Doi: 10.3390, ijerph18105232
- i. Palombo D., Rahmati M., Vignoletti F., Sanz-Esporrin J., Haugen H. J., Sanz M.,
**Hard and soft tissue healing around implants with a modified implant neck configuration:
An experimental in vivo preclinical investigation,**
Clin Oral Impl Res. 2021;00:1–15.
- j. Morón-Conejo B, Sanz-Sánchez I, Salido MP, Martínez-Rus F, Pradíes G,
The effect of a convergent transmucosal neck on soft tissues and radiographic outcomes: a 1-year follow-up randomized controlled trial,
Clin Oral Investig 2023 Feb 7. doi: 10.1007/s00784-023-04892-9.
- k. Mandillo-Alonso V, Cascos Sanchez R, Antonaya-Martin JL, Laguna-Martos M,
Evaluation of peri-implant soft and hard tissues behavior in screw-retained crowns by the biologically oriented preparation technique: Ambispective longitudinal analytical study,
J Clin Exp Dent. 2022 Jan; 14(1): e64–e71, doi: 10.4317/jced.58924
- l. Mandillo-Alonso V, Cascos Sanchez R, Antonaya-Martin JL, Laguna-Martos M,
Soft tissue thickness evaluation in screw-retained crowns by the biologically oriented preparation technique,
J Clin Exp Dent 2021 Dec 1;13(12): e1209-e1215, doi: 10.4317/jced.58952



rev. 06-26



Sweden & Martina S.p.A.
Via Veneto, 10
35020 Due Carrare (PD), Italy
Tel. +39.049.9124300
info@sweden-martina.com
www.sweden-martina.com

Sweden & Martina Ltd
Unit 1b
Amberley Court, Whitworth Road
Crawley, West Sussex, RH11 7XL
Toll free 0800 1123575
info.uk@sweden-martina.com

Authorized Distributor
Sweden & Martina Inc.
78 John Miller Way
Unit 1021 - Kearny, New Jersey 07032
Toll free 1-844-8MARTINA (1-844-862-7846)
info.us@sweden-martina.com

www.sweden-martinainc.com

Sweden & Martina Mediterranea S.L.
España - info.es@sweden-martina.com
Sweden & Martina Lda
Portugal - info.pt@sweden-martina.com
Sweden Martina Romania srl - România
info.ro@sweden-martina.com
Sweden Martina France Sarl -France
info.fr@sweden-martina.com

SOME PRODUCTS MAY NOT BE REGULATORY/RELEASED FOR SALE IN ALL MARKETS.

The implants, standard prosthetic components and surgical instruments contained in this issue are Medical devices and are manufactured by Sweden & Martina S.p.A. except for:

- the Locator Abutment and related components which are medical devices manufactured by Zest Anchors, Inc., 2875 Loker Avenue East, Carlsbad, CA 92010, USA. Locator is a trademark by Zest Anchors, Inc.
- The Screwdrivers for angled screw hole are medical devices manufactured by Talladium España S.L., Avenida Blondel, 54 3º, 25002 Lleida, Spin.
- The TWL torque wrench with control lever and its adapter are medical devices manufactured by Elos Medtech Pinol A/S, Engvej 33, 3330 Gorlose, Denmark.
- Caps for ball attachments, medical devices manufactured by Rhein83 S.r.l., Via E. Zago 10/ABC, 40128 Bologna, Italy.

They conform to the ISO 9001 and ISO 13485 standards and are certified with the CE in compliance with Regulation (EU) Medical Devices n. 2017/745.

The Sweden & Martina plant manufactures Medical Devices in compliance with the CGMPs in force in the USA and in other countries worldwide.



All trademarks herein are the property of Sweden & Martina S.p.A. unless otherwise indicated.
This material is intended for laboratories and clinicians and is not intended for patient distribution.
This material is not to be redistributed, duplicated, or disclosed without the express written consent of Sweden & Martina S.p.A.
For additional product information, including indications, contraindications, warnings, precautions, and potential adverse effects, see Sweden & Martina S.p.A. website.

The contents are updated at the time of publication. Check with the company for any subsequent updates.